



**STATE OF OKLAHOMA  
BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE**

Post Office Box 36307

Oklahoma City, Oklahoma 73136-2307  
Phone: (405) 425-2460 Fax: (405) 425-2490

www.bot.ok.gov

Dr. Jarrad Wagner, Ph.D. F-ABFT  
Chairman


Joshua Smith  
Director

**ATTESTATION**

I, the undersigned Joshua Smith (Director/Records Custodian), attest under oath pursuant to 12 O.S. §2902 (11)(a) this record was made, at or near the time of the occurrence of the matters set forth by or from information transmitted by a person having knowledge of those matters; was kept in the course of the regularly conducted business activity of the Board; and was made pursuant to the regularly conducted activity. The duplicate record provided is listed below:

- 1) The attached "Service Overview" is a true and accurate copy of maintenance performed on the Intoxilyzer 8000, serial number 80-002683, in compliance with the agency's Maintenance Policy and Procedures.

This packet consists of 16 pages.

Signed 

Name of Position: Director/Records Custodian

Date of Attestation: 07/13/2020

# Service Overview

Board of Tests for Alcohol and Drug Influence  
Enhancing public safety by ensuring the accuracy and scientific reliability  
of blood and breath testing



3600 N. Martin Luther King Bldg #9  
Oklahoma City, OK 73111  
405-425-2460

**BENCH CHECK DATE:** 07-13-2020

**SERIAL:** 80-002683

## **DOCUMENTATION WITHIN SERVICE PACKET:**

- BOT Certificate of Calibration and Operation
- Bench Check
- Bench Check Technician's Report
- BOT Dry Gas Certificate of Analysis
- Mock Subject Test
- BOT Certificates of Analysis on Certified Solutions (4)



Date:

Re: Request for Log of Test on Intoxilyzer 8000

Name:

Dear Sir:

The Log of Test information from a breath test conducted on a Intoxilyzer 8000 is printed in Section 2 of the Affidavit and is no longer retained as a separate form as with the Intoxilyzer 5000. The 8000 uses a nitrogen/ethanol dry gas for the external control test. The dry gas cylinder has the capacity to perform approximately 325 tests and does not require maintenance until its expiration date or can no longer provide the required pressure to perform a control test. The cylinder has an expiration date of two years from date filled at the manufacturer and is programmed in the 8000 upon installation, both dates are printed in Section 2 of the Affidavit. The Intoxilyzer is programmed to disable itself when the date of expiration is reached or the pressure monitor indicates 25 lbs psi.

Due to the reliability and stability of dry gas, maintenance on the 8000 is not required every 30 days or 25 tests, therefore maintenance is no longer performed by officers in the field. The Board of Tests for Alcohol and Drug Influence certifies each Intoxilyzer 8000 for field use by performing a bench check. The bench check is designed to check the operation and procedures programmed in the 8000 to ensure it is functioning properly for use in Oklahoma. A bench check is performed every time a new dry gas cylinder is installed or the instrument has returned from the manufacturer for repair. Once the instrument passes the bench check, all that is required to set it up for field use is to plug it in and turn it on and does not require a permit to set up. The Board of Tests has the ability to monitor each instrument once it is set up in the field.

Each Intoxilyzer has a database of valid permit numbers that it checks prior to allowing a breath test to be performed. The database is updated on a regular basis, if a permit number is invalid that officer will not be able to conduct a breath test on the 8000. If the 8000 prints an affidavit the breath test operator listed in Section 2 had a valid permit at the time of the test and the breath test is valid.

If you have any further questions please feel free to call me.

Sincerely,



J. Robert Blakeburn (OBA11059)  
State Director of Tests for Alcohol  
and Drug Influence

# Certificate of Calibration

This is to certify the calibration of **Intoxilyzer**® serial number 80-002683, manufactured by CMI Inc., a subsidiary of MPD, Inc. of Owensboro, Kentucky, was tested and found to conform to the National Highway Traffic Safety Administration (NHTSA) Standard for Devices to Measure Breath Alcohol (Federal Register, Vol.58 No.179, pp 48705-48710, Sept. 17, 1993) for accuracy and precision. Reference materials are traceable through the National Institute of Standards and Technology (NIST) to the International System of Units (SI).

Date July 6, 2020 Signed Alice Joler

Technician

**EMI** INC.

316 East 9<sup>th</sup> Street  
Owensboro, KY 42303  
USA

Part No. 650517 Rev.A

**STATE OF OKLAHOMA - BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE  
CMI INTOXILYZER BENCH CHECK PROCEDURE REPORT**

BENCH CHECK DATE 07/13/2020	TIME PROCEDURE BEGAN 14:08:0	INSTRUMENT LOCATION 80-002683
INSTRUMENT TYPE Intoxilyzer	INSTRUMENT MODEL I-8000	INSTRUMENT SERIAL NUMBER 80-002683

**TEST RESULTS**

**Diagnostics**

PASS

**Mouth Alcohol Test**

PASS

**RFI Test**

PASS

**Abort Test**

PASS

**Improper Sample Test**

PASS

**Interferent Detect Test**

PASS

**No Sample Given Test**

PASS

**Insufficient Sample Test**

PASS

**Air Blank Contamination Test**

PASS

**Range Exceeded Test**

PASS

**Duplicate Breath Sampling Results**

Value	Test 1	Test 2	Control Test
.02	0.020	0.018	0.019
.04	0.037	0.037	0.037
.10	0.097	0.096	0.097
.20	0.198	0.195	0.197

PASS

**CERTIFYING TECHNICIAN**

I swear under penalty of perjury that in regards to the above listed instrument, I have complied with the Intoxilyzer Bench Checking Procedure approved by the State of Oklahoma.

SIGNATURE 	DATE 7/13/2020
NAME PAPE, CHRISTOPHER	PERMIT # 0000023020



3600 N. Martin Luther King Bldg #9  
 Oklahoma City, OK 73111  
 405-425-2460

# BENCH CHECK TECHNICIAN'S REPORT

**Date:** 07/13/2020

**Start Time:** 14:00

**End Time:** 15:15

## INSTRUMENT

Intoxilyzer Model	Serial #
8000	80-002683

## GAS CANISTER

LOT #	EXP Date
14020080A1-001	07/05/2022

## REFERENCE

Simulator Model			
Guth 2100	Guth 2100	Guth 2100	Guth 2100
Serial #			
DR3591	DR3468	DR3667	DR3754
Concentration			
.02	.04	.10	.20
LOT #			
19050	19080	20190	20160
Manufactured Date			
02/04/2019	03/04/2019	04/06/2020	03/18/2020
Expiration Date			
02/04/2021	03/04/2021	04/06/2022	03/18/2022
Solution Commission Date			
7/7/2020	7/7/2020	7/7/2020	7/7/2020

## REASON FOR BENCH CHECK

- POST REPAIR   
  CYLINDER REPLACEMENT  
 TROUBLESHOOTING   
  ROUTINE MAINTENANCE

## COMMENTS


*Christopher Pape* 23020  
 Technician Signature / Permit #

Christopher Pape #23020  
 Printed Name

Certificate of Analysis - Dry Gas	
07/13/2020	0.080
Date of Analysis	Labelled target value (g/210L)
001	0.080
Cylinder #	Average test result
Christopher Pape <i>[Signature]</i> 23020	
BOT Technician Name and Signature	



8/16

ALPHA UNIT SUB  
Intoxilyzer - Alcohol Analyzer  
Model 8000 SN 80-003400  
07/13/2020 14:11

In accordance with BOT Policy and Procedure Statement BRT-2.2.0, Maintenance, the above referenced dry gas cylinder is suitable for use as an external reference with the Intoxilyzer 8000.  
7 Eastgate Dr., P.O. Box 790 • Jacksonville, IL 62651-0790  
217-245-2183 • Fax: 217-243-7634 • www.ilmoproducts.com

Test	g/210L	Time
Air Blank	0.000	14:11
Cal Check	0.080	14:12
Air Blank	0.000	14:12
Cal Check	0.081	14:12
Air Blank	0.000	14:13
Check	0.080	14:13
Blank	0.000	14:14
Check	0.080	14:14
Blank	0.000	14:15
Check	0.080	14:15
Blank	0.000	14:16
Check	0.080	14:16
Blank	0.000	14:16
Check	0.080	14:17
Blank	0.000	14:17
Check	0.080	14:18
Blank	0.000	14:18
Check	0.081	14:18
Blank	0.000	14:19
Check	0.080	14:19
Blank	0.000	14:20
Check Stats		
Average	0.0802	
Std Dev	0.0004	
Rel Std Dev(%)	0.5257	

## Certificate of Analysis

Certificate ID: 13020 ✓  
Part #: BAC105L080T  
Cylinder Size: 105L  
Lot Number: 14020080A1 ✓  
Expiration: 7/5/2022 ✓

*REV'D*  
*6/23/20*  
*CP*

0.080 BAC (For the calibration of instruments used to determine breath alcohol concentration)

Contents: 105 Liters @ 1000 psig 70°F (21°C)

Component:	Reported Concentration:	Analytical Accuracy (U, k=2):	Analytical Method:
Ethanol	208 ppm	+/-0.002 BAC(G/210L)	NDIR
Nitrogen	Balance	[5.2 ppm]	

Distributed by:  
CMI Inc.  
316 East Ninth Street  
Owensboro, KY 42303  
Phone 866-835-0690  
[www.alcoholtest.com](http://www.alcoholtest.com)

*[Signature]*  
Operator's Signature

\*Traceable to:  
Certified Reference Material - 262.4 µmol/mol  
Ethanol in Nitrogen - Serial No. GN0015026 Lot No. 050319E11

Store in dry area, away from sources of heat, ignition and direct sunlight. Do not allow storage area to exceed 52 °C (125 °F).

*[Signature]*  
Specialty Gas Lab Tech  
*[Signature]*

06-15-2020  
Issuance Date



The calibration results within this certificate were obtained using equipment and standards capable of producing analytical results traceable to NIST, and apply only to the items contained on this certificate. ILMO Products Company makes no warranty or representation as to the suitability of the use of any information provided for any particular purpose. The information use is at the sole discretion and risk of the user. Liability shall be limited to established replacement cost of this material or service.



**Service Evaluation Form**

**This form MUST be completed and enclosed with instrument to be serviced.  
Failure to complete and return this form may cause delays in service.**

**(Note: please ship items in their original shipping container(s) or a similar protective box.)**

1. Contact information: Customer Number **731362** (contact Customer Service)  
Name: **Board of Tests** Phone: **(405) 425 2460**  
Fax: **(405) 425 2490** Email: **christopher.pape@bot.ok.gov**

2. Bill to Address: **Board of Tests** Ship to Address: **Board of Tests**  
**3600 N ML King, Bldg 9** **3600 N ML King, Bldg 9**  
**OKC, OK 73111** **OKC, OK 73111**

3. Serial Number: **80-002683** Instrument Model: **Intoxilyzer 500 / 8000**

4. Detailed Description of Problem:  
**Flow sensor R value 100. Needs replaced.**

\*\*\*Hazardous Material Warning! - DO NOT return gas cylinders with instrument!\*\*\*

I Authorize Repairs Up To:  All  \$250  \$500  \$750  Other \$ \_\_\_\_\_  
Purchase Order Number (attach a copy of P.O. if applicable) \_\_\_\_\_

Authorized By:

Ship to:

**Chris Pape**

**Breath Testing Admin**

Name (Please Print)

Title

CMI, Inc



**06/23/20**

Attn: Service Dept.

Signature

Date

316 East Ninth Street

Call: 405-425-2460 or email

Owensboro, KY 42303

No, please send estimate before repairs are made.

Note: An estimate will be faxed before performing any repairs and may cause delays in service.  
An evaluation fee (\$79.00 or actual costs) will apply to estimates that are not repaired.



316 E 9th Street / Owensboro KY 42303 / USA  
 Phone: 866-835-0690 Fax: 270-685-6268

## SERVICE WORK ORDER

405072

DATE: 6/30/2020

**Bill To:**

Oklahoma Board Of Tests  
 Blvd Adcu Annex Bldg #9  
 3600 Martin Luther King  
 Oklahoma City, OK 73136-3000  
 USA  
 Customer #: 731362  
 Phone: 405-425-2424  
 Fax: 405-425-2490 fax

**Ship To:**

Oklahoma Board Of Tests  
 Blvd Adcu Annex Bldg #9  
 3600 Martin Luther King  
 Oklahoma City OK 73136-3000  
 USA  
 Phone: 405-425-2424

MODEL #: 002480OK

Serial Number: 80-002683

BILL CODE: Out of Warranty

**EXTRA PARTS RCVD:**

BT, SHELF

**DESCRIPTION OF PROBLEM**

LOW SENSOR R VALVE 100. NEEDS REPLACED.

**WORK PERFORMED:**

Flow sensor R: value 99. Breath hose o-ring missing. No other issues found. Replaced flow sensor and missing o-ring.

**PARTS USED**

Seq. No.	Part	Description	Quantity
50	402565	O-RING,NEOPRENE,,301IDx.06THK	1.00 EA
80	650517	CERTIFICATE OF CALIBRATION	2.00 EA
90	441169	COVER DUST,5/8in x 1/2in	1.00 EA
100	470145	CAP,PLASTIC,,25IDx.50LG,RED	3.00 EA
110	470154	CAP,PLSTC,,406IDx.438-.562,RED	1.00 EA
120	470218	SEAL,TAMPER EVIDENT,SERVICE	1.00 EA
130	210120	SENSOR,SMD,PRESSURE	1.00 EA

**LABOR / TESTING**

Misc Code	Description	Hours
LABR	Service Repair Labor	
LCAL	Service - Calibration Adjust	0.50
LFT	Service - Cal / Final Test	1.00

Technician Name: Gary Watts

SRV4050720001



# OFFICER'S AFFIDAVIT AND NOTICE OF REVOCATION/DISQUALIFICATION

11/16

Section 1

ARREST DATE 07/13/2020		ARREST TIME 1400		COUNTY # 55	CITY #	CITATION #		
ARREST LOCATION BOT				CITY OKLAHOMA CITY		COUNTY OKLAHOMA		
SUBJECT NAME NORTH, GAYLON LEVERNE				DATE OF BIRTH 08/09/1961	SEX M	RACE W	HEIGHT 510	WEIGHT 165
ADDRESS 3600 MARTIN LUTHER KING A			CITY OKLAHOMA CITY		STATE OK	ZIP CODE 73133		
DRIVER LICENSE # T083454871			EXPIRATION DATE 06/30/2019	STATE OK	CLASS A	RESTRICTIONS	ENDORSEMENTS M	
VEHICLE MAKE MOCK		MODEL MOCK	TAG # MOCK		STATE OK	CMV N	PLACARD N	

On the above date, time, and location, the above named person was arrested, and I had reasonable grounds to believe the person was driving or in actual physical control of a motor vehicle upon a public road, highway, street, turnpike, public parking lot, or other public place, or upon a private road, street, alley or lane which provides access to one or more single or multi-family dwelling within this state while under the influence of alcohol and/or other intoxicating substances as prohibited by law.  
(Describe driving behavior or circumstances):

Describe person's condition (odor, actions, etc.):

MOCK SUBJECT TEST, POST REPAIR BENCH CHECK, INSTALL NEW CYLINDER

THE PERSON WAS INFORMED OF THE IMPLIED CONSENT TEST REQUEST

BREATH-ALCOHOL ANALYSIS RECORD, REPORT OF THE PERSON NAMED IN SECTION 1 ABOVE, and LOG OF TEST AND MAINTENANCE RECORD.

Section 2

EQUIPMENT I-8000		SERIAL NUMBER 80-002683	SPECIALIST BOARD OF TESTS																																					
STANDARD Dry		TARGET VALUE 0.080	PERMIT NUMBER 23020																																					
MANUFACTURER ILMO		LOT # 14020080A1	DATE CYLINDER INSTALLED 07/13/2020	DATE CYLINDER EXPIRES 07/05/2022																																				
OBSERVATION START 1420	TEST DATE 07/13/2020	MOUTHPIECE G	<p style="text-align: center;">A 15-minute (or greater) deprivation/observation period was observed prior to and continuing to the time of test administration.</p> <p style="text-align: center;">I administered said Breath-Alcohol Test in compliance with the applicable operating procedure of the Board of Tests.</p> <p style="text-align: center;"><i>Christopher Pape</i> 23020 SIGNATURE OF OPERATOR</p>																																					
<table border="1" style="width: 100%; border-collapse: collapse;"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.00</td><td>14:56</td></tr> <tr><td>Subject Test</td><td>0.00</td><td>14:56</td></tr> <tr><td>Breath Volume</td><td>1.937 LITERS</td><td></td></tr> <tr><td>Air Blank</td><td>0.00</td><td>14:57</td></tr> <tr><td>Wait</td><td></td><td>14:59</td></tr> <tr><td>Air Blank</td><td>0.00</td><td>14:59</td></tr> <tr><td>Subject Test</td><td>0.00</td><td>14:59</td></tr> <tr><td>Breath Volume</td><td>2.128 LITERS</td><td></td></tr> <tr><td>Air Blank</td><td>0.00</td><td>15:00</td></tr> <tr><td>Cal Check</td><td>0.07</td><td>15:00</td></tr> <tr><td>Air Blank</td><td>0.00</td><td>15:01</td></tr> </tbody> </table>					Test	g/210L	Time	Air Blank	0.00	14:56	Subject Test	0.00	14:56	Breath Volume	1.937 LITERS		Air Blank	0.00	14:57	Wait		14:59	Air Blank	0.00	14:59	Subject Test	0.00	14:59	Breath Volume	2.128 LITERS		Air Blank	0.00	15:00	Cal Check	0.07	15:00	Air Blank	0.00	15:01
Test	g/210L	Time																																						
Air Blank	0.00	14:56																																						
Subject Test	0.00	14:56																																						
Breath Volume	1.937 LITERS																																							
Air Blank	0.00	14:57																																						
Wait		14:59																																						
Air Blank	0.00	14:59																																						
Subject Test	0.00	14:59																																						
Breath Volume	2.128 LITERS																																							
Air Blank	0.00	15:00																																						
Cal Check	0.07	15:00																																						
Air Blank	0.00	15:01																																						
<table border="1" style="width: 100%; border-collapse: collapse;"> <tr> <td colspan="2">NAME PAPE, CHRISTOPHER</td> </tr> <tr> <td>BADGE # 0000000484</td> <td>PERMIT # 0000023020</td> </tr> <tr> <td colspan="2">AGENCY BOARD OF TESTS</td> </tr> </table>			NAME PAPE, CHRISTOPHER		BADGE # 0000000484	PERMIT # 0000023020	AGENCY BOARD OF TESTS																																	
NAME PAPE, CHRISTOPHER																																								
BADGE # 0000000484	PERMIT # 0000023020																																							
AGENCY BOARD OF TESTS																																								

EFFECTIVE THIRTY (30) DAYS FROM DATE SERVED. Notice of Revocation is given that your driver's license is revoked or denied for a period of one hundred eighty (180) days or more, if you refused or failed the state's test(s). In addition, your commercial driving privilege in this state will be disqualified if you refused or failed the state's test(s) (operating a CMV or CDL holder) OR a test result of .04 or more while operating a commercial motor vehicle. You may appeal any Departmental action against your driver license as a result of this arrest by filing a petition in the District Court of the County in which you were arrested. Your petition must be filed within thirty (30) days of the date listed below. You may be eligible to participate in the Impaired Driver Accountability Program (IDAP) administered by the Department of Public Safety (DPS). Completion of IDAP may prevent a revocation as a result of this arrest from appearing on your driving record. Participation in IDAP may reduce the amount of fees you will be required to pay to the State. You must request IDAP within thirty (30) days of the date listed below. Mail the request for IDAP to the Department of Public Safety at P.O. Box 11415, Oklahoma City, OK 73136.

Section 3

**SERVING OF NOTICE:** I personally hand delivered a copy of the above containing the Notice of Revocation and the Notice of Department Impaired Driving Accountability Program to the person named in Section 1 above.

Date served: **07/13/2020** Signature of officer: \_\_\_\_\_ Badge #: \_\_\_\_\_

In accordance with Title 12 O.S. Section 426, "I state under penalty of perjury under the laws of Oklahoma that the foregoing is true and correct."

Section 4

Date: **07/13/2020** Place (location when signed): \_\_\_\_\_ Signature of arresting officer: \_\_\_\_\_

NAME PAPE, CHRISTOPHER		AGENCY BOARD OF TESTS	
BADGE 0000000484	AGENCY ADDRESS INTER-AGENCY	ZIP CODE 00000	PHONE 4054252460

OTHER WITNESSES:

Section 5

	NAME	TITLE	ADDRESS	PHONE
1				
2				

No Temporary License Issued:  
Test Below Legal Limit

# OFFICER'S AFFIDAVIT AND NOTICE OF REVOCATION/DISQUALIFICATION

12/16

Section 1

ARREST DATE <b>07/13/2020</b>		ARREST TIME <b>1400</b>		COUNTY # <b>55</b>	CITY #	CITATION #		
ARREST LOCATION <b>BOT</b>				CITY <b>OKLAHOMA CITY</b>		COUNTY <b>OKLAHOMA</b>		
SUBJECT NAME <b>NORTH, GAYLON LEVERNE</b>				DATE OF BIRTH <b>08/09/1961</b>	SEX <b>M</b>	RACE <b>W</b>	HEIGHT <b>510</b>	WEIGHT <b>165</b>
ADDRESS <b>3600 MARTIN LUTHER KING A</b>				CITY <b>OKLAHOMA CITY</b>		STATE <b>OK</b>	ZIP CODE <b>73133</b>	
DRIVER LICENSE # <b>T083454871</b>				EXPIRATION DATE <b>06/30/2019</b>	STATE <b>OK</b>	CLASS <b>A</b>	RESTRICTIONS	ENDORSEMENTS <b>M</b>
VEHICLE MAKE <b>MOCK</b>		MODEL <b>MOCK</b>		TAG # <b>MOCK</b>		STATE <b>OK</b>	CMV <b>N</b>	PLACARD <b>N</b>

On the above date, time, and location, the above named person was arrested, and I had reasonable grounds to believe the person was driving or in actual physical control of a motor vehicle upon a public road, highway, street, turnpike, public parking lot, or other public place, or upon a private road, street, alley or lane which provides access to one or more single or multi-family dwelling within this state while under the influence of alcohol and/or other intoxicating substances as prohibited by law.  
(Describe driving behavior or circumstances):

Describe person's condition (odor, actions, etc.):

**MOCK SUBJECT TEST, POST REPAIR BENCH CHECK, INSTALL NEW CYLINDER**

THE PERSON WAS INFORMED OF THE IMPLIED CONSENT TEST REQUEST

BREATH-ALCOHOL ANALYSIS RECORD, REPORT OF THE PERSON NAMED IN SECTION 1 ABOVE, and LOG OF TEST AND MAINTENANCE RECORD.

Section 2

EQUIPMENT <b>I-8000</b>		SERIAL NUMBER <b>80-002683</b>	SPECIALIST <b>BOARD OF TESTS</b>																																					
STANDARD <b>Dry</b>		TARGET VALUE <b>0.080</b>	PERMIT NUMBER <b>23020</b>																																					
MANUFACTURER <b>ILMO</b>		LOT # <b>14020080A1</b>	DATE CYLINDER INSTALLED <b>07/13/2020</b>	DATE CYLINDER EXPIRES <b>07/05/2022</b>																																				
OBSERVATION START <b>1420</b>	TEST DATE <b>07/13/2020</b>	MOUTHPIECE <b>G</b>	<p>A 15-minute (or greater) deprivation/observation period was observed prior to and continuing to the time of test administration.</p> <p>I administered said Breath-Alcohol Test in compliance with the applicable operating procedure of the Board of Tests.</p> <p style="text-align: center;"><i>[Signature]</i> <b>23020</b> SIGNATURE OF OPERATOR</p>																																					
<table border="1" style="width: 100%; border-collapse: collapse;"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.00</td><td>14:56</td></tr> <tr><td>Subject Test</td><td>0.00</td><td>14:56</td></tr> <tr><td>Breath Volume</td><td>1.937 LITERS</td><td></td></tr> <tr><td>Air Blank</td><td>0.00</td><td>14:57</td></tr> <tr><td>Wait</td><td></td><td>14:59</td></tr> <tr><td>Air Blank</td><td>0.00</td><td>14:59</td></tr> <tr><td>Subject Test</td><td>0.00</td><td>14:59</td></tr> <tr><td>Breath Volume</td><td>2.128 LITERS</td><td></td></tr> <tr><td>Air Blank</td><td>0.00</td><td>15:00</td></tr> <tr><td>Cal Check</td><td>0.07</td><td>15:00</td></tr> <tr><td>Air Blank</td><td>0.00</td><td>15:01</td></tr> </tbody> </table>					Test	g/210L	Time	Air Blank	0.00	14:56	Subject Test	0.00	14:56	Breath Volume	1.937 LITERS		Air Blank	0.00	14:57	Wait		14:59	Air Blank	0.00	14:59	Subject Test	0.00	14:59	Breath Volume	2.128 LITERS		Air Blank	0.00	15:00	Cal Check	0.07	15:00	Air Blank	0.00	15:01
Test	g/210L	Time																																						
Air Blank	0.00	14:56																																						
Subject Test	0.00	14:56																																						
Breath Volume	1.937 LITERS																																							
Air Blank	0.00	14:57																																						
Wait		14:59																																						
Air Blank	0.00	14:59																																						
Subject Test	0.00	14:59																																						
Breath Volume	2.128 LITERS																																							
Air Blank	0.00	15:00																																						
Cal Check	0.07	15:00																																						
Air Blank	0.00	15:01																																						
NAME <b>PAPE, CHRISTOPHER</b>			<div style="border: 1px solid black; padding: 5px;"> <p>Test Result BrAC g/210L <b>0.00</b></p> </div>																																					
BADGE # <b>0000000484</b>		PERMIT # <b>0000023020</b>																																						
AGENCY <b>BOARD OF TESTS</b>																																								

EFFECTIVE THIRTY (30) DAYS FROM DATE SERVED. Notice of Revocation is given that your driver's license is revoked or denied for a period of one hundred eighty (180) days or more, if you refused or failed the state's test(s). In addition, your commercial driving privilege in this state will be disqualified if you refused or failed the state's test(s) (operating a CMV or CDL holder) OR a test result of .04 or more while operating a commercial motor vehicle. You may appeal any Departmental action against your driver license as a result of this arrest by filing a petition in the District Court of the County in which you were arrested. Your petition must be filed within thirty (30) days of the date listed below. You may be eligible to participate in the Impaired Driver Accountability Program (IDAP) administered by the Department of Public Safety (DPS). Completion of IDAP may prevent a revocation as a result of this arrest from appearing on your driving record. Participation in IDAP may reduce the amount of fees you will be required to pay to the State. You must request IDAP within thirty (30) days of the date listed below. Mail the request for IDAP to the Department of Public Safety at P.O. Box 11415, Oklahoma City, OK 73136.

**SERVING OF NOTICE:** I personally hand delivered a copy of the above containing the Notice of Revocation and the Notice of Department Impaired Driving Accountability Program to the person named in Section 1 above.

Date served: **07/13/2020** Signature of officer: \_\_\_\_\_ Badge #: \_\_\_\_\_

In accordance with Title 12 O.S. Section 426, "I state under penalty of perjury under the laws of Oklahoma that the foregoing is true and correct." Section 4

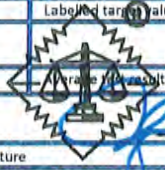
Date: **07/13/2020** Place (location when signed): \_\_\_\_\_ Signature of arresting officer: \_\_\_\_\_

NAME <b>PAPE, CHRISTOPHER</b>		AGENCY <b>BOARD OF TESTS</b>	
BADGE <b>0000000484</b>	AGENCY ADDRESS <b>INTER-AGENCY</b>	ZIP CODE <b>00000</b>	PHONE <b>4054252460</b>

No Temporary License Issued:  
Test Below Legal Limit

ARRESTEE'S COPY  
BOT Form AFF02 03/2020

**Certificate of Analysis - Wet Bath**

12-10-19	.020
Date of Analysis	Label target value (g/210L)
19050	.020
Lot Number	Average of Results (g/210L)
	
JOSH SMITH BOT Technician Name and Signature	



13/16

**GUTH LABORATORIES, INC.**  
 690 NORTH 67th STREET • HARRISBURG, PA 17111-4611 • TELEPHONE: 717-564-5470

BRT Form 6.0

In accordance with BRT - 2.2.0, Maintenance, the above referenced simulator solution is suitable for use as an external reference in maintenance of the Intoxilyzer 8000.

*Rec'd  
 JF  
 12/2019*

**CERTIFICATE OF ANALYSIS**

ALPHA INSTRUMENT 2  
 Intoxilyzer - Alcohol Analyzer  
 Model 8000  
 12/10/2019

SN 80-004923  
 14:52

**Certified Alcohol Reference Solution for Simulator**

Test	g/210L	Time
Air Blank	0.000	14:53
Cal Check	0.019	14:53
Air Blank	0.000	14:54
Cal Check	0.020	14:55
Air Blank	0.000	14:55
Cal Check	0.019	14:56
Air Blank	0.000	14:56
Cal Check	0.020	14:57
Air Blank	0.000	14:58
Cal Check	0.019	14:58
Air Blank	0.000	14:59
Cal Check	0.020	15:00
Air Blank	0.000	15:00
Cal Check	0.019	15:01
Air Blank	0.000	15:01
Cal Check	0.020	15:02
Air Blank	0.000	15:02
Cal Check	0.020	15:03
Air Blank	0.000	15:04
Cal Check	0.020	15:04
Air Blank	0.000	15:05
Cal Check Stats		
Average	0.0196	
Std Dev	0.0005	
Rel Std Dev(%)	2.6345	

Random Samples of Lot Number 19050 of Alcohol Reference Solution for Simulator were analyzed by gas chromatography on February 5, 2019, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain 0.0236% (w/vol) ethyl alcohol. The expiration date for this lot number is February 4, 2021 at 11:59 PM.

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of 0.020 g/210L +/- 3%.

The alcohol and water used in this solution were free of test interfering substances.

*JF*  
 10837  
 Operator's Signature


*Ted L. Pauley*  
 Ted L. Pauley, President  
 GUTH LABORATORIES, INC.

**NIST Traceability:**

Testing was conducted using Cerilliant Reference Standard lot number FN08031603 whose values are traceable to NIST.

All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.

**Certificate of Analysis – Wet Bath**

12-10-19	.040
Date of Analysis	Labelled tar(R) value (g/210L)
19080	.040
Lot Number	Alcohol Reference Solution (g/210L)
	



14/16

*Rec'd  
12/20/19*

**GUTH LABORATORIES, INC.**  
590 NORTH 67th STREET • HARRISBURG, PA 17111-4511 • TELEPHONE: 717-664-5470

JOSH SMITH  
BOT Technician Name and Signature BRT Form 6.0

In accordance with BRT – 2.2.0, Maintenance, the above referenced simulator solution is suitable for use as an external reference in maintenance of the Intoxilyzer 8000.

**CERTIFICATE OF ANALYSIS**

ALPHA INSTRUMENT 2  
Intoxilyzer - Alcohol Analyzer  
Model 8000 SN 80-004923  
12/10/2019 15:16

**Certified Alcohol Reference Solution for Simulator**

Test	g/210L	Time
Air Blank	0.000	15:17
Cal Check	0.040	15:17
Air Blank	0.000	15:18
Cal Check	0.040	15:19
Air Blank	0.000	15:19
Cal Check	0.040	15:20
Air Blank	0.000	15:20
Cal Check	0.040	15:21
Air Blank	0.000	15:22
Cal Check	0.040	15:22
Air Blank	0.000	15:23
Cal Check	0.040	15:24
Air Blank	0.000	15:24
Cal Check	0.040	15:25
Air Blank	0.000	15:25
Cal Check	0.039	15:26
Air Blank	0.000	15:27
Cal Check	0.039	15:27
Air Blank	0.000	15:28
Cal Check	0.039	15:29
Air Blank	0.000	15:29
Cal Check Stats		
Average	0.0397	
Std Dev	0.0004	
Rel Std Dev(%)	1.2167	

Random Samples of Lot Number **19080** of Alcohol Reference Solution for Simulator were analyzed by gas chromatography on **March 5, 2019**, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain **0.0481% (w/vol)** ethyl alcohol. The expiration date for this lot number is **March 4, 2021** at 11:59 PM.

When used in a calibrated Simulator, operating at **34°C +/- .2°C**, this solution will give a breath alcohol analysis instrument reading of **0.040 g/210L +/- 3%**.

The alcohol and water used in this solution were free of test interfering substances.

*[Signature]*  
16837  
Operator's Signature

*[Signature]*  
Ted L. Pauley, President  
GUTH LABORATORIES, INC.

**NIST Traceability:**  
Testing was conducted using Cerilliant Reference Standard lot number FN12181501 whose values are traceable to NIST.  
All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.



Date of Analysis	06/16/20	Labelled target value (g/210L)	.100
Lot Number	20190	Average test result (g/210L)	.098
BOT Technician Name and Signature Josh Smith			

**GUTH LABORATORIES, INC.**  
 BRT Form 6.0 RISE JRG, PA 17111-4511 • TELEPHONE: 717-564-5470

In accordance with BRT - 2.2.0, Maintenance, the above referenced simulator solution is suitable for use as an external reference in maintenance of the Intoxilyzer 8000.

Rec'd  
06/11/2020  
*[Signature]*

**CERTIFICATE OF ANALYSIS**

ALPHA INSTRUMENT  
 Intoxilyzer - Alcohol Analyzer  
 Model 8000  
 06/16/2020

SN 80-002591  
 16:18

Certified Alcohol Reference Solution for Simulator

Test	g/210L	Time
Air Blank	0.000	16:18
Cal Check	0.099	16:19
Air Blank	0.000	16:20
Cal Check	0.098	16:20
Air Blank	0.000	16:21
Cal Check	0.097	16:22
Air Blank	0.000	16:22
Cal Check	0.098	16:23
Air Blank	0.000	16:24
Cal Check	0.098	16:24
Air Blank	0.000	16:25
Cal Check	0.098	16:26
Air Blank	0.000	16:26
Cal Check	0.098	16:27
Air Blank	0.000	16:27
Cal Check	0.097	16:28
Air Blank	0.000	16:29
Cal Check	0.098	16:29
Air Blank	0.000	16:30
Cal Check	0.098	16:31
Air Blank	0.000	16:31
Cal Check Stats		
Average	0.0979	
Std Dev	0.0005	
Rel Std Dev(%)	0.5798	

Random Samples of Lot Number 20190 of Alcohol Reference Solution for Simulator were analyzed by gas chromatography on April 8, 2020, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain 0.1206% (w/vol) ethyl alcohol. The expiration date for this lot number is April 6, 2022 at 11:59 PM.

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of 0.100 g/210L +/- 3%.

The alcohol and water used in this solution were free of test interfering substances.

*[Signature]*  
 Operator's Signature 16837

*[Signature]*  
 Ted L. Pauley, President  
 GUTH LABORATORIES, INC.

**NIST Traceability:**  
 Testing was conducted using Cerilliant Reference Standard lot number FN02271802 whose values are traceable to NIST.  
 All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.

Certificate of Analysis - Wet Bath	
07/07/2020	0.200
Date of Analysis	Labelled Target Value (g/210L)
20160	0.197
Lot Number	Average Test Result (g/210L)
Christopher Pope	0.23030
ROT Technician Name and Signature	



**GUTH LABORATORIES, INC.**  
 690 NORTH 67th STREET • HARRISBURG, PA 17111-4511 • TELEPHONE: 717-664-5470

BRT Form 6.0

In accordance with BRT - 2.2.0, Maintenance, the above referenced simulator solution is suitable for use as an external reference in maintenance of the Intoxilyzer 8000.

### CERTIFICATE OF ANALYSIS

*RECD  
05/2020  
CP*

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number **20160** of Alcohol Reference Solution for Simulator were analyzed by gas chromatography on March 20, 2020, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain **0.2427%** (w/vol) ethyl alcohol. The expiration date for this lot number is March 18, 2022 at 11:59 PM.

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of 0.200 g/210L +/- 3%.

The alcohol and water used in this solution were free of test interfering substances.

LPRA UNIT SUB  
 Intoxilyzer - Alcohol Analyzer  
 Model 8000 SN 80-003400  
 7/07/2020 15:25

Test	g/210L	Time
ir Blank	0.000	15:26
al Check	0.196	15:26
ir Blank	0.000	15:27
al Check	0.196	15:28
ir Blank	0.000	15:28
al Check	0.197	15:29
ir Blank	0.000	15:29
al Check	0.196	15:30
ir Blank	0.000	15:31
al Check	0.198	15:31
ir Blank	0.000	15:32
al Check	0.198	15:33
ir Blank	0.000	15:33
al Check	0.197	15:34
ir Blank	0.000	15:35
al Check	0.198	15:35
ir Blank	0.000	15:36
al Check	0.198	15:37
ir Blank	0.000	15:37
al Check	0.197	15:38
ir Blank	0.000	15:38
al Check Stats		
Average	0.1971	
Std Dev	0.0008	
Rel Std Dev(%)	0.4442	

*Christopher Pope*  
 Operator's Signature

*Ted L. Pauley*  
 Ted L. Pauley, President  
 GUTH LABORATORIES, INC.

**NIST Traceability:**

Testing was conducted using Cerilliant Reference Standard lot number FN08101505 whose values are traceable to NIST.

All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.