



**STATE OF OKLAHOMA
BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE**

Dr. Jarrad Wagner, Ph.D. F-ABFT
Chairman

Post Office Box 36307

Oklahoma City, Oklahoma 73136-2307

Phone: (405) 425-2460 Fax: (405) 425-2490

www.bot.ok.gov

Joshua Smith
Director

ATTESTATION

I, the undersigned Joshua Smith (Director/Records Custodian), attest under oath pursuant to 12 O.S. §2902 (11)(a) this record was made, at or near the time of the occurrence of the matters set forth by or from information transmitted by a person having knowledge of those matters; was kept in the course of the regularly conducted business activity of the Board; and was made pursuant to the regularly conducted activity. The duplicate record provided is listed below:

- 1) The attached "Service Overview" is a true and accurate copy of maintenance performed on the Intoxilyzer 8000, serial number 80-002696, in compliance with the agency's Maintenance Policy and Procedures.

Signed _____
JS

Name of Position: Director/Records Custodian

Date of Attestation: 05/12/2023

AMENDED 09/19/2023



STATE OF OKLAHOMA
BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE

Dr. Jarrad Wagner, Ph.D., F-ABFT
Chairman

Post Office Box 36307
Oklahoma City, Oklahoma 73136-2307
Phone: (405) 425-2460 Fax: (405) 425-2490
www.bot.ok.gov

Joshua Smith
Director

9/19/2023

It was brought to our attention of a clerical error with the documentation of the 0.04 reference standard in the maintenance packet. The maintenance packet has been amended to include the correct documentation.

Christopher Pape
Breath Alcohol Testing Program Administrator



STATE OF OKLAHOMA
BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE

Dr. Jarrad Wagner, Ph.D., F-ABFT
 Chairman

Post Office Box 36307
 Oklahoma City, Oklahoma 73136-2307
 Phone: (405) 425-2460 Fax: (405) 425-2490
 www.bot.ok.gov

Joshua Smith
 Director

Certificate of Analysis – Reference Solution


ALPHA INSTRUMENT
 Intoxilyzer - Alcohol Analyzer
 Model 8000
 10/24/2022

SN 80-002591
 15:45

| Test | g/210L | Time |
|-----------------|--------|-------|
| Air Blank | 0.000 | 15:45 |
| Cal Check | 0.039 | 15:46 |
| Air Blank | 0.000 | 15:47 |
| Cal Check | 0.039 | 15:47 |
| Air Blank | 0.000 | 15:48 |
| Cal Check | 0.039 | 15:49 |
| Air Blank | 0.000 | 15:49 |
| Cal Check | 0.039 | 15:50 |
| Air Blank | 0.000 | 15:50 |
| Cal Check | 0.039 | 15:51 |
| Air Blank | 0.000 | 15:52 |
| Cal Check | 0.039 | 15:52 |
| Air Blank | 0.000 | 15:53 |
| Cal Check | 0.039 | 15:54 |
| Air Blank | 0.000 | 15:54 |
| Cal Check | 0.039 | 15:55 |
| Air Blank | 0.000 | 15:56 |
| Cal Check | 0.039 | 15:56 |
| Air Blank | 0.000 | 15:57 |
| Cal Check | 0.040 | 15:57 |
| Air Blank | 0.000 | 15:58 |
| Cal Check Stats | | |
| Average | 0.0391 | |
| Std Dev | 0.0003 | |
| Rel Std Dev(%) | 0.8087 | |

NIST Traceable Breath-Alcohol Equivalent (g/210L +/- 3%) Ratio of ethanol and water free of test interfering substances.

| | | |
|--|------------------------------|----------------------------|
| Manufacturer: Guth | Manufactured Date: 10/4/2021 | Concentration: 0.040 |
| Lot No.: 21410 | Expiration Date: 9/30/2023 | Equipment Model: Guth 2100 |
| The above solution has been analyzed via Infrared Spectroscopy utilizing the Intoxilyzer 8000 and was performed in accordance with the procedures prescribed by the State Director of Tests pursuant to 40 OAC 30-1-3. The above NIST traceable reference solution meets the requirements for use pursuant to 40 OAC 25-1-2.1. Approved breath-alcohol measurement equipment and reference solutions. Equipment used appears on the current or supplemented Conforming Products List of Calibrating Units for Breath Alcohol Testers (72 FR 34747), published by the National Highway Traffic Safety Administration. | | |
| Signature <i>Christopher Pape</i> | Date: 10/24/2022 | |
| Name Christopher Pape | Permit No.: 23020 | |

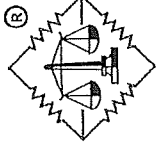

GUTH LABORATORIES, INC.

590 NORTH 67th STREET • HARRISBURG, PA 17111-4511 • TELEPHONE: 717-564-5470

RCVD
5/9/2022
CP

CERTIFICATE OF ANALYSIS

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number **21410** of Alcohol Reference Solution for Simulator were analyzed by gas chromatography on **October 4, 2021**, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain **0.0485%** (w/vol) ethyl alcohol. The expiration date for this lot number is **September 30, 2023** at 11:59 PM.

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of **0.040 g/210L +/- 3%**.

The alcohol and water used in this solution were free of test interfering substances.

Ted L. Pauley, President
GUTH LABORATORIES, INC.

NIST Traceability:
Testing was conducted using Cerilliant Reference Standard lot number **FN08211802** whose values are traceable to NIST.
All balances are calibrated annually by an outside agency using NIST traceable weights.
Calibration verification is done prior to each use utilizing NIST traceable weights.

Service Overview

5/20



Board of Tests for Alcohol and Drug Influence
Enhancing public safety by ensuring the accuracy and scientific reliability
of blood and breath testing.

3600 N. Martin Luther King Bldg #9
Oklahoma City, OK 73111
405-425-2460

BENCH CHECK DATE: 5/12/2023

SERIAL: 80-002696

DOCUMENTATION WITHIN SERVICE PACKET:

- BOT Certificate of Maintenance, Calibration and Operation
- Bench Check
- BOT Dry Gas Certificate of Analysis
- Ilmo Specialty Gases Certificate of Analysis
- Mock Subject Test
- BOT Certificate of Analysis on Certified Solutions (4)
- Guth Laboratories Certificate of Analysis (4)



6/20

STATE OF OKLAHOMA BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE

Post Office Box 36307
Oklahoma City, Oklahoma 73136-2307
Phone: (405) 425-2460 Fax: (405) 425-2490
www.bot.ok.gov

CERTIFICATE of MAINTENANCE, CALIBRATION, and OPERATION

This certifies the Intoxilyzer 8000 referenced below, manufactured by CMI, Inc., a subsidiary of MPD, Inc. of Owensboro, Kentucky, is maintained by the Board of Tests for Alcohol and Drug Influence in accordance with 40 O.A.C. 30-1-3. The instrument reported below is hereby deemed to be in compliance with Board rules via linearity accuracy examination and operational/functional procedure evaluation also referred to as a "Bench Check". The Intoxilyzer 8000 appears on the National Highway Traffic Safety Administration's Conforming Products List of Evidential Breath Alcohol Measurement Devices (82 FR 50940 *et seq*) as an approved device. In accordance with 40 O.A.C. 25-1-2, this instrument is approved and certified for use in the State of Oklahoma as an automated evidential breath-alcohol measurement device by the Board of Tests pursuant to 47 O.S. § 759.

All analyses conducted by this instrument are considered precise, valid, and in accordance with 40 O.A.C. 30-1-3(b) Operating Procedure(s). A printed breath test operator affidavit containing the seal of the Board by this instrument confirms a calibration check successfully verified the instrument's calibration at the time of test. The printed breath test affidavit also serves as proof that safeguards in place to protect the breath samples, including but not limited to Ambient Air Interference, Radio Frequency Interference, Residual Mouth Alcohol, Deficient Samples, etc. were not triggered and the test is a true and accurate scientific measurement of the person's alcohol concentration.

- Maintenance was performed on the Intoxilyzer 8000, Serial No. 80-002696 on 5/12/2023 in accordance with 40 O.A.C. 30-1-3 and internal BOT maintenance policy and procedures.
- A dry gas canister appearing on the National Highway Traffic Safety Administration's Conforming Products List of Calibrating Units for Breath Alcohol Testers (77 FR 64588 et seq) and approved in accordance with 40 O.A.C. 25-1-3 and BOT maintenance policy was installed.
 - Manufacturer: ILMO Unit: 208 ppm / 105 L
 - Cylinder Lot No. - Canister No.: 06723080A4-036
 - BOT analysis average test result: 0.081
 - Cylinder expiration date: 4/5/2025
- In accordance with 40 O.A.C. 25-1-2.1, the approved reference standards and equipment were used:

| | | | | |
|---------------------------------|-----------|-----------|-----------|-----------|
| Simulator Model | Guth 2100 | Guth 2100 | Guth 2100 | Guth 2100 |
| Serial # | | | | |
| | DR2834 | DR3694 | DR2455 | DR3470 |
| Concentration | .02 | .04 | .10 | .20 |
| LOT # | 21340 | 21410 | 22080 | 22050 |
| Manufactured Date | 8/24/2021 | 9/30/2021 | 3/7/2022 | 2/7/2022 |
| Expiration Date | 8/24/2023 | 9/30/2023 | 3/7/2024 | 2/7/2024 |
| Solution Commission Date | 5/12/2023 | 5/12/2023 | 5/12/2023 | 5/12/2023 |

Completed by BOT Personnel:


Signature/Name Christopher Pape

Breath-Alcohol Testing Program Administrator
Title


Certified by: Joshua Smith
Signature Name

State Director
Title

**STATE OF OKLAHOMA - BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE
CMI INTOXILYZER BENCH CHECK PROCEDURE REPORT**

| | | |
|--------------------------------|---------------------------------|---------------------------------------|
| BENCH CHECK DATE 05/12/2023 | TIME PROCEDURE BEGAN 12:16:4 | INSTRUMENT LOCATION 80-002696 |
| INSTRUMENT TYPE Intoxilyzer | INSTRUMENT MODEL I-8000 | INSTRUMENT SERIAL NUMBER 80-002696 |

TEST RESULTS

Diagnostics
PASS

Mouth Alcohol Test
PASS

RFI Test
PASS

Abort Test
PASS

Improper Sample Test
PASS

Interferent Detect Test
PASS

No Sample Given Test
PASS

Insufficient Sample Test
PASS

Air Blank Contamination Test
PASS

Range Exceeded Test
PASS


Duplicate Breath Sampling Results

| Value | Test 1 | Test 2 | Control Test |
|-------|--------|--------|--------------|
| .02 | 0.017 | 0.017 | 0.018 |
| .04 | 0.037 | 0.036 | 0.038 |
| .10 | 0.096 | 0.096 | 0.098 |
| .20 | 0.191 | 0.191 | 0.193 |

PASS

CERTIFYING TECHNICIAN

I swear under penalty of perjury that in regards to the above listed instrument, I have complied with the Intoxilyzer Bench Checking Procedure approved by the State of Oklahoma.

| | |
|---|------------------------|
| SIGNATURE  | DATE 5/12/23 |
| NAME PAPE, CHRISTOPHER | PERMIT # 0000023020 |



STATE OF OKLAHOMA
BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE

Post Office Box 36307

Dr. Jarrad Wagner, Ph.D., F-ABFT
 Chairman

Oklahoma City, Oklahoma 73136-2307
 Phone: (405) 425-2460 Fax: (405) 425-2490
 www.bot.ok.gov


Joshua Smith
 Director

Certificate of Analysis – Dry Gas

ALPHA INSTRUMENT

Intoxilyzer - Alcohol Analyzer SN 80-002591
 Model 8000
 04/21/2023 13:51

| Test | g/210L | Time |
|-----------------|--------|-------|
| Air Blank | 0.000 | 13:52 |
| Cal Check | 0.081 | 13:52 |
| Air Blank | 0.000 | 13:53 |
| Cal Check | 0.081 | 13:53 |
| Air Blank | 0.000 | 13:54 |
| Cal Check | 0.081 | 13:54 |
| Air Blank | 0.000 | 13:54 |
| Cal Check | 0.082 | 13:55 |
| Air Blank | 0.000 | 13:55 |
| Cal Check | 0.081 | 13:56 |
| Air Blank | 0.000 | 13:56 |
| Cal Check | 0.081 | 13:56 |
| Air Blank | 0.000 | 13:57 |
| Cal Check | 0.081 | 13:57 |
| Air Blank | 0.000 | 13:58 |
| Cal Check | 0.081 | 13:58 |
| Air Blank | 0.000 | 13:59 |
| Cal Check | 0.081 | 13:59 |
| Air Blank | 0.000 | 14:00 |
| Cal Check | 0.000 | 14:00 |
| Cal Check Stats | | |
| Average | 0.0811 | |
| Std Dev | 0.0003 | |
| Rel Std Dev(%) | 0.3899 | |

| Known Breath-Alcohol Equivalent (g/210L) Ratio of Nitrogen and Ethanol Gas (208 ppm/105 L) $0.08 \text{ g/210L} = 208 \text{ PPM} \div 2605^*$ | |
|---|--------------------------------------|
| Manufacturer: ILMO | Manufacturer Cert. No.: 15566 |
| Lot No.: 06723080A4 | Cylinder No.: 036 |
| The above 105 L cylinder has been analyzed via Infrared Spectroscopy utilizing the Intoxilyzer 8000 and was performed in accordance with the procedures prescribed by the State Director of Tests pursuant to 40 OAC 30-1-3. The above NIST traceable reference standard meets the requirements for use as a calibrating unit in accordance with 47 O.S. 759, 40 OAC 25-1-3, and appears on the current conforming products list of calibrating units for breath alcohol testers published by the U.S. Department of Transportation in the Federal Register [77 FR 64588 et seq]. * See Dubowski, K.M. and Essary N., <i>Journal of Anal. Toxicol.</i> (1996) 20(6): pp. 484-491. [g/210L = PPM \div 2605]. | |
| Signature  | Date: 04/21/2023 |
| Name Christopher Pape | Permit No.: 23020 |



7 Eastgate Dr. • P.O. Box 790 • Jacksonville, IL 62651-0790
217-245-2183 • Fax: 217-243-7634 • www.ilmoproducts.com

Certificate of Analysis

Certificate ID: 15566
Part #: BAC105L080T
Cylinder Size: 105L
Lot Number: 06723080A4
Expiration: 4/5/2025

RQVD
4/6/23
CP

0.080 BAC (For the calibration of instruments used to determine breath alcohol concentration)

Contents: 105 Liters @ 1000 psig 70°F (21°C)

| Component: | Reported Concentration: | Analytical Accuracy | |
|------------|-------------------------|----------------------|---------|
| | | (U, k=2): | Method: |
| Ethanol | 208 ppm | +/-0.002 BAC(g/210L) | NDIR |
| Nitrogen | Balance | [5.2 ppm] | |


Distributed by:
CMI Inc.
316 East Ninth Street
Owensboro, KY 42303
Phone 866-835-0690
www.alcoholtest.com

*Traceable to:

Certified Reference Material - 261.0 µmol/mol

Ethanol in Nitrogen - Serial No. ND7017 Lot No. 080722E2

Store in dry area, away from sources of heat, ignition and direct sunlight. Do not allow storage area to exceed 52 °C (125 °F).



Specialty Gas Lab Tech

03-21-2023

Issuance Date



The calibration results within this certificate were obtained at the facility listed above using equipment and standards capable of producing analytical results traceable to NIST, and apply only to the items contained on this certificate. ILMO Products Company makes no warranty or representation as to the suitability of the use of any information provided for any particular purpose. The information use is at the sole discretion and risk of the user. Liability shall be limited to established replacement cost of this material or service. This certificate applies only to the items described and shall not be reproduced other than in full, without written approval from the issuing facility.


J. Banovic

OFFICER'S IMPAIRED DRIVING AFFIDAVIT

10/20

| | | | | | | |
|-----------------------------------|--|-------------------------------|---------------|--------------------|----------------------------|----------|
| ARREST DATE 05/12/2023 | | COUNTY # 55 | CITY # | CITATION # | Section 1 | |
| ARREST LOCATION BOT | | CITY | | COUNTY OKLAHOMA | | |
| SUBJECT NAME BREATHTEST, DELTA | | DATE OF BIRTH 01/01/1985 | SEX M | HEIGHT 510 | WEIGHT 200 | |
| ADDRESS 3600 MLK | | CITY OKLAHOMA CITY | STATE OK | ZIP CODE 73111 | RESTRICTIONS /ENDORSEMENTS | |
| DRIVER LICENSE # B083463042 | | EXPIRATION DATE 10/31/2022 | STATE OK | CLASS D | STATE OK | CMV N |
| VEHICLE MAKE MOCK | | MODEL MOCK | TAG # MOCK | PLACARD N | N | |

On the above date, time, and location, the above named person was arrested, and I had reasonable grounds to believe the person was driving or in actual physical control of a motor vehicle upon a public road, highway, street, turnpike, public parking lot, or other public place, or upon a private road, street, alley or lane which provides access to one or more single or multi-family dwelling within this state while under the influence of alcohol and/or other intoxicating substances as prohibited by law.
(Describe driving behavior or circumstances):

MOCK SUBJECT TEST

Describe person's condition (odor, actions, etc.):

ROUTINE MAINTENANCE, INSTALLED NEW CYLINDER

THE PERSON WAS INFORMED OF THE IMPLIED CONSENT TEST REQUEST

| | | | | | |
|----------------------|--|----------------------------|---------------------------------------|-------------------------------------|--|
| EQUIPMENT I-8000 | | SERIAL NUMBER 80-002696 | SPECIALIST BOARD OF TESTS | Section 2 | |
| STANDARD Dry | | TARGET VALUE 0.080 | PERMIT NUMBER 23020 | | |
| MANUFACTURER ILMO | | LOT # 06723080A4 | DATE CYLINDER INSTALLED 05/12/2023 | DATE CYLINDER EXPIRES 04/05/2025 | |

OBSERVATION START TEST DATE
1305 05/12/2023

RESULT TABLE

| Test | g/210L | Time |
|---------------|--------------|-------|
| Air Blank | 0.00 | 13:41 |
| Subject Test | 0.00 | 13:42 |
| Breath Volume | 2.128 LITERS | |
| Air Blank | 0.00 | 13:42 |
| Wait | | 13:44 |
| Air Blank | 0.00 | 13:44 |
| Subject Test | 0.00 | 13:45 |
| Breath Volume | 2.140 LITERS | |
| Air Blank | 0.00 | 13:45 |
| Cal Check | 0.08 | 13:45 |
| Air Blank | 0.00 | 13:46 |

A 15-minute (or greater) deprivation/observation period was observed prior to and continuing to the time of test administration.

I administered said Breath-Alcohol Test in compliance with the applicable operating procedure of the Board of Tests.

[Signature]

SIGNATURE OF OPERATOR

| | |
|---------------------------|-----------------------|
| NAME PAPE, CHRISTOPHER | PERMIT # 000000484 |
| BADGE # 000000484 | 0000023020 |
| AGENCY BOARD OF TESTS | |



| |
|---|
| Test Result BrAC g/210L 0.00 |
|---|

What happens next?

- The Arresting Officer is required to submit the Officer's Impaired Driving Affidavit to the State of Oklahoma.
- Upon receipt and review of the Officer's Impaired Driving Affidavit and any associated test results, if applicable, the State of Oklahoma may revoke your Oklahoma driving privileges. If you hold a commercial driver license or were driving a commercial motor vehicle, the State of Oklahoma may disqualify your privileges to operate a commercial motor vehicle. You will continue to have driving privileges, if otherwise eligible, unless or until an Order of Revocation and/or Disqualification is issued by the State of Oklahoma.
- The State of Oklahoma will send you the Order of Revocation and/or Disqualification to the address you have on file with Service Oklahoma. To confirm or change your address on file, contact Service Oklahoma.
- The Revocation and/or Disqualification will be effective 55 days from the date of the Order, including 10 days allowed for mailing, 47 O.S. §§2-116, 753, 754.
- The Order of Revocation and/or Disqualification will contain important information about your rights to appeal the revocation and/or disqualification of your driving privileges. To appeal the applicable Order of Revocation and/or Disqualification your petition must be filed with the District Court of the County in which you were arrested within 40 days of the date of the order, including 10 days allowed for mailing, 47 O.S. §§2-116, 6-211.

In accordance with Title 12 O.S. Section 426, "I state under penalty of perjury under the laws of Oklahoma that the foregoing is true and correct." Section 4

Date: _____ Place (location when signed): _____ Signature of arresting officer: _____

| | |
|---------------------------|--------------------------|
| NAME PAPE, CHRISTOPHER | AGENCY BOARD OF TESTS |
| BADGE 000000484 | INTER-AGENCY |
| AGENCY ADDRESS | ZIP CODE 00000 |
| PHONE 4054252460 | PHONE 4054252460 |

OTHER WITNESSES:

| | | | |
|------|-------|---------|-------|
| NAME | TITLE | ADDRESS | PHONE |
| 1 | | | |
| 2 | | | |

OFFICER'S IMPAIRED DRIVING AFFIDAVIT

11/20

Section 1

| | | | | | |
|-----------------------------------|-------------------------------|-----------------------------|--------------|-------------------|-----------------------------|
| ARREST DATE 05/12/2023 | ARREST TIME 1300 | COUNTY # 55 | CITY # | CITATION # | COUNTY OKLAHOMA |
| ARREST LOCATION BOT | | | | | |
| SUBJECT NAME BREATHTEST, DELTA | | DATE OF BIRTH 01/01/1985 | SEX M | HEIGHT 510 | WEIGHT 200 |
| ADDRESS 3600 MLK | | CITY OKLAHOMA CITY | STATE OK | ZIP CODE 73111 | RESTRICTIONS / ENDORSEMENTS |
| DRIVER LICENSE # B083463042 | EXPIRATION DATE 10/31/2022 | STATE OK | CLASS D | STATE OK | CMV N |
| VEHICLE MAKE MOCK | MODEL MOCK | TAG # MOCK | PLACARD N | | |

On the above date, time, and location, the above named person was arrested, and I had reasonable grounds to believe the person was driving or in actual physical control of a motor vehicle upon a public road, highway, street, turnpike, public parking lot, or other public place, or upon a private road, street, alley or lane which provides access to one or more single or multi-family dwelling within this state while under the influence of alcohol and/or other intoxicating substances as prohibited by law.

(Describe driving behavior or circumstances):
MOCK SUBJECT TEST

Describe person's condition (odor, actions, etc.):

ROUTINE MAINTENANCE, INSTALLED NEW CYLINDER

THE PERSON WAS INFORMED OF THE IMPLIED CONSENT TEST REQUEST
BREATH-ALCOHOL ANALYSIS RECORD, REPORT OF THE PERSON NAMED IN SECTION 1 ABOVE, and LOG OF TEST AND MAINTENANCE RECORD.

Section 2

| EQUIPMENT I-8000 | SERIAL NUMBER 80-002696 | SPECIALIST BOARD OF TESTS | TARGET VALUE 0.080 | PERMIT NUMBER 23020 | DATE CYLINDER EXPIRES 04/05/2025 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
|---|----------------------------|------------------------------|-----------------------|---------------------------------------|-------------------------------------|------|------|--------|-----------|------|-------|--------------|------|-------|---------------|--------------|--|-----------|------|-------|------|--|-------|-----------|------|-------|--------------|------|-------|---------------|--------------|--|-----------|------|-------|-----------|------|-------|-----------|------|-------|
| STANDARD Dry | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| MANUFACTURER ILMO | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| OBSERVATION START 1305 | | TEST DATE 05/12/2023 | | DATE CYLINDER INSTALLED 05/12/2023 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| <p>RESULT TABLE</p> <table style="width: 100%; border-collapse: collapse;"> <thead> <tr> <th style="width: 60%;">Test</th> <th style="width: 10%;">Time</th> <th style="width: 30%;">Result</th> </tr> </thead> <tbody> <tr> <td>Air Blank</td> <td>0.00</td> <td>13:41</td> </tr> <tr> <td>Subject Test</td> <td>0.00</td> <td>13:42</td> </tr> <tr> <td>Breath Volume</td> <td>2.128 LITERS</td> <td></td> </tr> <tr> <td>Air Blank</td> <td>0.00</td> <td>13:42</td> </tr> <tr> <td>Wait</td> <td></td> <td>13:44</td> </tr> <tr> <td>Air Blank</td> <td>0.00</td> <td>13:44</td> </tr> <tr> <td>Subject Test</td> <td>0.00</td> <td>13:45</td> </tr> <tr> <td>Breath Volume</td> <td>2.140 LITERS</td> <td></td> </tr> <tr> <td>Air Blank</td> <td>0.00</td> <td>13:45</td> </tr> <tr> <td>Cal Check</td> <td>0.08</td> <td>13:45</td> </tr> <tr> <td>Air Blank</td> <td>0.00</td> <td>13:46</td> </tr> </tbody> </table> <p>A 15-minute (or greater) deprivation/observation period was observed prior to and continuing to the time of test administration.</p> <p>I administered said Breath-Alcohol Test in compliance with the applicable operating procedure of the Board of Tests.</p> <p style="text-align: center;"><i>[Signature]</i> SIGNATURE OF OPERATOR</p> | | | | | | Test | Time | Result | Air Blank | 0.00 | 13:41 | Subject Test | 0.00 | 13:42 | Breath Volume | 2.128 LITERS | | Air Blank | 0.00 | 13:42 | Wait | | 13:44 | Air Blank | 0.00 | 13:44 | Subject Test | 0.00 | 13:45 | Breath Volume | 2.140 LITERS | | Air Blank | 0.00 | 13:45 | Cal Check | 0.08 | 13:45 | Air Blank | 0.00 | 13:46 |
| Test | Time | Result | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| Air Blank | 0.00 | 13:41 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| Subject Test | 0.00 | 13:42 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| Breath Volume | 2.128 LITERS | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| Air Blank | 0.00 | 13:42 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| Wait | | 13:44 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| Air Blank | 0.00 | 13:44 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| Subject Test | 0.00 | 13:45 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| Breath Volume | 2.140 LITERS | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| Air Blank | 0.00 | 13:45 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| Cal Check | 0.08 | 13:45 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| Air Blank | 0.00 | 13:46 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| NAME PAPE, CHRISTOPHER | | PERMIT # 000000484 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| BADGE # 000000484 | | AGENCY BOARD OF TESTS | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |
| AGENCY BOARD OF TESTS | | PERMIT # 0000023020 | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | | |



| |
|---|
| Test Result BrAC g/210L 0.00 |
|---|

Section 3

What happens next?

1. The Arresting Officer is required to submit the Officer's Impaired Driving Affidavit to the State of Oklahoma.
2. Upon receipt and review of the Officer's Impaired Driving Affidavit and any associated test results, if applicable, the State of Oklahoma may revoke your Oklahoma driving privileges. If you hold a commercial driver license or were driving a commercial motor vehicle, the State of Oklahoma may disqualify your privileges to operate a commercial motor vehicle. You will continue to have driving privileges, if otherwise eligible, unless or until an Order of Revocation and/or Disqualification is issued by the State of Oklahoma.
3. The State of Oklahoma will send you the Order of Revocation and/or Disqualification to the address you have on file with Service Oklahoma. To confirm or change your address on file, contact Service Oklahoma.
4. The Revocation and/or Disqualification will be effective 55 days from the date of the Order, including 10 days allowed for mailing, 47 O.S. §52-116, 753, 754.
5. The Order of Revocation and/or Disqualification will contain important information about your rights to appeal the revocation and/or disqualification of your driving privileges. To appeal the applicable Order of Revocation and/or Disqualification your petition must be filed with the District Court of the County in which you were arrested within 40 days of the date of the order, including 10 days allowed for mailing, 47 O.S. §52-116, 6-211.

In accordance with Title 12 O.S. Section 426, "I state under penalty of perjury under the laws of Oklahoma that the foregoing is true and correct."

Section 4

Date: _____ Place (location when signed): _____ Signature of arresting officer: _____

| | | | |
|---------------------------|--------------------------|-------------------|---------------------|
| NAME PAPE, CHRISTOPHER | AGENCY BOARD OF TESTS | ZIP CODE 00000 | PHONE 4054252460 |
| BADGE 000000484 | | INTER-AGENCY | |



**STATE OF OKLAHOMA
BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE**

Dr. Jarrad Wagner, Ph.D., F-ABFT
Chairman
Post Office Box 36307
Oklahoma City, Oklahoma 73136-2307
Phone: (405) 425-2460 Fax: (405) 425-2490
www.bot.ok.gov

Joshua Smith
Director


Certificate of Analysis – Reference Solution

ALPHA INSTRUMENT
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-002591
09/23/2022 09:37

| Test | g/210L | Time |
|-----------------|--------|-------|
| Air Blank | 0.000 | 09:37 |
| Cal Check | 0.019 | 09:38 |
| Air Blank | 0.000 | 09:39 |
| Cal Check | 0.019 | 09:39 |
| Air Blank | 0.000 | 09:40 |
| Cal Check | 0.019 | 09:41 |
| Air Blank | 0.000 | 09:41 |
| Cal Check | 0.019 | 09:42 |
| Air Blank | 0.000 | 09:42 |
| Cal Check | 0.019 | 09:43 |
| Air Blank | 0.000 | 09:44 |
| Cal Check | 0.019 | 09:44 |
| Air Blank | 0.000 | 09:45 |
| Cal Check | 0.019 | 09:46 |
| Air Blank | 0.000 | 09:46 |
| Cal Check | 0.019 | 09:47 |
| Air Blank | 0.000 | 09:47 |
| Cal Check | 0.019 | 09:48 |
| Air Blank | 0.000 | 09:49 |
| Cal Check | 0.019 | 09:49 |
| Air Blank | 0.000 | 09:50 |
| Cal Check Stats | | |
| Average | 0.0190 | |
| Std Dev | 0.0000 | |
| Rel Std Dev (%) | 0.0000 | |

NIST Traceable Breath-Alcohol Equivalent (g/210L +/- 3%) Ratio of ethanol and water free of test interfering substances.

| | | |
|---|------------------------------|----------------------------|
| Manufacturer: Guth | Manufactured Date: 8/24/2021 | Concentration: 0.02 |
| Lot No.: 21340 | Expiration Date: 8/24/2023 | Equipment Model: Guth 2100 |
| <p>The above solution has been analyzed via Infrared Spectroscopy utilizing the Intoxilyzer 8000 and was performed in accordance with the procedures prescribed by the State Director of Tests pursuant to 40 OAC 30-1-3. The above NIST traceable reference solution meets the requirements for use pursuant to 40 OAC 25-1-2.1. Approved breath-alcohol measurement equipment and reference solutions. Equipment used appears on the current or supplemented Conforming Products List of Calibrating Units for Breath Alcohol Testers (72 FR 34747), published by the National Highway Traffic Safety Administration.</p> | | |
| Signature | Date: 9/23/2022 | |
| Name Christopher Pape | Permit No.: 23020 | |

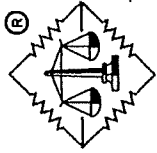

GUTH LABORATORIES, INC.

590 NORTH 67th STREET • HARRISBURG, PA 17111-4611 • TELEPHONE: 717-564-5470

RCVD
5/9/2022
CP

CERTIFICATE OF ANALYSIS

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number **21340** ✓ of Alcohol Reference Solution for Simulator were analyzed by gas chromatography on **August 25, 2021**, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain **0.0242%** (w/vol) ethyl alcohol. The expiration date for this lot number is **August 24, 2023** at 11:59 PM. ✓

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of **0.020 g/210L +/- 3%**.

The alcohol and water used in this solution were free of test interfering substances.

Ted L. Pauley, President
GUTH LABORATORIES, INC.

NIST Traceability:
Testing was conducted using Cerilliant Reference Standard lot number **FN05122004** whose values are traceable to NIST.
All balances are calibrated annually by an outside agency using NIST traceable weights.
Calibration verification is done prior to each use utilizing NIST traceable weights.



14/20

STATE OF OKLAHOMA BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE

Post Office Box 36307

Dr. Jarrad Wagner, Ph.D., F-ABFT
Chairman

Oklahoma City, Oklahoma 73136-2307
Phone: (405) 425-2460 Fax: (405) 425-2490
www.bot.ok.gov


Joshua Smith
Director

Certificate of Analysis – Reference Solution

ALPHA INSTRUMENT 2
Intoxilyzer – Alcohol Analyzer SN 80-004923
Model 8000
03/08/2022 11:53

| Test | g/210L | Time |
|-----------------|--------|-------|
| Air Blank | 0.000 | 11:53 |
| Cal Check | 0.041 | 11:54 |
| Air Blank | 0.000 | 11:54 |
| Cal Check | 0.041 | 11:55 |
| Air Blank | 0.000 | 11:56 |
| Cal Check | 0.042 | 11:56 |
| Air Blank | 0.000 | 11:57 |
| Cal Check | 0.042 | 11:58 |
| Air Blank | 0.000 | 11:58 |
| Cal Check | 0.042 | 11:59 |
| Air Blank | 0.000 | 11:59 |
| Cal Check | 0.041 | 12:00 |
| Air Blank | 0.000 | 12:01 |
| Cal Check | 0.041 | 12:01 |
| Air Blank | 0.000 | 12:02 |
| Cal Check | 0.041 | 12:02 |
| Air Blank | 0.000 | 12:03 |
| Cal Check | 0.041 | 12:04 |
| Air Blank | 0.000 | 12:04 |
| Cal Check | 0.042 | 12:05 |
| Air Blank | 0.000 | 12:05 |
| Cal Check Stats | | |
| Average | 0.0414 | |
| Std Dev | 0.0005 | |
| Rel Std Dev(%) | 1.2473 | |

| | |
|--|-------------------------------|
| NIST Traceable Breath-Alcohol Equivalent (g/210L +/- 3%) Ratio of ethanol and water free of test interfering substances. | |
| Manufacturer: Guth | Manufactured Date: 03/01/2021 |
| Lot No.: 21070 | Expiration Date: 03/01/2023 |
| The above solution has been analyzed via Infrared Spectroscopy utilizing the Intoxilyzer 8000 and was performed in accordance with the procedures prescribed by the State Director of Tests pursuant to 40 OAC 30-1-3. The above NIST traceable reference solution meets the requirements for use pursuant to 40 OAC 25-1-2.1. Approved breath-alcohol measurement equipment and reference solutions. Equipment used appears on the current or supplemented Conforming Products List of Calibrating Units for Breath Alcohol Testers (72 FR 34747), published by the National Highway Traffic Safety Administration. | |
| Signature  | Date: 03/08/2022 |
| Name Christopher Pape | Permit No.: 23020 |

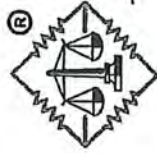**GUTH LABORATORIES, INC.**

590 NORTH 67th STREET • HARRISBURG, PA 17111-4511 • TELEPHONE: 717-564-5470

RCVD
8/26/21
CP**CERTIFICATE OF ANALYSIS**

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number **21070** of Alcohol Reference Solution for Simulator were analyzed by gas chromatography on **March 2, 2021**, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain **0.0488%** (w/vol) ethyl alcohol. The expiration date for this lot number is **March 1, 2023** at 11:59 PM.

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of **0.040 g/210L +/- 3%**.

The alcohol and water used in this solution were free of test interfering substances.

Ted L. Pauley, President
GUTH LABORATORIES, INC.

NIST Traceability:

Testing was conducted using Cerilliant Reference Standard lot number **FN08211802** whose values are traceable to NIST.
All balances are calibrated annually by an outside agency using NIST traceable weights.
Calibration verification is done prior to each use utilizing NIST traceable weights.



**STATE OF OKLAHOMA
BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE**

Dr. Jarrad Wagner, Ph.D., F-ABFT
Chairman

Post Office Box 36307
Oklahoma City, Oklahoma 73136-2307
Phone: (405) 425-2460 Fax: (405) 425-2490
www.bot.ok.gov

Joshua Smith
Director


Certificate of Analysis – Reference Solution

ALPHA INSTRUMENT
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-002591
06/27/2022 11:40

| Test | g/210L | Time |
|-----------------|--------|-------|
| Air Blank | 0.000 | 11:40 |
| Cal Check | 0.100 | 11:41 |
| Air Blank | 0.000 | 11:41 |
| Cal Check | 0.100 | 11:42 |
| Air Blank | 0.000 | 11:43 |
| Cal Check | 0.100 | 11:43 |
| Air Blank | 0.000 | 11:44 |
| Cal Check | 0.100 | 11:45 |
| Air Blank | 0.000 | 11:45 |
| Cal Check | 0.100 | 11:46 |
| Air Blank | 0.000 | 11:47 |
| Cal Check | 0.101 | 11:47 |
| Air Blank | 0.000 | 11:48 |
| Cal Check | 0.101 | 11:49 |
| Air Blank | 0.000 | 11:49 |
| Cal Check | 0.101 | 11:50 |
| Air Blank | 0.000 | 11:50 |
| Cal Check | 0.100 | 11:51 |
| Air Blank | 0.000 | 11:52 |
| Cal Check | 0.100 | 11:52 |
| Air Blank | 0.000 | 11:53 |
| Cal Check Stats | | |
| Average | 0.1003 | |
| Std Dev | 0.0004 | |
| Rel Std Dev(%) | 0.4815 | |

NIST Traceable Breath-Alcohol Equivalent (g/210L +/- 3%) Ratio of ethanol and water free of test interfering substances.

| | | |
|---|-------------------------------|----------------------------|
| Manufacturer: Guth | Manufactured Date: 03/07/2022 | Concentration: 0.10 |
| Lot No.: 22080 | Expiration Date: 03/07/2024 | Equipment Model: Guth 2100 |
| <p>The above solution has been analyzed via Infrared Spectroscopy utilizing the Intoxilyzer 8000 and was performed in accordance with the procedures prescribed by the State Director of Tests pursuant to 40 OAC 30-1-3. The above NIST traceable reference solution meets the requirements for use pursuant to 40 OAC 25-1-2.1. Approved breath-alcohol measurement equipment and reference solutions. Equipment used appears on the current or supplemented Conforming Products List of Calibrating Units for Breath Alcohol Testers (72 FR 34747), published by the National Highway Traffic Safety Administration.</p> | | |
| Signature | Date: 06/27/2022 | |
| Name Christopher Pape | | Permit No.: 23020 |

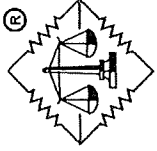

GUTH LABORATORIES, INC.

590 NORTH 67th STREET • HARRISBURG, PA 17111-4511 • TELEPHONE: 717-564-5470

CERTIFICATE OF ANALYSIS

RCVD
5/9/2022
CP

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number **22080** of

Alcohol Reference Solution for Simulator were analyzed by gas chromatography on **March 9, 2022**, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain **0.1214%** (w/vol) ethyl alcohol. The expiration date for this lot number is **March 7, 2024** at 11:59 PM.

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of **0.100 g/210L +/- 3%**.

The alcohol and water used in this solution were free of test interfering substances.

Ted L. Pauley, President
GUTH LABORATORIES, INC.

NIST Traceability:

Testing was conducted using Cerilliant Reference Standard lot number **FN03052002** whose values are traceable to NIST.

All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.



**STATE OF OKLAHOMA
BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE**

Dr. Jarrad Wagner, Ph.D., F-ABFT
Chairman

Post Office Box 36307
Oklahoma City, Oklahoma 73136-2307
Phone: (405) 425-2460 Fax: (405) 425-2490
www.bot.ok.gov


Joshua Smith
Director

Certificate of Analysis – Reference Solution

ALPHA INSTRUMENT
Intoxilyzer - Alcohol Analyzer SN 80-002591
Model 8000
04/19/2023 10:57

| Test | g/210L | Time |
|-----------------|--------|-------|
| Air Blank | 0.000 | 10:58 |
| Cal Check | 0.198 | 10:58 |
| Air Blank | 0.000 | 10:59 |
| Cal Check | 0.198 | 11:00 |
| Air Blank | 0.000 | 11:00 |
| Cal Check | 0.199 | 11:01 |
| Air Blank | 0.000 | 11:02 |
| Cal Check | 0.199 | 11:02 |
| Air Blank | 0.000 | 11:03 |
| Cal Check | 0.198 | 11:04 |
| Air Blank | 0.000 | 11:04 |
| Cal Check | 0.199 | 11:05 |
| Air Blank | 0.000 | 11:06 |
| Cal Check | 0.199 | 11:06 |
| Air Blank | 0.000 | 11:07 |
| Cal Check | 0.199 | 11:08 |
| Air Blank | 0.000 | 11:08 |
| Cal Check | 0.199 | 11:09 |
| Air Blank | 0.000 | 11:10 |
| Air Blank | 0.000 | 11:11 |
| Cal Check Stats | | |
| Average | 0.1987 | |
| Std Dev | 0.0004 | |
| Rel Std Dev(%) | 0.2431 | |

| | |
|---|-----------------------------|
| NIST Traceable Breath-Alcohol Equivalent (g/210L +/- 3%) Ratio of ethanol and water free of test interfering substances. | |
| Manufacturer: Guth | Manufactured Date: 2/7/2022 |
| Lot No.: 22050 | Expiration Date: 2/7/2024 |
| <p>The above solution has been analyzed via Infrared Spectroscopy utilizing the Intoxilyzer 8000 and was performed in accordance with the procedures prescribed by the State Director of Tests pursuant to 40 OAC 30-1-3. The above NIST traceable reference solution meets the requirements for use pursuant to 40 OAC 25-1-2.1. Approved breath-alcohol measurement equipment and reference solutions. Equipment used appears on the current or supplemented Conforming Products List of Calibrating Units for Breath Alcohol Testers (72 FR 34747), published by the National Highway Traffic Safety Administration.</p> | |
| Signature | Date: 4/19/2023 |
| Name Christopher Pape | Permit No.: 23020 |

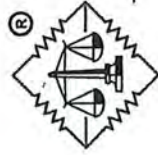

GUTH LABORATORIES, INC.

590 NORTH 67th STREET • HARRISBURG, PA 17111-4511 • TELEPHONE: 717-564-5470

CERTIFICATE OF ANALYSIS

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number **22050** of Alcohol Reference Solution for Simulator were analyzed by gas chromatography on **February 9, 2022**, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain **0.2437%** (w/vol) ethyl alcohol. The expiration date for this lot number is **February 7, 2024** at 11:59 PM.

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of **0.200 g/210L +/- 3%**.

The alcohol and water used in this solution were free of test interfering substances.

RCVD
2/23/23
cp

Ted L. Pauley, President
GUTH LABORATORIES, INC.

NIST Traceability:
Testing was conducted using Cerilliant Reference Standard lot number **FN06231703** whose values are traceable to NIST.
All balances are calibrated annually by an outside agency using NIST traceable weights.
Calibration verification is done prior to each use utilizing NIST traceable weights.

Location Returned From

Ellis Co

Unit Serial #

80-002696

Date: 20/20

5/1/23

Inventory Process

- Unit checked in to BOT
- New Unit Checked out

- Unlink Gas Canister from Unit
- Remove old cylinder from unit (After Lab Evaluation)

Old Cylinder #: 11821080A2-008

- Old cylinder disposed in inventory
- Old cylinder archived

Check in Procedure

- Connect unit to ethernet connection 172.019.004. 141
- Download any stored records: 2
- Change unit location to Serial number
- Update Unit Date and time

Previous time: 5/10/23 1115

New Time: 5/10/23 1116

- Execute a new Forms load on the unit

Notes from the Field

- Routine Swap
- Issues in the Field

Lab Evaluation

Breath Hose Hold Vacuum?:

Test Chamber Hold Vacuum?:

DVM

3µm 12975 Chamber Temp 47.00 Snap Test ok?

9µm 13802 Breath Hose Temp 45 RFI Detected?

Peripheral Tests

Speaker (F5) Power LED Green (F6) Red (F7) Fan (F9)

Gas Regulator (F11) Pump (F10) Receipt paper

Atmospheric Sensor
 Pressure 975 Correction Factor 1.03
 Flow Sensor Resistance 100 Tank Sensor Value 983

Is regulator close to Tank Sensor Value?
 Yes / No

Sent to CMI (If needed)
 Date Sent _____
 Invoice Sent _____
 Date Returned _____

Reason Sent

Repairs needed:

Repairs Made:

Ready for Bench Check?

Inventory Process

- Linked Cylinder to Unit
- Changed Unit Status
- Bench Check Report Completed

Bench Check Technician Report

- Mock Subject Test Completed
- Certificate of Calibration and Operation
- Service Overview Completed

To complete if unit sent to CMI

- CMI Workorder
- CMI Certificate of Calibration
- CMI Service Evaluation Form