

## STATE OF OKLAHOMA BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE

Dr. Jarrad Wagner, Ph.D. F-ABFT Chairman

Post Office Box 36307 Oklahoma City, Oklahoma 73136-2307 Phone: (405) 425-2460 Fax: (405) 425-2490 www.bot.ok.gov

Joshua Smith Director

#### ATTESTATION

I, the undersigned *Joshua Smith* (Director/Records Custodian), attest under oath pursuant to 12 O.S. §2902 (11)(a) this record was made, at or near the time of the occurrence of the matters set forth by or from information transmitted by a person having knowledge of those matters; was kept in the course of the regularly conducted business activity of the Board; and was made pursuant to the regularly conducted activity. The duplicate record provided is listed below:

1) The attached "Service Overview" is a true and accurate copy of maintenance performed on the Intoxilyzer 8000, serial number 80-003290 , in compliance with the agency's Maintenance Policy and Procedures.

This packet consists of 13 pages.

Signed //

Name of Position: Director/Records Custodian

Date of Attestation: 01/11/2021

## **Service Overview**



3600 N. Martin Luther King Bldg #9 Oklahoma City, OK 73111 405-425-2460

BENCH CHECK DATE: 01-11-2021

SERIAL: 80-003290

#### DOCUMENTATION WITHIN SERVICE PACKET:

- **BOT Certificate of Calibration and Operation**
- Bench Check
- Bench Check Technician's Report
- **BOT Dry Gas Certificate of Analysis**
- Mock Subject Test
- BOT Certificates of Analysis on Certified Solutions (4)



## STATE OF OKLAHOMA BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE

Administrative Office at:
Board of Tests for Alcohol and Drug Influence
Post Office Box 36307
Oklahoma City, Oklahoma 73136-2307

### **Certificate of Calibration and Operation**

This is to certify that the calibration of INTOXILYZER 8000, serial number 80-003290, manufactured by CMI Inc., a subsidiary of MPD, Inc. of Owensboro, Kentucky, was tested and found to be in compliance with National Highway Traffic Safety Administration Standard for Devices to Measure Breath Alcohol (F.R., No. 179 48705-48710 Sept. 17, 1993) and with International Recommendation OIML R 126 Edition 1998 (E). Calibration solutions are traceable to NIST standard material 1828.

In addition the operational procedures and functions of this Intoxilyzer have also been tested and found to be in compliance with the Oklahoma Board of Tests for Alcohol and Drug Influence - Rules and Regulations OAC Title 40: 25-1-2 and 30-1-3 and is hereby certified for use as an evidential breath testing device within this state.

CERTIFICATION TYPE	DATE
Manufacturer Calibration	07-16-08
Board of Tests Bench Check	08-22-08
Board of Tests Bench Check	09-02-09
Board of Tests Bench Check	07-12-11
Manufacturer Calibration	10-13-11
Board of Tests Bench Check	10-24-11
Board of Tests Bench Check	10-10-13
Board of Tests Bench Check	10-05-15
Board of Tests Bench Check	01-20-17
Board of Tests Bench Check	10-02-18
Manufacturer Calibration	08-26-20
Board of Tests Bench Check	08-31-20
Board of Tests Bench Check	01-11-21



Presentation of this form certifies that the Intoxilyzer, SN: 80-003290, functioned properly at the time of the breath test and hereby certifies the breath test result as valid.

Date:

Re: Request for Log of Test on Intoxilyzer 8000

Name:

#### Dear Sir:

The Log of Test information from a breath test conducted on a Intoxilyzer 8000 is printed in Section 2 of the Affidavit and is no longer retained as a separate form as with the Intoxilyzer 5000. The 8000 uses a nitrogen/ethanol dry gas for the external control test. The dry gas cylinder has the capacity to perform approximately 325 tests and does not require maintenance until its expiration date or can no longer provide the required pressure to perform a control test. The cylinder has an expiration date of two years from date filled at the manufacturer and is programmed in the 8000 upon installation, both dates are printed in Section 2 of the Affidavit. The Intoxilyzer is programmed to disable itself when the date of expiration is reached or the pressure monitor indicates 25 lbs psi.

Due to the reliability and stability of dry gas, maintenance on the 8000 is not required every 30 days or 25 tests, therefore maintenance is no longer performed by officers in the field. The Board of Tests for Alcohol and Drug Influence certifies each Intoxilyzer 8000 for field use by performing a bench check. The bench check is designed to check the operation and procedures programmed in the 8000 to ensure it is functioning properly for use in Oklahoma. A bench check is performed every time a new dry gas cylinder is installed or the instrument has returned from the manufacturer for repair. Once the instrument passes the bench check, all that is required to set it up for field use is to plug it in and turn it on and does not require a permit to set up. The Board of Tests has the ability to monitor each instrument once it is set up in the field.

Each Intoxilyzer has a database of valid permit numbers that it checks prior to allowing a breath test to be performed. The database is updated on a regular basis, if a permit number is invalid that officer will not be able to conduct a breath test on the 8000. If the 8000 prints an affidavit the breath test operator listed in Section 2 had a valid permit at the time of the test and the breath test is valid.

If you have any further questions please feel free to call me.

Sincerely,

J. Robert Blakeburn (OBA11059) State Director of Tests for Alcohol and Drug Influence

poher Blokelung

## STATE OF OKLAHOMA - BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE

1/11/2021	LO:53:4		NSTRUMENT LOCATION 80-00329	0	
RUMENT TYPE	LU:53:4	EL	80-00329	INSTRUMENT SE	RIAL NUMBER
ntoxilyzer	I-8000			80-00	
		TEST RE	ESULTS		
Diagnostics					
PASS					
Mouth Alcohol Test					
PASS					
RFI Test					
PASS		Du	plicate Breath	Sampling Re	eulte
Abort Test		Du	phoate breath	Sampling Re	Tis in The
PASS		Value	Test 1	Test 2	Control Test
Improper Sample Tes	st _				1501
PASS		.02	0.019	0.019	0.020
Interferent Detect Tes	st	.04	0.038	0.037	0.039
PASS		.10	0.098	0.098	0.099
No Sample Given Tes	st		1 10 10 10 10 10 10 10 10 10 10 10 10 10	172-172-172-17	1 1 1 1 1 1 1 1 1 1
PASS		.20	0.197	0.194	0.198
Insufficient Sample Te	est		PA	SS	
PASS					
Air Blank Contamination	Test				
PASS					
Range Exceeded Tes	t				
PASS					

PAPE,


3600 N. Martin Luther King Bldg #9 Oklahoma City, OK 73111 405-425-2460

## **BENCH CHECK TECHNICIAN'S REPORT**

INSTRUMENT ntoxilyzer Model 8000 GAS CANISTER OT# 20520080A1-026 REFERENCE	Simulator Mo	80-003  Date  09/05/2  del  Guth 2100	
8000 <b>GAS CANISTER</b> OT # 20520080A1-026	EXP Simulator Mo	80-003  Date  09/05/2  del  Guth 2100	022
GAS CANISTER OT# 20520080A1-026	_Simulator Mo	Date 09/05/2 del Guth 2100	022
OT# 20520080A1-026	_Simulator Mo	09/05/2 odel Guth 2100	
	_Simulator Mo	09/05/2 odel Guth 2100	
REFERENCE	100	Guth 2100	C.:#- 0400
	100	Guth 2100	0.45 0400
Guth 2100 Guth 2		Annual Control of the	(a)(iii) (711)()
Suit 2100 Suit 2		ŧ	Guii 2 100
DR2834 DR370	00	DR2455	DR3470
	Concentra		
.02 .04		.10	.20
00070	LOT#		20100
20070 2006	Carrier and the second	20190	20160
00/40/0000	Manufacture		1
02/13/2020 02/10/2		04/06/2020	03/18/2020
00/40/0000	Expiration I		
02/13/2022 02/10/2	75343	04/06/2022	03/18/2022
	lution Commis	sion Date	
1/4/2021 1/4/20	21	1/4/2021	1/4/2021
REASON FOR BENCH CHEC POST REPAIR CYLINDER TROUBLESHOOTING RO COMMENTS	R REPLACE	EMENT INTENANCE	
Unit came in with multiple RFI dete	ects in 1 day.	After evaluation, ur	nit was not the issue
111			
With offer 230	20		
echnician Signature / Permit #			

Christopher Pape #23020

Printed Name

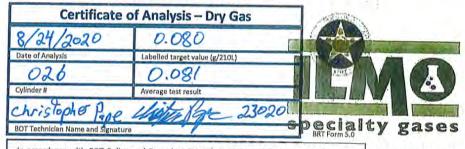

in accordance with BOT Policy and Procedure Statement இது இரும் இரும் இரும் இரும் Jacksonville, IL 62651-0790 referenced dry gas cylinder is suitable for use as an external reference with the intexilvzer 8000.
217-245-2183 • Fax: 217-243-7634 • www.ilmoproducts.com

## Certificate of Analysis

Certificate ID:

13202

Part #:

BAC105L080T

Cylinder Size:

105L

Lot Number:

Expiration:

20520080A1 9/5/2022

RCV'D 8/19/2020 CP

0.080 BAC (For the calibration of instruments used to determine breath alcohol concentration)

Contents:

105 Liters @ 1000 psig 70°F (21°C)

Analytical

Component:

Reported Concentration:

Accuracy (U, k=2):

Analytical Method:

Ethanol Nitrogen 208 ppm Balance

+/-0.002 BAC(G/210L)

[5.2 ppm]

Distributed by:

CMI Inc.

316 East Ninth Street

Owensboro, KY 42303 -

Phone 866-835-0690 www.alcoholtest.com Operator's Signature

ALPHA INSTRUMENT

Model 8000

08/24/2020

Test

Air Blank

Cal Check

Air Blank

Std Deu

Cal Check Stats Average

Rel Std Dev(%)

Intoxilyzer - Alcohol Analyzer

g/210L

0.000

0.081

0.000

0.081

0.000

0.081

0.000

0.081

0.000

0.081

0.000

0.082

0.000

0.081

0.000

0.081

0.000

0.082

0.000

0.081

0.000

0.0812

0.0004

0.5192

SN 80-002591

09:20

Time

09:21

19:21

09:22

09:22

09:23

09:23

09:23

09:24

09:24

09:25

09:25

09:25

09:26

09:26

09:27

09:27

09:28

09:28

09:29

09:29

09:29

\*Traceable to:

Certified Reference Material - 262.4 µmol/mol

Ethanol in Nitrogen - Serial No. GN0015026 Lot No. 050319E11

Store in dry area, away from sources of heat, ignition and direct sunlight. Do not allow storage area to exceed 52 °C (125 °F).

Specialty/Gas Lab Tech

Issuance Date

Accreditation #61895

he calibration results within this certificate were obtained using equipment and standards capable of producing analytical results traceable to NIST, and apply only to the items ontained on this certificate. ILMO Products Company makes no warranty or representation as to the suitability of the use of any information provided for any particular urpose. The information use is at the sole discretion and risk of the user. Liability shall be limited to established replacement cost of this material or service.

#### 8/13

# OFFICER'S AFFIDAVIT AND NOTICE OF REVOCATION/DISQUALIFICATION

ARREST DATE		ARREST TIME	2.16/17/17	COUNTY#	CIT	Y#	CITATION	,	_	Section
01/11/20: ARREST LOCATION	21	1100		55	CITY	Y2.	1 1 1 1 1 1 1 1 1 1 1 1 1 1 1 1 1 1 1 1	COUNTY	,	
BOT					City				LAHO	MA
NORTH, G	AYLON LEVE	RNE			0.8 /	09/19	61 M	RACE	510	165
ADDRESS	rin Luther		CITY	ALIONA C			STATE	1.77	ZIP COD	E
DRIVER LICENSE #		KING A	OKL	AHOMA C EXPIRATION D	ATE	STATE 0	OK	RESTRIC	731	5 3 ENDORSEMENTS
T0834548	71	Incorn		06/30	A COLUMN TO THE PARTY OF THE PA	OK	A		Total	M
MOCK		MOCE	7		TAC M	OCK		OK.	N CMV	N N
nighway, street, turnpike,		ublic place, or upon a prival	d i nad reasonable gro le road, street, alley or	lane which provides	access to one or	g of in actual pr more single or r	nulti-family dwel	a motor ver ling within th	nis state w	a public road, nile under the
escribe person's condition	n (odor, actions, etc.):	Count / Was								
	OTING BENCH	THE PERSON	WAS INFORMED OF							
REATH-ALCOHOL ANA	LYSIS RECORD, REPORT		IN SECTION 1 ABON	E, and LOG OF TE	ST AND MAINTE	NANCE RECOF	RD.			Section 2
I-8000			-003290	BOARD	OF TES					
Dry			0.080		23020	R				
MANUFACTURER			LOT#	00071	DATE CYLINDE 01/11	RINSTALLED	DATI	CYLINDE	REXPIRE	S
ILMO DBSERVATION START	TEST DATE	MOUTHPIECE	1					9/05	/202	2
1110 RESULT TABLE Test	01/11/202 g/210L	1 G	period v	ite (or greater) vas observed p the time of tes	orior to and	continuing	on	OF .	HE S	TATE
Air Blank Subject Test Breath Volu Air Blank Wait Air Blank Subject Test Breath Volu Air Blank	0.00 0.00 0.00	11:40 11:42 11:42 11:42	complian procession	ered said B ce with the edure of the	applicabl	e operati	7.00	SHEAT BE	1907	
Cal Check	0.08	11:43	PAPE, CH	IRISTOPHEI	p				est Re	
Air Blank	0.00	11:44	BADGE # 0000000	484	PERMIT # 00000	23020		100	AC g/:	
ays or more, if you refu operating a CMV or CDIs is a result of this arrest to ou may be eligible to evocation as a result of equest IDAP within this EERVING OF NOTI accountability Programmers	b) DAYS FROM DATE SEI sed or failed the state's test. holder) OR a test result op yfiling a petition in the Di participate in the Impair of this arrest from appear rty (30) days of the date CE: I personally hand am to the person nan	st(s). In addition, your co of ,04 or more while oper, strict Court of the County ad Driver Accountability ring on your driving red listed below. Mall the r d delivered a copy of ned in Section 1 abo	mmercial driving privating a commercial reliable to the which you were a program (IDAP) actord. Participation request for IDAP to the above contained.	our driver's license illege in this state w notor vehicle. You rrested. Your peliti dministered by the n IDAP may reduc the Department of	vill be disqualifie may appeal any ion must be filed e Department o ce the amount of f Public Safety	d if you refused Departmental I within thirty (3 If Public Safety of fees you with at P.O. Box 11	d or failed the s action against 0) days of the o y (DPS). Com Il be required 415, Oklahom otice of Dep	tate's test(s your driver date listed to pletion of l to pay to to la City, OK	ilicense below. DAP may he State. 73136.	You must
	itle 12 O.S. Section 4			under the laws	of Oklahoma	that the for	egoing is tru	e and co	rrect."	Section 4
	2021 Place (locat			IAGENCY		of arresting o				
PAPE, CH	IRISTOPHER	GENCY		BOARD	OF TE	STS ZIP CODE 00000	PHONE 4054	2524	60	
		TOTALCT				55500	14054	4724	50	
000000048										Continu 6
000000048			TITLE	ADDRESS			PHONE			Section 5
THER WITNESSE			TITLE	ADDRESS			PHONE			Section 5

## OFFICER'S AFFIDAVIT AND NOTICE OF

9/13

M IV PLACARD N IN	COUNTY OKL RACE W	STATE OK ASS A	.9 OK F TAG # MOCK s driving or in actual phys	AHOMA CITY EXPIRATION DATE 06/30/20	₹ nd I had reasonable gr	MODEL MOCK	LON LEVERN	ADDRESS 3600 MARTIN DRIVER LICENSE # T083454871 VEHICLE MAKE MOCK
IGHT WEIGHT 510 165 P CODE 73133 DNS ENDORSEME M IN PLACARD N Upon a public road, tate while under the	RESTRIC TATE OK	STATE OK  ASS A	STATE CL. OR FAGE MOCK Striving or in actual phys	AHOMA CITY EXPIRATION DATE 06/30/203	OKI	MODEL MOCK	I LUTHER K.	BOT SUBJECT NAME NORTH, GAYI NORTH, GAYI NORTH, GAYI NORTH 3600 MARTIN ORIVER LICENSE # T083454871 VEHICLE MAKE MOCK
IGHT WEIGHT 510 165 P CODE 73133 DNS ENDORSEME M IN PLACARD N Upon a public road, tate while under the	RACE W RESTRIC	STATE OK  ASS A	STATE CL. OR FAGE MOCK Striving or in actual phys	AHOMA CITY EXPIRATION DATE 06/30/20	OKI	MODEL MOCK	I LUTHER K.	SUBJECT NAME NORTH, GAYI NORTH, GAYI NORTH, GAYI 3600 MARTIN ORIVER LICENSE # T083454871 VEHICLE MAKE MOCK
510   165 P CODE 73133 INS ENDORSEME M IN PLACARD N IN PLACARD N In upon a public road, tate while under the	RESTRIC	STATE OK  ASS A	STATE CL. OR FAGE MOCK Striving or in actual phys	AHOMA CITY EXPIRATION DATE 06/30/20	OKI	MODEL MOCK	I LUTHER K.	NORTH, GAYI DDRESS 3600 MARTIN RIVER LICENSE # T083454871 EHICLE MAKE MOCK
P CODE 73133  ONS ENDORSEME M PLACARD N PLACARD N Pupon a public road, tale while under the	RESTRIC TATE OK	STATE OK ASS A	STATE CL.  9 OK A  TAG# MOCK  6 driving or in actual phys	AHOMA CITY EXPIRATION DATE 06/30/20	OKI	MODEL MOCK	I LUTHER K.	DDRESS 3600 MARTIN RIVER LICENSE # T083454871 EHICLE MAKE MOCK
DNS ENDORSEME M  IV PLACARD  N upon a public road, tate while under the	TATE OK	ASS A	.9 OK F TAG # MOCK s driving or in actual phys	EXPIRATION DATE 06/30/20	₹ nd I had reasonable gr	MODEL MOCK son was arrested, and	ation, the above named pers	RIVER LICENSE # FO 8 3 4 5 4 8 7 1 EHICLE MAKE MOCK
M PLACARD N PLAC	TATE OK	A sical control	.9 OK F TAG # MOCK s driving or in actual phys	06/30/203	nd I had reasonable gr	MOCK son was arrested, and ace, or upon a private	ation, the above named pers	T083454871 HICLE MAKE MOCK
N PLACARD N Upon a public road, tate while under the	OK a motor vehi	sical control	TAG # MOCK	ounds to believe the person w	nd I had reasonable gr	MOCK son was arrested, and ace, or upon a private	ation, the above named pers	HICLE MAKE MOCK
N N upon a public road, tate while under the	OK a motor vehi	sical control	MOCK s driving or in actual phys	ounds to believe the person w. lane which provides access to	nd I had reasonable gr	MOCK son was arrested, and ace, or upon a private	ation, the above named pers	MOCK
tale while under the	a motor veh ling within thi	sical control ulti-family dw	s driving or in actual phys one or more single or mul	ounds to believe the person w. lane which provides access to	nd I had reasonable gr	son was arrested, and	ation, the above named pers	of the second se
Capit						rohibited by law.	parking lot, or other public pla intoxicating substances as pound umstances):	thway, street, turnpike, public p luence of alcohol and/or other le escribe driving behavior or circu MOCK SUBJECT
Carif							, actions, etc.):	scribe person's condition (odor
Carit						FCV	NG DENCH CUI	TROUBLESHOOTI
Carit						ECK	NG BENCH CHI	.ROUBEBHOOTI
Cast								
Cant				THE IMPLIED CONSENT TE				
Section		),	AINTENANCE RECORD.				RECORD, REPORT OF TH	EATH-ALCOHOL ANALYSIS
			TESTS	BOARD OF	NUMBER -003290			UIPMENT [ - 8000
					TARGET VALUE	180		ANDARD
					0.080			Dry
	E CYLINDER		LINDER INSTALLED	DATE CY	LOT#			NUFACTURER
2022	9/05/	(	11/2021	0080A1 01/	20520	Notes that work has		LMO
		1	ation/observation	ite (or greater) depriv	A 15-min	MOUTHPIECE G	/11/2021	SERVATION START TEST
E Ch			and continuing	vas observed prior to	period	1 G	./11/2021	SULT TABLE
ATATE	193		nistration.	the time of test admi	to	44.5	-/2105	
** * (%)	194	_	032-1-1011		F 15 15 15	Time	g/210L	Test
· * * /0	3/ * *	in	Alcohol Test i	ered said Breath	Ladminis	11:39	0.00	Air Blank
1 1 1 1 1 1 1 1 1 1 1 1 1 1 1 1 1 1 1	2 × 00			ce with the appli		11:40	0.00	Subject Test
图51.8(+ * * )?	5 ***	В			The second secon		1.757 LITERS	Breath Volume
** 0	EA ++		d of Tests.	edure of the Boar	proc	11:40	0.00	
	19/1	0	22020	to New	1/1/2:		0.00	
	1		23020	19.47	0710			
1907				F OPERATOR	SIGNATURE C	22.72		Breath Volume
					NAME OF THE OWNER OWNER OF THE OWNER OWNE	11:43	0.00	
t Result	T			ID TOMODUED	TO STATE OF THE ST	11:43	0.08	
Cg/210L	Br		to the same of	49 440 4 41 4 4 4 4 4 4 4 4 4 4 4 4 4 4 4		11:44	0.00	Air Blank
		- 11						
. () (1		_		111111111111111111	AGENCY			
	U			F TESTS	BOARD O			a barrier State and
	Br		d of Tests.	FOPERATOR  HRISTOPHER 0484  F TESTS	SIGNATURE CONTROL OF SIGNATURE	11:42 11:42 11:42 11:43 11:43 11:44	0.00 0.00 0.00 1,683 LITERS 0.00 0.08	Air Blank Wait Air Blank Subject Test Breath Volume Air Blank Cal Check Air Blank

Certifica	ate of Analysis - Wet Bath	- Contraction	
8/11/200	20 0.020		
Date of Analysis	Labellad tan (g/210L)		
20070	0.020	Some of the second	
Lot Number		ABORATORIES,	INC
christophes	I Che The NORTH Can 3 P	RESTO - HARRISBURG, PA 17111-4511	
BOT Technician Name and S	ignature	BRT Form 5.0	

In accordance with BRT - 2.2.0, Maintenance, the above referenced simulator solution is suitable for use as an external reference in maintenance of the intoxilyzer 8000.

ALPHA INSTRUMENT 2

Intoxilyzer - Alcohol Analyzer

Model 8000 08/11/2020

SN 80-004923 09:42 OF ANALYSIS

g/210L	Time
0,000	09: 42
0.020	09:43
0,000	09:43
0.020	09:44
0.000	09:45
0.020	09:45
0.000	09:46
0.020	09:47
0.000	09:47
0.021	09:48
0.000	09:48
0.020 .	09:49
0.000	09:49
0.020	09:50
0.000	09:51
0.020	09:51
0.000	09:52
0.020	09:53
0.000	09:53
0.019	09:54
0.000	09:54
0.0200	
0.0004	
2,3569	
	0.000 0.020 0.000 0.020 0.000 0.020 0.000 0.021 0.000 0.020 0.000 0.020 0.000 0.020 0.000 0.020 0.000

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number 20070 of Alcohol Reference Solution for Simulator were analyzed by chromatography on February 17, 2020, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain 0.0242% (w/vol) ethyl alcohol. The expiration date for this lot number is February 13, 2022 at 11:59 PM.

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of 0.020 g/210L +/- 3%.

The alcohol and water used in this solution were free of test interfering substances.

perator's Signature

Ted L. Pauley, President GUTH LABORATORIES, INC.

NIST Traceability:

Testing was conducted using Cerilliant Reference Standard lot number FN08031603 whose values are traceable to NIST.

All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.



#### CERTIFICATE OF ANALYSIS

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number 20060 of

Alcohol Reference Solution for Simulator were analyzed by
gas chromatography on February 11, 2020, using a Perkin Elmer Gas
Chromatograph Autosystem XL S/N: 610N9030209, and found to contain
0.0483% (w/vol) ethyl alcohol. The expiration date for this lot
number is February 10, 2022 at 11:59 PM.

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of 0.040 g/210L +/- 3%.

The alcohol and water used in this solution were free of test interfering substances.

Operator S Signature

ALPHA UNIT SUB

Model 8000

07/22/2020

Test

Air Blank

Cal Check

Air Blank

Std Deu

Cal Check Stats Average

Rel Std Dev(%)

Intoxilyzer - Alcohol Analyzer

g/210L

0.000

0.040

0.000

0.039

0.000

0.039

0.000

0.039

0.000

0.039

0.000

0.039

0.000

0.040

0.000

0.040

0.000

0.040

0.000

0.040

0.000

0.0395

0.0005

1.3342

SN 80-003400

14:57

Time

14:57

14:58

14:59

14:59

15:00

15:01

15:01

15:02

15:02

15:03

15:04

15:04

15:05

15:06

15:06

15:07

15:07

15:08

15:09

15:09

15:10

Ted L. Pauley, President GUTH LABORATORIES, INC.

NIST Traceability:

Testing was conducted using Cerilliant Reference Standard lot number FN12181501 whose values are traceable to NIST.

All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.

06 16 28   100 Date of Analysis   Labelled target value (g/2101) 20 90   .098	CEITHICATE (
Date of Analysis Labelled target value (g/210L)  20196 .098	06/16/20
100	
The second secon	20190
Lot Number Auglage text result (g/210L)	Lot Number

LIBONATORIES, INC. NGET BR CHOME OF THE URG, PA 17111- 4511 . TELEPHONE: 717-564-5470

In accordance with BRT - 2.2.0, Maintenance, the above referenced simulator solution is suitable for use as an external reference in maintenance of the Intoxilyzer 8000.

16:18

#### CERTIFICATE OF ANALYSIS

ALPHA INSTRUMENT Intoxilyzer - Alcohol Analyzer SN 80-002591 Certified Alcohol Reference Solution for Simulator Model 8000

Test 9/210L Time Air Blank 0.000 16:18 Cal Check 0.099 16:19 Air Blank 0.000 16:20 Cal Check 0.098 16:20 Air Blank 0.000 16:21 Cal Check 0.097 16:22 Air Blank 0.000 16:22 Cal Check 0.098 16:23 Air Blank 0.000 16:24 Cal Check 0.098 16:24 Air Blank 0.000 16:25 Cal Check 0.098 16:26 Air Blank 0.000 16:26 Cal Check 0.098 16:27 Air Blank 0.000 16:27 Cal Check 0.097 16:28 Rir Blank 0.000 16:29 tal Check 0.098 16:29 lir Blank 0.000 16:30 'al Check 0.098 15:31 Hir Blank 0.000 16:31 'al Check Stats Average 0.0979

06/16/2020

Std Deu

Rel Std Deu[%]

Random Samples of Lot Number 20190 Alcohol Reference Solution for Simulator were analyzed by gas chromatography on April 8, 2020, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain 0.1206% (w/vol) ethyl alcohol. The expiration date for this lot number is April 6, 2022 at 11:59 PM.

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of 0.100 g/210L +/- 3%.

The alcohol and water used in this solution were free of test interfering substances.

Operato Signature

0.0005


0.5798

Ted L. Pauley, President GUTH LABORATORIES, INC.

NIST Traceability:

Testing was conducted using Cerilliant Reference Standard lot number FN02271802 whose values are traceable to NIST.

All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.

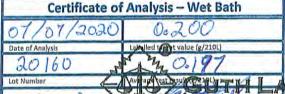

SN 80-003400

15:25

Time

15:26

15:26

15:27

15:28

15:28

15:29

15:30

15:31

15:31

15:32

15:33

15:33

15:34

15:35

15:35

15:36

15:37

15:37

15:38

15:38

SUTH LABORATORIES, INC.

CLASUSU BRT Form 6.0

In accordance with BRT – 2.2.0, *Maintenance*, the above referenced simulator solution is suitable for use as an external reference in maintenance of the Intoxilyzer 8000.

#### CERTIFICATE OF ANALYSIS

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number 20160 of
Alcohol Reference Solution for Simulator were analyzed by
gas chromatography on March 20, 2020, using a Perkin Elmer Gas
Chromatograph Autosystem XL S/N: 610N9030209, and found to contain
0.2427% (w/vol) ethyl alcohol/ The expiration date for this lot
number is March 18, 2022 at 11:59 PM.

When used in a calibrated Simulator, operating at 34°C +/- 2°C, this solution will give a breath alcohol analysis instrument reading of 0.200 g/210L +/- 3%.

The alcohol and water used in this solution were free of test interfering substances.

Operator is Segnature

LPHR UNIT SUB atoxilyzer - Alcohol Analyzer

g/210L

0,196

0.000

0,196

0.000

0,196

0.000

0.000

0.198

0.000

0.197

0.000

0.198

0.000

0.198

0.000

0.197

0.000

0.1971

0.0008

0.4442

0008 lebc

7/07/2020

ir Blank

al Check

ir Blank

al Check

ir Blank

al Check

ir Blank al Check

ir Blank

31 Check

ir Blank

al Check

ir Blank

luerage .

itd Dev le) Std Dev(2)

1) Check Stats

Ted L. Pauley, President GUTH LABORATORIES, INC.

NIST Traceability:

Testing was conducted using Cerilliant Reference Standard lot number FN08101505 whose values are traceable to NIST.

All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.