

STATE OF OKLAHOMA BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE

Dr. Jarrad Wagner, Ph.D. F-ABFT Chairman

Post Office Box 36307 Oklahoma City, Oklahoma 73136-2307 Phone: (405) 425-2460 Fax: (405) 425-2490 www.bot.ok.gov

Joshua Smith Director

ATTESTATION

I, the undersigned <u>Joshua Smith</u> (Director/Records Custodian), attest under oath pursuant to 12 O.S. §2902 (11)(a) this record was made, at or near the time of the occurrence of the matters set forth by or from information transmitted by a person having knowledge of those matters; was kept in the course of the regularly conducted business activity of the Board; and was made pursuant to the regularly conducted activity. The duplicate record provided is listed below:

1) The attached "Service Overview" is a true and accurate copy of maintenance performed on the Intoxilyzer 8000, serial number 80-003292, in compliance with the agency's Maintenance Policy and Procedures.

This packet consists of 16 pages.

Name of Position: Director/Records Custodian

Date of Attestation: 05/18/2022

Service Overview



3600 N. Martin Luther King Bldg #9 Oklahoma City, OK 73111 405-425-2460

BENCH CHECK DATE: 05-13-2019

SERIAL: 80-003292

DOCUMENTATION WITHIN SERVICE PACKET:

- BOT Certificate of Calibration and Operation
- CMI Certificate of Calibration 04-30-2018
- Bench Check
- Bench Check Technician's Report
- BOT Dry Gas Certificate of Analysis
- CMI Service Evaluation Form
- CMI Work Order
- Mock Subject Test
- BOT Certificates of Analysis on Certified Solutions (4)



STATE OF OKLAHOMA BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE

Administrative Office at:
Board of Tests for Alcohol and Drug Influence
Post Office Box 36307
Oklahoma City, Oklahoma 73136-2307

Certificate of Calibration and Operation

This is to certify that the calibration of INTOXILYZER 8000, serial number 80-003292, manufactured by CMI Inc., a subsidiary of MPD, Inc. of Owensboro, Kentucky, was tested and found to be in compliance with National Highway Traffic Safety Administration Standard for Devices to Measure Breath Alcohol (F.R., No. 179 48705-48710 Sept. 17, 1993) and with International Recommendation OIML R 126 Edition 1998 (E). Calibration solutions are traceable to NIST standard material 1828.

In addition the operational procedures and functions of this Intoxilyzer have also been tested and found to be in compliance with the Oklahoma Board of Tests for Alcohol and Drug Influence - Rules and Regulations OAC Title 40: 25-1-2 and 30-1-3 and is hereby certified for use as an evidential breath testing device within this state.

CERTIFICATION TYPE	DATE
Manufacturer Calibration	07-16-08
Board of Tests Bench Check	08-22-08
Board of Tests Bench Check	09-04-09
Manufacturer Calibration	08-12-11
Board of Tests Bench Check	08-24-11
Board of Tests Bench Check	08-08-13
Board of Tests Bench Check	05-19-15
Board of Tests Bench Check	08-02-17
Manufacturer Calibration	04-30-18
Board of Tests Bench Check	05-13-19



Presentation of this form certifies that the Intoxilyzer, SN: 80-003292, functioned properly at the time of the breath test and hereby certifies the breath test result as valid.

Date:

Re: Request for Log of Test on Intoxilyzer 8000

Name:

Dear Sir:

The Log of Test information from a breath test conducted on a Intoxilyzer 8000 is printed in Section 2 of the Affidavit and is no longer retained as a separate form as with the Intoxilyzer 5000. The 8000 uses a nitrogen/ethanol dry gas for the external control test. The dry gas cylinder has the capacity to perform approximately 325 tests and does not require maintenance until its expiration date or can no longer provide the required pressure to perform a control test. The cylinder has an expiration date of two years from date filled at the manufacturer and is programmed in the 8000 upon installation, both dates are printed in Section 2 of the Affidavit. The Intoxilyzer is programmed to disable itself when the date of expiration is reached or the pressure monitor indicates 25 lbs psi.

Due to the reliability and stability of dry gas, maintenance on the 8000 is not required every 30 days or 25 tests, therefore maintenance is no longer performed by officers in the field. The Board of Tests for Alcohol and Drug Influence certifies each Intoxilyzer 8000 for field use by performing a bench check. The bench check is designed to check the operation and procedures programmed in the 8000 to ensure it is functioning properly for use in Oklahoma. A bench check is performed every time a new dry gas cylinder is installed or the instrument has returned from the manufacturer for repair. Once the instrument passes the bench check, all that is required to set it up for field use is to plug it in and turn it on and does not require a permit to set up. The Board of Tests has the ability to monitor each instrument once it is set up in the field.

Each Intoxilyzer has a database of valid permit numbers that it checks prior to allowing a breath test to be performed. The database is updated on a regular basis, if a permit number is invalid that officer will not be able to conduct a breath test on the 8000. If the 8000 prints an affidavit the breath test operator listed in Section 2 had a valid permit at the time of the test and the breath test is valid.

If you have any further questions please feel free to call me.

Sincerely,

J. Robert Blakeburn (OBA11059) State Director of Tests for Alcohol and Drug Influence

O. pohr & Blakeline

Certificate of Calibration

This is to certify the calibration of **MONIVEY** serial number manufactured by CMI Inc., a subsidiary of MPD, Inc. of Owensboro, Kentucky, was tested and found to conform to the National Highway Traffic Safety Administration (NHTSA) Standard for Devices to Measure Breath Alcohol (Federal Register, Vol.58 No.179, pp 48705-48710, Sept. 17, 1993) for accuracy and precision. Reference materials are traceable through the National Institute of Standards and Technology (NIST) to the International System of Units (SI).

Date April 30, 2018 Signed Man Phill
Technician
INC.
316
Owe

316 East 9th Street Owensboro, KY 42303

Part No. 650517 Rev.A

STATE OF OKLAHOMA - BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE

ENCH CHECK DATE 05/13/2019	TIME PROCEDURE BEGAN 09:19:5		NSTRUMENT LOCATION 80-00329	2	
ntoxilyzer	INSTRUMENT MO			80-00	
4		TEST RE	SULTS		
Diagnostic PASS	es				
Mouth Alcoho	l Test				
RFI Test PASS					
Abort Tes		Du	plicate Breath	Sampling Re	esults
PASS		Value	Test 1	Test 2	Control Test
Improper Samp	le Test				
PASS		.02	0.019	0.019	0.018
Interferent Dete	ct Test	.04	0.039	0.038	0.037
PASS		.10	0.096	0.096	0.096
No Sample Give	en Test				
PASS		.20	0.192	0.192	0.193
Insufficient Samp	ple Test		PA	SS	
Air Blank Contamin	nation Test				
Range Exceede PASS	d Test				
	CE	ERTIFYING T	ECHNICIAN		
swear under penalty he Intoxliyzer Bench					
TURE				DATE	5-13-19
MITH, JOSHUA				PERMIT#	016837

End Time: 1030

3600 N. Martin Luther King Bldg #9 Oklahoma City, OK 73111 405-425-2460

BENCH CHECK TECHNICIAN'S REPORT

Start Time: 0900

Date: 05 / 13 / 19

000		Serial #				
000	80-003292					
SAS CANISTER						
OT#	- A - 1-	EXP Date				
99919080A1	# 040	08-05-2	021			
REFERENCE						
0 11 0400	And the second s	tor Model				
Guth 2100	Guth 2100	Guth 2100	Guth 2100			
DR3591	DR3753	erial # DR3594	DR3754			
DINOSSI		entration	DR3734			
.02	.04	.10	.20			
		T #				
18020	17410	18070	17340			
7.000	Manufac	tured Date				
01-09-2018	12-06-2017	02-26-2018	10-09-2017			
	Expira	tion Date				
01-09-2020	12-06-2019	02-26-2020	10-09-2019			
	Solution Cor	mmission Date				
05-06-19	05-06-19	15-06-19	05-06-19			
05-06-19	05-06-19	05-06-19	05-06-19			
05-06-19 REASON FOR BEN POST REPAIR TROUBLESHOOT COMMENTS	ICH CHECK CYLINDER REPL	ACEMENT	05-06-19			
REASON FOR BEN POST REPAIR DESTROUBLESHOOT	ICH CHECK CYLINDER REPL	ACEMENT	05-06-1			

Certificate	of Analysis – Dry Gas	
05/13/19	,080	
Date of Analysis	Labelled target value (g/210L)	
# 640	960	
Cylinder #	Average test result	
BOT Technician Name and Signatur	. K	specialty gases

In accordance with BOT Policy and Procedure Statement BRIE as Control of the Procedure Procedure 200 - Jacksonville, IL 62651-0790 referenced dry gas cylinder is suitable for use as an external reference with the stoward 2997-243-7634。www.ilmoproducts.com

Certificate of Analysis

Certificate ID:

11721

Part #:

BAC105L080T

Cylinder Size:

105L

Lot Number:

00919080A1

Expiration:

3/5/2021 4

0.080 BAC (For the calibration of instruments used to determine breath alcohol concentration)

Contents:

105 Liters @ 1000 psig 70°F (21°C)

Analytical

Reported

Accuracy

Analytical

Component:

Concentration:

(U, k=2):

Method:

Ethanol

208 ppm

+/-0.002 BAC(G/210L) NDIR

[5.2 ppm]

Nitrogen

Balance

CMI Inc. 316 East Ninth Stn Owensboro, KY 4 Phone 866-835-06 www.alcoholtest.c

Distributed by:

Operator's Signature

ALPHA INSTRUMENT

Model 8000

05/13/2019

Test

Air Blank

Cal Check

Air Blank

Std Deu

Cal Check Stats Average

Rel Std Deu(%)

Intoxilyzer - Alcohol Analyzer

g/210L

0.000

0.079

0.000

0.079

0.000

0.079

0.000

0.079

0.000

0.079

0.000

0.080

0.000

0.081

0.000

0.080

0.000

0.080

0.000

0.079

0.000

0.0795

0.0007

0.8894

SN 80-002591

09:28

Time

09:28

09:29

09:29

09:30

09:30

19:31

09:31

09:31

09:32

09:32

09:33

09:33

09:34

09:34

09:35

09:35

09:35

09:36

09:36

09:37

09:37

*NIST Traceable Reference Material Cylinder No. CC274523 / Job No. 09160306 Certified 362.2 µmol/mol Ethanol in Nitrogen Store in dry area, away from sources of heat, ignition and direct sunlight. Do not allow storage area to exceed 52 °C (125 °F).

Specialty Gas Lab Tech

Accreditation #61895

CMI INC. Service Evaluation Form

This form MUST be completed and enclosed with instrument to be serviced. Failure to complete and return this form may cause delays in service.

	•	ng container	r(s) or a similar protective box.)				
1. Contact information: Custom	ner Number 7	31362	(contact Customer Service)				
Name: Board of Tests		Phone: (405) 425 2460					
Fax: (405) 425 2490	Email: bot @	@dps.state.o	k.us				
2. Bill to Address: Board of Tests 3600 N ML King, OKC, OK 73111	•	to Address:	Board of Tests 3600 N ML King, Bldg 9 OKC, OK 73111				
3. Serial Number: 80-003292	Instru	ument Model:	Intoxilyzer 500 / 8000				
4. Detailed Description of Problem: DVM;s registering equally. Un	it experiencing o	cal check out	t of tolerance.				
Hazardous Material Warning! - <u>DO NOT</u> return gas cylinders with instrument!							
i azaradad material vyariini	g! - <u>DO NO I</u> retui	rn gas cylinde	ers with instrument!***				
☐ I Authorize Repairs Up To: ☐A		_					
_	II <u></u> \$250 ∑\$5	00	Other \$				
☐ I Authorize Repairs Up To: ☐A	II <u></u> \$250 ∑\$5	00	Other \$				
☐ I Authorize Repairs Up To: ☐A	II <u></u> \$250 ∑\$5	00	Other \$				
☑ I Authorize Repairs Up To: ☐A Purchase Order Number (attach	II <u></u> \$250 ∑\$5	00	Other \$				
✓ I Authorize Repairs Up To: ☐A Purchase Order Number (attach Authorized By:	II \$250 \$50 a copy of P.O. if a	00	Other \$				
 ✓ I Authorize Repairs Up To: ☐ A Purchase Order Number (attach Authorized By: Joshua Smith 	II \$250 \$50 a copy of P.O. if a Breath Alcohol Testing Program Admin	00	Other \$				
 ✓ I Authorize Repairs Up To: ☐ A Purchase Order Number (attach Authorized By: Joshua Smith 	II \$250 \$50 a copy of P.O. if a Breath Alcohol Testing Program Admin Title	00	Other \$				
I Authorize Repairs Up To: A Purchase Order Number (attach Authorized By: Joshua Smith Name (Please Print) Signature Call: 405-425-2460 Fax 405-42	Breath Alcohol Testing Program Admin Title 02/28/2018 Date	00 \$750 applicable) Ship CMI Attr 316	Other \$ to: I, Inc n: Service Dept.				
I Authorize Repairs Up To: A Purchase Order Number (attach Authorized By: Joshua Smith Name (Please Print) Signature Call: 405-425-2460 Fax 405-42 No, please send estimate before	Breath Alcohol Testing Program Admin Title 02/28/2018 Date 5-2490 repairs are made	00	Other \$ to: I, Inc n: Service Dept. S East Ninth Street				



SERVICE WORK ORDER

DATE 3 - 5 - 18

RCVD VIA 11PS

316 E. 9th St., Owensboro KY 42303 Phone 866-835-0690 Fax: 270-685-6268

401733

BILL TO:		SHIP TO	D:0K	Board of	LFL9T '	
ATTN:		ATTN:				
		3/000	N	me King		
		Bldo	19			
CITY STATE ZIP		CITY KC)	STATE	731	ZIP
CUSTOMER NO. 731362		PHONE	15-1	125-2460		
PHONE		FAX	05-	425-2490		
MODEL # 8000 SERIAL #80-003292		000	2481	OOK	BILL CODE	
EXTRA PARTS RCVD						
DESCRIPTION OF PROBLEM DV MS (*GISKring) equally, Pal charles	-\ -\	2 tolore	140.			
				PARTS	USED	
Infrared / Fuel Cell PRELIMINARY TEST INFORMATION DVM's both 13.5k, source OK. Breath O-ring cracked. Need CR25 update. Cell indet filter clogged.	_V; IN	05 CV	QTY	PART NO.	DESCRI	PTION
oring cracked. Need CR25 update.	VERIF	IED PROBLEM?	1	923002	Diode SA	no
WORK PERFORMED:		res	/	440904	Filter S	
Transfer to the second		a 9	1	402565	0-ring	
Replaced cracked o-ring and filt	er .	screen.			-	
Updated CRZS.						
Reset toreamo gains,						
Performed system and update ch	ec	£5.				
Alan Can At const	1		-			
Autocal performed to correct concect out of tolerance.						
Eneck out of Tolerance,						
PN						
Miscellaneous						
PRELIMINARY TEST INFORMATION						
WORK PERFORMED:						

REPAIR
CAL. ADJ.
CALIBRATON
FINAL

1/2 HOURS
HOURS
HOURS
DATE
DATE
DATE
DATE
DATE
DATE
DATE

SHIPPING STICKER HERE

SPECIAL SHIPPING INSTRUCTIONS

OFFICER'S AFFIDAVIT AND NOTICE OF

1		11	
	l /	-	h
	L/	1	v

05/13/2019	/	ARREST TIME 1000		COUNTY #	#	81		0	MOC			
RREST LOCATION MOCK					CITY	CKM				COUN	TY OCK	
BJECT NAME					DA.	TE OF B			SEX	RACE		WEIGH
OCK, SUBJE			CITY		C	08/0	9/1		. M STATE	M	510 ZIP CODE	
600 MARTIN	I LUTHER	KING A	OKLA	HOMA C		Ist	TATE	CLAS	OK	IRESTE	731	33 ENDORSEM
083454871 ICLE MAKE		MODEL		06/30		9 (OK	A		- [M
OCK		MOCE MOCE				TAG #				OK	N	PLACARD N
ee above date, time, and loci vay, street, turnpike, public p nce of alcohol and/or other it ribe driving behavior or circl DCK SUBJECT ibe person's condition (odor	parking lot, or other public intoxicating substances a umstances): TEST POST	c place, or upon a priva is prohibited by law.	le road, street, alley or la	ne which provides	access to o	one or mo	ore single	or multi-	family dwo	elling within	this state wh	ile under the
		THE PERSON	WAS INFORMED OF TI	HE IMPLIED CON	SENT TES	T REQUI	EST					
ATH-ALCOHOL ANALYSIS	RECORD, REPORT OF				ST AND MA	INTENA	NCE REC	CORD.				Secti
PMENT - 8000			NUMBER 5	SPECIALIST BOARD	OF T	EST	S					
IDARD LY			TARGET VALUE		PERMIT N	UMBER						
JFACTURER			0.080 LOT#		1683	INDER IN	NSTALLE	D	DA.	TE CYLIND	ER EXPIRES	<u> </u>
JMO			009190	80A1	05/1	13/2	2019)			5/202	
RVATION START TEST	DATE 5/13/2019	MOUTHPIECE G	A 15-minute	e (or greater)					T'			
LT TABLE	7/13/2019	_ J G		s observed p							75 0	
est	g/210L	Time	to th	ne time of tes	st admin	istratio	on.			- F.	THE	ATE
										4/*		1.161
r Blank	0.00	10:31	I administe								AND	
ıbject Test Breath Volume	0.00 1.691 LITERS	10:31	compliance	e with the	applica	able (opera	ting		8 **		
r Blank	0.00	10:32	proced	dure of the	Board	of T	ests.			[4] **		#(* * ¥)
it		10:34	1							13/1		
r Blank	0.00	10:34								1,		シゾノ
bject Test Breath Volume	0.00 1.613 LITERS	10:35	SIGNATURE OF C	OPERATOR							1907	
r Blank	0.00	10:35										
1 Check	0.08	10:36	NAME		-				Γ	,	Test Res	sult
r Blank	0.00	10:36	SMITH, JO BADGE#	OSHUA	PERMIT#				4		rAC g/2	
			00000004		0000	016	837	,			_	
			AGENCY	mpama					 	(0.0	U
CTIVE THIRTY (30) DAYS	S EROM DATE CERVE	ID Notice of Dec		TESTS			16.		<u> </u>			
or more, if you refused or t ating a CMV or CDL holde y must receive your writter	railed the state's test(s) r) OR a test result of .0 n request for an administ t form is available at Safety, P.O. Box 11415	 In addition, your co 4 or more while opera strative hearing on th https://www.ok.gov/ , OKC, OK 73136. It 	mmercial driving privile ating a commercial mot e Department approve dps/documents/ReqH shall be a misdemeand and imprisonment, for	ge in this state w or vehicle. To ap ed form, within fi learing.pdf. Sub or punishable by i	rill be disqua opeal the dr ifteen (15) o omit your he imprisonme	alified if river's lic days from earing re	you refusense revenue the da equest for the set the se	sed or fa ocationate you re rm in pe	ailed the disqualifeceived reson at 3	state's test ication, the lotice of re 600 N. ML	t(s) e Departmen vocation or K Ave. or by	disqualificati mail to Leg
approved hearing reques on, Department of Public S of not more than Five Hur n's license, permit, or othe	r evidence of driving p	rivilege is in the custo		nt officer or the D	epartment.							
pproved hearing reques in, Department of Public S of not more than Five Hur 's license, permit, or othe VING OF NOTICE: I	r evidence of driving pr personally hand-d	rivilege is in the custo elivered a copy of	the above containi	nt officer or the Di ing the Notice	epartment. of Revoc	cation t	o the p	erson	named	in Sectio	on 1 above	€.
oproved hearing reques n, Department of Public S of not more than Five Hur t's license, permit, or othe	r evidence of driving pr personally hand-d	rivilege is in the custo elivered a copy of	the above containi	nt officer or the Di ing the Notice	epartment. of Revoc	cation t	o the p	erson	named Ba	in Sectionadge #:_	on 1 above	
oproved hearing reques in, Department of Public S of not more than Five Hur i's license, permit, or othe //ING OF NOTICE: I	personally hand-dispersion of driving purposes of the personal purposes	rivilege is in the custo elivered a copy of ignature of officer , "I state under pe	the above containi	nt officer or the Di ing the Notice	of Revoc	cation t	o the p	erson	Ba	in Sectionadge #: ue and co	on 1 above	Section
pproved hearing reques in, Department of Public S of not more than Five Hur i's license, permit, or othe VING OF NOTICE: I served: 05/13 cordance with Title 12	personally hand-diagrams of 2019 S 2 O.S. Section 426 Place (location	rivilege is in the custo elivered a copy of ignature of officer , "I state under pe	the above containi	ing the Notice Inder the laws	of Revoc of Oklaho Signatu	oma thure of a	at the f	erson	Ba	in Section adge #:_ ue and co	on 1 above	
pproved hearing reques in, Department of Public S of not more than Five Hur i's license, permit, or othe VING OF NOTICE: I served: 05/13 cordance with Title 12 05/13/2019	personally hand-diagrams of 2019 Section 426 Place (location HUA	rivilege is in the custo elivered a copy of ignature of officer , "I state under pe	the above containi	nt officer or the Ding the Notice	of Revoc of Oklaho Signatu	oma thure of a	at the f	erson oregoi	Ba	in Sectional Section Section 19 and control of the section 19 and	on 1 above	
proved hearing reques in, Department of Public S of not more than Five Hur i's license, permit, or othe //ING OF NOTICE: I served: 05/13 cordance with Title 12 05/13/201	personally hand-diagrams of 2019 S 2 O.S. Section 426 Place (location	ivilege is in the custo elivered a copy of ignature of officer , "I state under pe when signed):	the above containi	ing the Notice Inder the laws	of Revoc of Oklaho Signatu	oma thure of a	at the f	oregoi	Ba	adge #:_ ue and co	on 1 above	
proved hearing reques n, Department of Public S of not more than Five Hur 's license, permit, or othe //ING OF NOTICE: I served: 05/13 cordance with Title 12 05/13/201 CITH, JOS: 00000435 ER WITNESSES:	personally hand-diagrams of 2019 S 2 O.S. Section 426 9 Place (location HUA AGENCY ADDRESS	ivilege is in the custo elivered a copy of ignature of officer , "I state under pe when signed):	the above containi	ing the Notice Inder the laws	of Revoc of Oklaho Signatu	oma thure of a	at the f	oregoi	Ba	adge #:_ ue and co	orrect."	
pproved hearing reques in, Department of Public S of not more than Five Hur n's license, permit, or othe VING OF NOTICE: I served: 05/13 Cordance with Title 12 05/13/2019	personally hand-diagrams of 2019 S 2 O.S. Section 426 9 Place (location HUA AGENCY ADDRESS	ivilege is in the custo elivered a copy of ignature of officer , "I state under pe when signed):	the above containi	ing the Notice Inder the laws	of Revoc of Oklaho Signatu	oma thure of a	at the f	oregoi	Ba	adge #:_ ue and co	orrect."	Section

12/16

OFFICER'S AFFIDAVIT AND NOTICE OF REVOCATION/DISOUALIFICATION

1	E / OCMITO	TANDID		/A.A. A		LOI					
			_								Section 1
ARREST DATE	ARREST TIME		COUNTY#		CITY#	CI	TATION#				
05/13/2019	1000		55		81	ľ	10CF	ζ			
ARREST LOCATION				CITY				COUNT	Y		
MOCK				MOC	CKM			MO	CK		
SUBJECT NAME					E OF BIRTH		SEX	RACE	HEIGHT	W	/EIGHT
MOCK, SUBJECT TEST				0	8/09/:	1961	M	W	510) :	165
ADDRESS		CITY				ST	ATE		ZIP COL	E	
3600 MARTIN LUTHER	KING A	OKLAHO	MA CI	TY			ΣK		731	L33	
DRIVER LICENSE #			RATION DATE		STATE	CLASS		RESTRI	CTIONS	ENDO	RSEMENTS
T083454871		0.	6/30/:	2019	9 OK	A				M	
VEHICLE MAKE	MODEL				TAG#		ST	ATE	CMV	PLACA	\RD
MOCK	MOCK				MOKC			ΣK	N	N	
On the above date, time, and location, the above named highway, street, turnpike, public parking lot, or other publinfluence of alcohol and/or other intoxicating substances (Describe driving behavior or circumstances):	lic place, or upon a private road, stre	esonable grounds to t eet, alley or lane whic	pelieve the per ch provides acc	son was o	driving or in actu ne or more singl	ual physical c e or multi-far	ontrol of nily dwelli	a motor ve ng within t	hicle upon his state w	a public	c road, der the

MOCK SUBJECT TEST. POST REPAIR. NEW GAS CANISTER INSTALL

Describe person's condition (odor, actions, etc.):

NONE

THE PERSON WAS INFORMED OF THE IMPLIED CONSENT TEST REQUEST

		THETERSON	VAS INFORMED OF THE IMPLIE	CONSENT TEST REQUEST	
BREATH-ALCOHOL ANALYSIS	RECORD, REPORT OF TH	IE PERSON NAMED	IN SECTION 1 ABOVE, and LOG	OF TEST AND MAINTENANCE RECORD.	Section 2
EQUIPMENT		SERIAL	NUMBER SPECIALIST		
I-8000		80-	-003292 BOAR	D OF TESTS	
STANDARD			TARGET VALUE	PERMIT NUMBER	
Dry			0.080	16837	
MANUFACTURER			LOT#	DATE CYLINDER INSTALLED	DATE CYLINDER EXPIRES
ILMO			00919080A1	05/13/2019	03/05/2021
	DATE	MOUTHPIECE	A 15-minute (or gre	ater) deprivation/observation	T
	5/13/2019	G			
RESULT TABLE				ed prior to and continuing	THE STAR
Test	g/210L	Time	to the time of	f test administration.	0° ** \^* * * 0
Air Blank	0.00	10:31	I administered sai	d Breath-Alcohol Test in	V ** 10
Subject Test	0.00	10:31	compliance with	the applicable operating	10 + 10 × 10 × 10
Breath Volume	1.691 LITERS		_		日 ** 20
Air Blank	0.00	10:32	procedure of	the Board of Tests.	A + 15
Wait		10:34	_		C) */
Air Blank	0.00	10:34			
Subject Test	0.00	10:35			1907
Breath Volume	1.613 LITERS		SIGNATURE OF OPERATOR		1.507
Air Blank	0.00	10:35			
Cal Check	0.08	10:36	NAME		Test Result
Air Blank	0.00	10:36	SMITH, JOSHUA		BrAC g/210L
			BADGE #	PERMIT#	DIAC 9/210L
			0000000435	0000016837	
			AGENCY	7.0	
			BOARD OF TEST	.'S	

EFFECTIVE THIRTY (30) DAYS FROM DATE SERVED. Notice of Revocation is given that your driver's license is revoked or denied for a period of one hundred eighty (180)
days or more, if you refused or failed the state's test(s). In addition, your commercial driving privilege in this state will be disqualified if you refused or failed the state's test(s) (operating a CMV or CDL holder) OR a test result of .04 or more while operating a commercial motor vehicle. To appeal the driver's license revocation/disqualification, the Department of Public Safety must receive your written request for an administrative hearing on the Department approved form, within fifteen (15) days from the date you received notice of revocation or disqualification. The approved hearing request form is available at https://www.ok.gov/dps/documents/ReqHearing.pdf. Submit your hearing request form in person at 3600 N. MLK Ave. or by mail to Legal Division, Department of Public Safety, P.O. Box 11415, OKC, OK 73136. It shall be a misdemeanor punishable by imprisonment for not less than seven (7) days nor more than six (6) months, or by a fine of not more than Five Hundred Dollars (\$500.00), or by both such fine and imprisonment, for any person to apply for a renewal or a replacement license to operate a motor vehicle while the person's license, permit, or other evidence of driving privilege is in the custody of a law enforcement officer or the Department.

SERVING OF NOTICE: I personally hand-delivered a copy of the above containing the Notice of Revocation to the person named in Section 1 above.

Date served: 05/13/2019 Signature of officer:	Badge #:
In accordance with Title 12 O.S. Section 426, "I state under penalty of perjury to	under the laws of Oklahoma that the foregoing is true and correct." Section 4
Date: 05/13/2019 Place (location when signed):	Signature of arresting officer:
SMITH, JOSHUA	BOARD OF TESTS
BADGE O 0 0 0 0 0 4 3 5 AGENCY ADDRESS INTER-AGENCY	ZIP CODE PHONE 00000 "4054252460"



Certificate of Analysis – Wet Bath

09-2b-2018 , 82 9 3/16

Date of Analysis Labelled target value (g/210L)

18029 0.02 9 4

S. Nurph C. Average test result-(g/210L)

2.5

BOT Technician Name and Signature

In accordance with BRT - 2.1.0, Maintenance, the above referenced simulator solution is

CERTIFICATE OF ANALYSIS

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number 18020 of

Alcohol Reference Solution for Simulator were analyzed by
gas chromatography on January 11, 2018, using a Perkin Elmer Gas
Chromatograph Autosystem XL S/N: 610N9030209, and found to contain
0.0238% (w/vol) ethyl alcohol. The expiration date for this lot
number is January 9, 2020 at 11:59 PM.

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of 0.020 g/210L +/- 3%.

The alcohol and water used in this solution were free of test interfering substances.

Signature

RLPHR INSTRUMENT

Model 8000

09/26/2018

Rir Blank

Cal Check Air Blank

Cal Check

Air Blank

Cal Check

Rir Blank Cal Check

Air Blank

Cal Check

Air Blank

Cal Check

Air Blank

Tal Check

Rir Blank

Cal Check

Pir Blank

'al Check

Hr Blank

al Check

ir Blank

5td Deu

al Check Stats Average

Rel Std Deu(%)

Intoxilyzer - Alcohol Analyzer

g/210L

.0.000

0.020

0.000

0.020

0.000

0.020

0.020

0.000

0.020

0.000

0.020

0.000

0.021

0.000

0.021

0.000

0.021

0.000

0.021

0.000

0.0204

0.0005

2.5312

SN 80-002591

15:36

Time

15:36

15:37

15:37

15:38 15:39

15:39

15:40

15:41

15:41

15:42

15:42

15:43

15:44

15:44

15:45

15:46

15:46

15:47

15:47

15:48

Ted L. Pauley, President GUTH LABORATORIES, INC.

NIST Traceability:

Testing was conducted using Cerilliant Reference Standard lot number FN08031603 whose values are traceable to NIST.

All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.



Certificate of Analysis - Wet Bath 09-26-18 Labelled target value (g/210L) Date of Analysis 0.0396



In accordance with BRT - 2.1.0, Maintenance, the above referenced simulator solution is A sittle for use of the intoxilyzer 8000.

Average test result (g/210L)

CERTIFICATE OF

Certified Alcohol Reference Solution for Simulator

ALPHA INSTRUMENT Intoxilyzer - Alcohol Analyzer Model 8000 SN 80-002591 15:54 g/210L Time

09/26/2018 Test Air Blank 0.000 15:54 Cal Check 15:55 Air Blank 0.000 15:56 Cal Check 0.039 15:56 Air Blank 0.000 15:57 Cal Check 0.03915:58 Rir Blank 0.000 15:58 Cal Check 15:59 Air Blank 0.000 15:59 Cal Check 0.040 16:00 Air Blank 0.000 16:01 Cal Check 0.040 16:01 Air Blank 0.000 16:02 Cal Check 0.039 16:03 Air Blank 0.000 16:03 Cal Check 0.040 16:04 Air Blank 0.000 16:04 Cal Check 0.040 16:05 Air Blank 0.000 16:06 Cal Check 0.039 16:06 Air Blank 0.000 Cal Check Stats Average Std Dev 0.0005 Rel Std Deu(2) 1.3040

Random Samples of Lot Number Alcohol Reference Solution for Simulator were analyzed by chromatography on December 7, 2017, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain 0.0481% (w/vol) ethyl alcohol. The expiration date for this lot number is December 6, 2019 at 11:59 PM.

When used in a calibrated Simulator, operating at +/- .2°C, this solution will give a breath alcohol analysis instrument reading of 0.040 g/210L +/- 3%.

The alcohol and water used in this solution were free of test interfering substances.


Operator's Signature

Ted L. Pauley, President GUTH LABORATORIES, INC.

NIST Traceability:

Testing was conducted using Cerilliant Reference Standard lot number FN12181501 whose values are traceable to NIST.

All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.

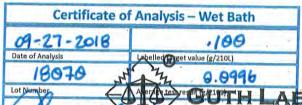

BOT Technician Name and Signatur

ALPHR INSTRUMENT

Model 8000

09/27/2018

Test

Air Blank

Cal Check

Air Blank

-Cal Check

Air Blank

Cal Check

Air Blank

Call Check

Rir Blank

Std Deu Rel Std Deu(%)

Cal Check Stats Average

Operator

Intoxilyzer - Alcohol Analyzer

g/210L

0,000

0.100

0.000

0.000

0.099

0.000

0.100

0.000

0.100

0,000

0.100

0.000

0.099

0.000

0.099

0.000

0.099

0.000

0.100

0.000

0.0996

0.5184

SN 80-002591

08:28

Time

08:28

08:29

08:29

08:30

08:31

08:31

08:32

08:33

08:33

08:34

08:34

08:35

08:36

08:36

08:37

08:38

08:38

08:39

08:39

08:40

08:41

590 NORTH 67th STREET A HARRISBURG, PA 17111-4511 TELEPHONE: 717-564-5470

BRT Form 6.0

In actordance with BRT - 2.1.0, Maintenance, the above referenced simulator solution is suitable for use as an external reference in maintenance of the Intoxilyzer 8000.

ERTIFICATE OF ANALYSIS

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number 18070 Alcohol Reference Solution for Simulator were analyzed by chromatography on February 28, 2018, using a Perkin Elmer Gas gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain 0.1219% (w/vol) ethyl alcohol. The expiration date for this lot number is February 26, 2020 at 11:59 PM.

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of 0.100 g/210L +/- 3%.


The alcohol and water used in this solution were free of test interfering substances.

Ted L. Pauley, President GUTH LABORATORIES, INC.

NIST Traceability:

Testing was conducted using Cerilliant Reference Standard lot number FN04271602 whose values are traceable to NIST.

All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.

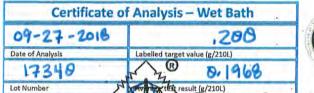

GUTH LABORATORIES, INC

HARRISBURG, PA 17111- 4511 ® TELEPHONE: 717-564-5470

BRT Form 6.0

In accordance with BRT -2.1.0, Majolenance, the above referenced simulator solution is suitable for use as an external reference in maintenance of the Intoxilyzer 8000.

CERTIFICATE OF ANALYSIS

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number 17340 of
Alcohol Reference Solution for Simulator were analyzed by
gas chromatography on October 10, 2017, using a Perkin Elmer Gas
Chromatograph Autosystem XL S/N: 610N9030209, and found to contain
0.2410% (w/vol) ethyl alcohol. The expiration date for this lot
number is October 9, 2019 at 11:59 PM.

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of 0.200 g/210L +/- 3%.

The alcohol and water used in this solution were free of test interfering substances.

Operator's Xignature

ALPHA INSTRUMENT

Model 8000

09/27/2018

Air Blank

Cal Check

Rir Blank

:al Check

ir Blank

al Check

lir Blank

Std Deu

ial Check Stats Average

Rel Std Dev(%)

Intoxilyzer - Alcohol Analyzer

g/210L

0.000

0.196

0.00

0,196

0.000

0.197

0.000

0.197

0.000

0.197

0.000

0.197

0.000

0.197

0.000

0.197

0.000

0.197

0.000

0, 197

0.000

0.1968

0.0004

0.2142

SN 80-002591

15:27

Time

15:27

15:28

15:29

15:29

15:30

15:31

15:31

15:32

15:33

15:33

15:34

15:35

15:35

15:36

15:36

15:37

15:38

15:38

15:39

15:40

Ted L. Pauley, President GUTH LABORATORIES, INC.

NIST Traceability:

Testing was conducted using Cerilliant Reference Standard lot number FN08101505 whose values are traceable to NIST.

All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.