

STATE OF OKLAHOMA BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE

Dr. Jarrad Wagner, Ph.D. F-ABFT Chairman

Post Office Box 36307 Oklahoma City, Oklahoma 73136-2307 Phone: (405) 425-2460 Fax: (405) 425-2490 www.bot.ok.gov

Joshua Smith Director

ATTESTATION

I, the undersigned <u>Joshua Smith</u> (Director/Records Custodian), attest under oath pursuant to 12 O.S. §2902 (11)(a) this record was made, at or near the time of the occurrence of the matters set forth by or from information transmitted by a person having knowledge of those matters; was kept in the course of the regularly conducted business activity of the Board; and was made pursuant to the regularly conducted activity. The duplicate record provided is listed below:

1) The attached "Service Overview" is a true and accurate copy of maintenance performed on the Intoxilyzer 8000, serial number 80-003323, in compliance with the agency's Maintenance Policy and Procedures.

This packet consists of 13 pages.

Name of Position: Director/Records Custodian

Date of Attestation: 05/18/2022

Service Overview



3600 N. Martin Luther King Bldg #9 Oklahoma City, OK 73111 405-425-2460

BENCH CHECK DATE: 05-06-2019

SERIAL: <u>80-003323</u>

DOCUMENTATION WITHIN SERVICE PACKET:

- BOT Certificate of Calibration and Operation
- Bench Check
- Bench Check Technician's Report
- BOT Dry Gas Certificate of Analysis
- Mock Subject Test
- BOT Certificates of Analysis on Certified Solutions (4)



STATE OF OKLAHOMA BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE

Administrative Office at:
Board of Tests for Alcohol and Drug Influence
Post Office Box 36307
Oklahoma City, Oklahoma 73136-2307

Certificate of Calibration and Operation

This is to certify that the calibration of INTOXILYZER 8000, serial number 80-003323, manufactured by CMI Inc., a subsidiary of MPD, Inc. of Owensboro, Kentucky, was tested and found to be in compliance with National Highway Traffic Safety Administration Standard for Devices to Measure Breath Alcohol (F.R., No. 179 48705-48710) Sept. 17, 1993) and with International Recommendation OIML R 126 Edition 1998 (E). Calibration solutions are traceable to NIST standard material 1828.

In addition the operational procedures and functions of this Intoxilyzer have also been tested and found to be in compliance with the Oklahoma Board of Tests for Alcohol and Drug Influence - Rules and Regulations OAC Title 40: 25-1-2 and 30-1-3 and is hereby certified for use as an evidential breath testing device within this state.

CERTIFICATION TYPE	DATE
Manufacturer Calibration	07-21-08
Board of Tests Bench Check	10-10-08
Board of Tests Bench Check	03-16-10
Manufacturer Calibration	03-28-12
Board of Tests Bench Check	05-23-12
Manufacturer Calibration	09-13-13
Board of Tests Bench Check	09-18-13
Board of Tests Bench Check	12-09-14
Board of Tests Bench Check	09-11-15
Board of Tests Bench Check	06-16-17
Board of Tests Bench Check	05-06-19
	L



Presentation of this form certifies that the Intoxilyzer, SN: 80-003323, functioned properly at the time of the breath test and hereby certifies the breath test result as valid.

Date:

Re: Request for Log of Test on Intoxilyzer 8000

Name:

Dear Sir:

The Log of Test information from a breath test conducted on a Intoxilyzer 8000 is printed in Section 2 of the Affidavit and is no longer retained as a separate form as with the Intoxilyzer 5000. The 8000 uses a nitrogen/ethanol dry gas for the external control test. The dry gas cylinder has the capacity to perform approximately 325 tests and does not require maintenance until its expiration date or can no longer provide the required pressure to perform a control test. The cylinder has an expiration date of two years from date filled at the manufacturer and is programmed in the 8000 upon installation, both dates are printed in Section 2 of the Affidavit. The Intoxilyzer is programmed to disable itself when the date of expiration is reached or the pressure monitor indicates 25 lbs psi.

Due to the reliability and stability of dry gas, maintenance on the 8000 is not required every 30 days or 25 tests, therefore maintenance is no longer performed by officers in the field. The Board of Tests for Alcohol and Drug Influence certifies each Intoxilyzer 8000 for field use by performing a bench check. The bench check is designed to check the operation and procedures programmed in the 8000 to ensure it is functioning properly for use in Oklahoma. A bench check is performed every time a new dry gas cylinder is installed or the instrument has returned from the manufacturer for repair. Once the instrument passes the bench check, all that is required to set it up for field use is to plug it in and turn it on and does not require a permit to set up. The Board of Tests has the ability to monitor each instrument once it is set up in the field.

Each Intoxilyzer has a database of valid permit numbers that it checks prior to allowing a breath test to be performed. The database is updated on a regular basis, if a permit number is invalid that officer will not be able to conduct a breath test on the 8000. If the 8000 prints an affidavit the breath test operator listed in Section 2 had a valid permit at the time of the test and the breath test is valid.

If you have any further questions please feel free to call me.


Sincerely,

J. Robert Blakeburn (OBA11059) State Director of Tests for Alcohol and Drug Influence

O. pohr &Blokelune

STATE OF OKLAHOMA - BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE CMI INTOXILYZER BENCH CHECK PROCEDURE REPORT


	WI INTOXILYZER E			EDURE RE	PORT
15/06/2019	TIME PROCEDURE BEGAN 12:51:4		NSTRUMENT LOCATION 80-00332	3	
RUMENT TYPE ntoxilyzer	INSTRUMENT MO			80-00	
		TEST RE	ESULTS		
Diagnosti PASS					
Mouth Alcoho					
RFI Tes PASS		D.	unlineste Duneste	Canadia a Da	
Abort Te	est i	Du	plicate Breath	Sampling Re	suits
PASS		Value	Test 1	Test 2	Control Test
Improper Sam	ple Test				
PASS		.02	0.020	0.020	0.020
Interferent Det	ect Test	.04	0.039	0.039	0.040
PASS		.10	0.097	0.098	0.100
No Sample Giv	ven Test				
PASS		.20	0.195	0.196	0.199
Insufficient Sam PASS			PA	SS	
Air Blank Contami PASS					
Range Exceed PASS					
			TECHNICIAN		
he Intoxliyzer Benc	y of perjury that in re h Checking Procedur				
NTURE	J			DATE	15/06/2019
MITTIL TOCKET				PERMIT #	01.602.7
MITH, JOSHUA	<i>A</i>				016837

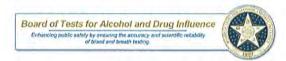

3600 N. Martin Luther King Bldg #9 Oklahoma City, OK 73111 405-425-2460

BENCH CHECK TECHNICIAN'S REPORT

00		Serial # 80-00 3323	
		600000000	
S CANISTER		VD Date	
3519000A5	# 019	OH-06-2	07.1
FERENCE	30011	04-09- 2	
FERENCE	Simulato	Model	
Guth 2100	Guth 2100	Guth 2100	Guth 2100
Guil 2100	Seria		Guil 2100
DR3591	DR3753	DR3594	DR3754
	Concen		
.02	.04 LO	.10	.20
18020	17410	18070	17340
10020	Manufactu	7 - 1 - 1 - 1 - 1 - 1 - 1 - 1 - 1 - 1 -	17040
01-09-2018	12-06-2017	02-26-2018	10-09-2017
707.70.20.10	Expiration		18.55 55 11
01-09-2020	12-06-2019	02-26-2020	10-09-2019
V(4) (V(2) (2) (1) (1)	Solution Com	mission Date	
	05-060 19	05/06/19	06/11/19
me - M - 19		09100117	05/06/19
05-06-19			
ASON FOR BEI		CEMENT	
ASON FOR BEI	NCH CHECK CYLINDER REPLA	CEMENT	

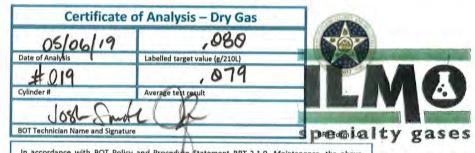

In accordance with BOT Policy and Procedure Statement BRT-2.1.0, Maintenance, the above referenced dry gas cylinder is suitable for use as an external/rearing who மிப்பால் இடு a Boom 790 + Jacksonville, IL 62651-0790 217-245-2183 • Fax: 217-243-7634 • www.ilmoproducts.com

Certificate of Analysis

PECN 22/1

Certificate ID:

11790

Part #:

BAC105L080T

Cylinder Size:

105L

Lot Number: Expiration:

03519080A5

4/5/2021

().080 BAC (For the calibration of instruments used to determine breath alcohol concentration)

Contents:

105 Liters @ 1000 psig 70°F (21°C)

Analytical

Reported

Analytical Accuracy

Component:

Concentration:

(U, k=2):

Method:

Ethanol Nitrogen 208 ppm Balance

*NIST Traceable Reference Material

Cylinder No. CC274523 / Job No. 09160306

Certified 362.2 µmol/mol Ethanol in Nitrogen

+/-0.002 BAC(G/210L)

[5.2 ppm]

Distributed by:

CMI Inc.

316 East Ninth Street Owensboro, KY 4230 Phone 866-835-0690

www.alcoholtest.com

Store in dry area, away from sources of heat, ignition and direct sunlight. Do not allow storage area to exceed 52 °C (125 °F).

Specialty Gas Lab Tech

Accreditation #61895

The calibration results within this certificate were obtained using equipment and standards capable of producing analytical results traceable to NIST, and apply only to the items contained on this certificate. ILMO Products Company makes no warranty or representation as to the suitability of the use of any information provided for any particular purpose. The information use is at the sole discretion and risk of the user. Liability shall be limited to established replacement cost of this material or service.

05/06/2019 Test q/210L Time Air Blank 0.000 12:52 Cal Check 0.079 12:53 0.000 Air Blank 12:53 Cal Check 0.080 12:54 Air Blank 0.000 12:54 Cal Check 0.080 12:54 Air Blank 0.000 12:55 Cal Check 0.079 12:55 Air Blank 0.000 12:56 Cal Check 0.079 12:56 Air Blank 0.000 12:57 0.079 Cal Check 12:57 0.000 12:57 Air Blank 0.080 12:58 Cal Check 0.000 Air Blank 12:58 0.080 12:59 Cal Check 0.000 Air Blank 12:59 Cal Check 0.079 13:00 Air Blank 0.000 13:00 Cal Check 0.079 13:00 Air Blank 0.000 13:01 Cal Check Stats 0.0794 Average Std Deu 0.0005

ALPHA INSTRUMENT

Model 8000

Rel Std Deu(%)

Intoxilyzer - Alcohol Analyzer

SN 80-002591

12:52

berator's Signatur

0.6503

8/13

OFFICER'S AFFIDAVIT AND NOTICE OF REVOCATION/DISQUALIFICATION

				_						Section
05/06/2019		ARREST TIME 1230		COUNTY 55		81	CITATIO			
ARREST LOCATION		1230		33	CITY		11100	COUNT		
MOCK			A		MOC			MO		
MOCK, SUBJE	CT TEST					OF BIRTH 3/09/196	1 N		HEIGHT 510	WEIGHT 165
ADDRESS 3600 MARTIN		KING A	CITY OI	KLAHOMA C	CITY	<u> </u>	STATE		ZIP COD	
DRIVER LICENSE # T083454871			l	EXPIRATION D)/2019	OK CI	ASS 7	RESTRI	CTIONS	ENDORSEMENTS
VEHICLE MAKE		MODEL		106/30	<u> </u>	AG#	Α	STATE	CMV	M PLACARD
MOCK		MOKC			1	MOCK		OK	N	N
On the above date, time, and loca highway, street, turnpike, public p. influence of alcohol and/or other ir (Describe driving behavior or circu MOCK SUBJECT	arking lot, or other pub ntoxicating substances imstances): TEST. ROU'	lic place, or upon a private as prohibited by law.	road, street, alle	ey or lane which provides	s access to one	or more single or m	sical contro ulti-family d	or a motor ve	nicie upon his state wl	a public road, nile under the
Describe person's condition (odor,	, actions, etc.):									
NONE										
		THE BEDGON I	AC INCODMED	OF THE IMBLIED CO	NEEDIT TEST I	POURCE				
BREATH-ALCOHOL ANALYSIS	RECORD REPORT O			OF THE IMPLIED CO			`			Section 2
EQUIPMENT		SERIAL	NUMBER	SPECIALIST						Section 2
I-8000		80-	003323		OF TE					
STANDARD Dry			O.08		PERMIT NUM					
MANUFACTURER			LOT#	1000075	DATE CYLINI	DER INSTALLED		ATE CYLINDE		
OBSERVATION START TEST	DATE	MOUTHPIECE		19080A5		6/2019		04/05	/202	1
	/06/2019			inute (or greater			n			_
RESULT TABLE			perio	d was observed to the time of te	•				THE S	TATA
Test	g/210L	Time		to the time of te	st aummins			(St.)		
Air Blank	0.00	13:39	I admin	istered said E	Breath-A	Icohol Test	in			
Subject Test	0.00	13:41		ance with the			- 1	5		多門
Breath Volume	2.125 LITER		-	ocedure of the		-	6	[4] ***		#(*** A
Air Blank Wait	0.00	13:41 13:43	Pι	occurre or the	c Doard	or rests.				
Air Blank	0.00	13:43						10%	$\mathcal{L}^{(3)}$	
Subject Test	0.00 1.750 LITER	13:44	SIGNATUR	RE OF OPERATOR			-		1907	
Breath Volume Air Blank	0.00	13:44								
Cal Check	0.07	13:44	NAME	TOCHILA				7	Γest Re	sult
Air Blank	0.00	13:45	SMITH, BADGE#	JOSHUA	PERMIT#			Bı	rAC g/	210L
			00000	00435	0000	016837		(0.0	0 (
			AGENCY BOARD	OF TESTS				•	<i>y</i> . c	, 0
EFFECTIVE THIRTY (30) DAY: days or more, if you refused or	S FROM DATE SER'	VED. Notice of Revoca	tion is given th	at your driver's license	is revoked or will be disqual	denied for a period	d of one hu or failed th	ndred eighty e state's test((180) (s)	Section 3
(operating a CMV or CDL holde Safety must receive your writter The approved hearing reques Division, Department of Public 8 a fine of not more than Five Hur person's license, permit, or other	r) OR a test result of a request for an admi at form is available a Safety, P.O. Box 114 adred Dollars (\$500.0	.04 or more while opera nistrative hearing on the at https://www.ok.gov/o 15, OKC, OK 73136. It s 00), or by both such fine	ing a commerce Department a ps/documents hall be a misde and imprisonme	ial motor vehicle. To a pproved form, within a large per selection of the l	appeal the driv fifteen (15) da abmit your hea apply for a ren	rer's license revoca ays from the date youring request form in t for not less than s	tion/disqua ou received o person al even (7) d	lification, the d notice of rev : 3600 N. MLI ays nor more	Departme ocation of Ave. or the	disqualification. by mail to Legal b) months, or by
SERVING OF NOTICE:				ontaining the Notic	e of Revoca	ation to the pers	on name	d in Sectio	n 1 abov	e.
Date served: 05/06	5/2019	Signature of officer						Badge #:		
In accordance with Title 1:	2 O.S. Section 42	26, "I state under pe	nalty of perj	ury under the law	s of Oklaho	ma that the fore	going is	true and co	orrect."	Section 4
Date: 05/06/201	9 Place (location	on when signed):			_ Signatur	e of arresting of	ficer:			
NAME SMITH, JOS				AGENCY BOARD	OF T	ESTS				
BADGE 000000435	AGENCY ADDRESS INTER-A	GENCY		BOARD	OF I	ZIP CODE 00000	PHONE	54252	14601	1
	7141 D.	<u> </u>				100000	1 +0	J4434	. 100	
OTHER WITNESSES:			TITLE	ADDRESS			IPHC	NF		Section 5
1							1,10	F 4 No.		
NAME 2			TITLE	ADDRESS			PHC	NE		
4							1			

9/13

OFFICER'S AFFIDAVIT AND NOTICE OF REVOCATION/DISQUALIFICATION

	1	E V OCITI		, 1 , 1	3Q0 1 1				011				S	ection 1
05/06/2019		ARREST TIME 1230			COUNTY #	!	CITY 8			MOC				
ARREST LOCATION		1230				CITY				1100	COUN			
MOCK SUBJECT NAME							OCK	BIRTH		SEX		OCK HEIGHT	WEIG	3HT
MOCK, SUBJEC	CT TEST	***************************************						09/1		M		510) 10	65
3600 MARTIN	LUTHER	KING A		OKLA:	HOMA C				(oK		731	.33	
DRIVER LICENSE # T083454871					06/30			STATE OK	CLASS	5	REST	RICTIONS	ENDORSI M	EMENTS
VEHICLE MAKE MOCK		MOKC		***************************************		•	TAG	# OCK			STATE OK	CMV N	PLACARE N	
On the above date, time, and locatic highway, street, turnpike, public part influence of alcohol and/or other into (Describe driving behavior or circum	king lot, or other pub exicating substances estances):	l person was arrested, and lic place, or upon a private as prohibited by law.	road, street,	alley or lan	e which provides	access to	is driving	g or in actua more single	or multi-fa	control imily dw	of a motor	vehicle upon	a public ro	pad, the
MOCK SUBJECT T	EST. ROU	TINE MAINTE	NANCE	. NEW	GAS CA	NISTI	ER .	LNSTA	ידיו					
Describe person's condition (odor, a	ctions, etc.):													
NONE														
DDEATH ALCOHOL ANALYSIS DE	CORD REPORT O				E IMPLIED CON				2000				C.	ti 2
BREATH-ALCOHOL ANALYSIS RE EQUIPMENT			NUMBER		PECIALIST	51 AND M	1AIN I EI	MANCE REC	JUND.					ection 2
I-8000		80-	0033		BOARD		ГES							
STANDARD Dry				T VALUE		PERMIT N		R						
MANUFACTURER			LOT#	······································		DATE CY	LINDEF	RINSTALLE				DER EXPIRE		
OBSERVATION START TEST DA	ATF	MOUTHPIECE			80A5	<u> </u>		/2019		()4/0.	5/202	21	
	⁾ 06/2019				(or greater)									
RESULT TABLE			pe		s observed p e time of tes				ıg			THE S	TAT	
Test	g/210L	Time		to the	e time of tes	st aumm	IIISHIa			_	19	***		
Air Blank	0.00	13:39	I adm	inister	ed said B	reath-	Alco	ohol Te	est in				/* * * *	0
_	0.00	13:41	1		with the						8 1		AF	
Breath Volume Air Blank	2.125 LITER 0.00	.S 13:41			lure of the			-	8		14 +		**	
Wait	0.00	13:43		p10 000	ure or m	Dom	u 01	r Coto.					划 :	
Air Blank	0.00	13:43									12			'
1	0.00 1.750 LITER	13:44	SIGNAT	URE OF O	PERATOR							190		
Air Blank	0.00	13:44												
1	0.07	13:44	NAME SMITH	.TO	SHUA							Test Re		
Air Blank	0.00	13:45	BADGE #			PERMIT#				┨ ┠	I	3rAC g/	210L	
			0000 AGENCY	0004	:35	000	001	L6837	, 	_		0.0	0 (
			BOAR	D OF	TESTS					L				
EFFECTIVE THIRTY (30) DAYS days or more, if you refused or fai (operating a CMV or CDL holder) Safety must receive your written r The approved hearing request Division, Department of Public Sa a fine of not more than Five Hund person's license, permit, or other	iled the state's testi OR a test result of equest for an admi form is available a fety, P.O. Box 114' red Dollars (\$500.0	is). In addition, your cor .04 or more while opera .nistrative hearing on the .th https://www.ok.gov/c/.5, OKC, OK 73136. It solo), or by both such fine	nmercial driv ting a comm e Departmer Ips/docume thall be a mis and imprisor	ring privileg ercial moto at approve ants/ReqHo sdemeanor ament, for a	ge in this state vor vehicle. To a sed form, within the earing.pdf. Sulr punishable by any person to a	vill be disq ppeal the fifteen (15 bmit your l imprisonm pply for a	qualified driver's i) days hearing nent for renewa	d if you refu- s license rev from the da g request fo r not less th	sed or fai vocation/o te you re rm in per an seven	iled the disqual ceived son at (7) da	state's testification, the notice of r 3600 N. M ys nor mo	st(s) ne Departme evocation of ILK Ave. or re than six (ent of Pub r disqualifi by mail to 6) months	ication. Legal or bv
SERVING OF NOTICE: 1 p			-			•		n to the n	erenn r	namer	l in Soct	ion 1 abov	<i>1</i> 0	
Date served: 05/06	12010	Signature of officer			ing the Notice						Badge #:	on rabo	vu.	
In accordance with Title 12	O.S. Section 42	6, "I state under pe										correct."	Se	ction 4
Date: 05/06/2019			- •					f arrestin	•	•				
NAME SMITH, JOSH					AGENCY BOARD	OF		STS						
	GENCY ADDRESS INTER – A						7	OODE		ONE	5425	0.4.6.0		



Certificate of Analysis – Wet Bath

09-26-2018 , 929-0/13

Date of Analysis Labelled target value (g/210L)

18029 0.0204

Average test result (g/210L)

BOT Technician Name and Signature

BRT Form 6.0

In accordance with BRT – 2.1.0, Maintenance, the above referenced simulator solution is suitable for use as an external reference in maintenance of the Intoxilyzer 8000.

CERTIFICATE OF ANALYSIS

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number 18020 of

Alcohol Reference Solution for Simulator were analyzed by
gas chromatography on January 11, 2018, using a Perkin Elmer Gas
Chromatograph Autosystem XL S/N: 610N9030209, and found to contain
0.0238% (w/vol) ethyl alcohol. The expiration date for this lot
number is January 9, 2020 at 11:59 PM.

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of 0.020 g/210L +/- 3%.

The alcohol and water used in this solution were free of test interfering substances.

Ted L. Pauley, President GUTH LABORATORIES, INC.

Operator's Signature

ALPHA INSTRUMENT

Model 8000

09/26/2018

Test

Air Blank

Cal Check

Air Blank

Tal Check

air Blank

Cal Check

Air Blank

al Check

IIr Blank

al Check

ir Blank

al Check

ir Blank

average atd Dev

al Check Stats

tel Std Deu(%)

Intoxilyzer - Alcohol Analyzer

g/210L

0.000

0.020

0.000

0.020

0.000

0.000

0.020

0.000

0.020

0.000

0.020

0.000

0.000

0.021

0.000

0.021

0.000

0.021

0.000

0.0204

0.0005

2.5312

SN 80-002591

15:36

Time

15:36 15:37

15:37

15:38 15:39

15:39

15:40

15:41

15:41

15:42

15:43

15:44

15:44

15:45

15:46

15:47

15:47

15:48

15:49

NIST Traceability:

Testing was conducted using Cerilliant Reference Standard lot number FN08031603 whose values are traceable to NIST.

All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.



Certificate of Analysis - Wet Bath 09-26-18 Labelled target value (g/210L) Date of Analysis 0.0396

BOT Technician Name and Signature

In accordance with BRT - 2.1.0, Maintenance, the above referenced simulator solution is Ashitable for use as external reference in maintenance of the intoxilyzer 8000.

CERTIFICATE OF

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number 17410 Alcohol Reference Solution for Simulator were analyzed by chromatography on December 7, 2017, using a Perkin Elmer Gas gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain 0.0481% (w/vol) ethyl alcohol. The expiration date for this lot number is December 6, 2019 at 11:59 PM.

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of 0.040 g/210L +/- 3%.

The alcohol and water used in this solution were free of test interfering substances.

Operator Signature

ALPHA INSTRUMENT Intoxilyzer - Alcohol Analyzer

g/210L

0.000

0.040

0.000

0.039

0.000

0.039

0,000

0.040

0.000

0.040

0,000

0.040

0,000

0.039

0.000

0.040

0.000

0.040

0.000

0.039

0.000

0.0396

0.0005

SN 80-002591

15:54

Time

15:54

15:55

15:56

15:56

15:57

15:58

15:58

15:59

15:59

16:00

16:01

16:01

16:02

16:03

16:03

16:04

16:04

16:05

16:06

16:06

Model 8000

09/26/2018

Test

Air Blank

Cal Check

Air Blank

Ruerage Std Dev

Cal Check Stats


Rel Std Dev(%)

Ted L. Pauley, President GUTH LABORATORIES, INC.

NIST Traceability:

Testing was conducted using Cerilliant Reference Standard lot number FN12181501 whose values are traceable to NIST.

All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.

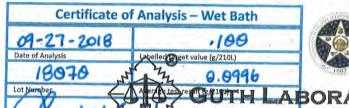

HARRISBURG, PA 17111-4511
 TELEPHONE: 717-564-5470

In accordance with BRT - 2.1.0, Maintenance, the above referenced simulator solution is suitable for use as an external reference in maintenance of the Intoxilyzer 8000.

CERTIFICATE OF ANALYSIS

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number 18070 Alcohol Reference Solution for Simulator were analyzed by chromatography on February 28, 2018, using a Perkin Elmer Gas gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain 0.1219% (w/vol) ethyl alcohol. The expiration date for this lot number is February 26, 2020 at 11:59 PM.

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of 0.100 g/210L +/- 3%.

The alcohol and water used in this solution were free of test interfering substances.

Ted L. Pauley, President GUTH LABORATORIES, INC.

ALPHA INSTRUMENT Intoxilyzer - Alcohol Analyzer

g/210L

0.100

0.000

0.100

0.000

0.099

0.000

0.100

0.000

0.100

0.000

0.100

0.000

0.099

0.000

0.099

0.000

0.099

0.000

0.100

0.0996

0.0005

0.5184

Signature

SN 80-002591

Time

08:28

08:29

08:29

08:30

08:31

08:31

08:32

08:33

08:33

08:34

08:34

08:35

08:36

08:36

08:37

08:38

08:38

08:39

18:39

08:40

Model 8000

09/27/2018

Test

Air Blank

Cal Check

Air Blank

Cal Check

Rir Blank

Cal Check


Air Blank Cal Check Stats

Average

Std Deu Rel Std Deu(%)

> NIST Traceability: Testing was conducted using Cerilliant Reference Standard lot number FN04271602 whose values are traceable to NIST.

> All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.

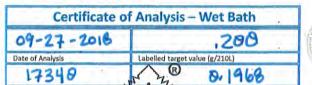

Lot Number

ALPHA INSTRUMENT

Model 8000

09/27/2018

Test

Air Blank

Cal Check

Pir Blank

:al Check

tir Blank

al Check

iir Blank

al Check

ifr Blank

Average Std Deu

al Check Stats

Rel Std Dev(%)

Intoxilyzer - Alcohol Analyzer

g/210L

0,000

0.196

0.000

0.196

0.000

0.197

0.000

0.197

0.000

0.197

0.000

0.197

0.000

0.197

0.000

0.197

0.000

0.197

0.000

0.197

0.000

0.1968

0.0004

0.2142

SN 80-002591

15:27

Time

15:27

15:28

15:29

15:29

15:30

15:31

15:31

15:32

15:33

15:33

15:34

15:35

15:35

15:36

15:36

15:37

15:38

15:38

15:39

15:40

BOT Technician Name and Signature

GUTH LABORATORIES, INC.

NORTH 67th STREET A HARRISBURG, PA 17111-4511 TELEPHONE: 717-564-5470

BRT Form 6.0

In accordance with BRT – 2.1.0, Maintenance, the above referenced simulator solution is suitable for use as an external reference in maintenance of the Intoxilyzer 8000.

CERTIFICATE OF ANALYSIS

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number 17340 of

Alcohol Reference Solution for Simulator were analyzed by
gas chromatography on October 10, 2017, using a Perkin Elmer Gas
Chromatograph Autosystem XL S/N: 610N9030209, and found to contain
0.2410% (w/vol) ethyl alcohol. The expiration date for this lot
number is October 9, 2019 at 11:59 PM.

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of 0.200 g/210L +/- 3%.

The alcohol and water used in this solution were free of test interfering substances.

Operator's Signature

Ted L. Pauley, President GUTH LABORATORIES, INC.

NIST Traceability:

Testing was conducted using Cerilliant Reference Standard lot number FN08101505 whose values are traceable to NIST.

All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.