

Service Overview

Board of Tests for Alcohol and Drug Influence
Enhancing public safety by ensuring the accuracy and scientific reliability
of blood and breath testing



3600 N. Martin Luther King Bldg #9
Oklahoma City, OK 73111
405-425-2460

BENCH CHECK DATE: 06-12-2017

SERIAL: 80-006119

DOCUMENTATION WITHIN SERVICE PACKET:

- BOT Certificate of Calibration and Operation
- Bench Check
- Bench Check Technician's Report
- BOT Dry Gas Certificate of Analysis
- Mock Subject Test
- BOT Certificates of Analysis on Certified Solutions (4)



**STATE OF OKLAHOMA
BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE**

Administrative Office at:
Board of Tests for Alcohol and Drug Influence
Post Office Box 36307
Oklahoma City, Oklahoma 73136-2307

Certificate of Calibration and Operation

This is to certify that the calibration of **INTOXILYZER 8000**, serial number 80-006119, manufactured by CMI Inc., a subsidiary of MPD, Inc. of Owensboro, Kentucky, was tested and found to be in compliance with National Highway Traffic Safety Administration Standard for Devices to Measure Breath Alcohol (F.R., No. 179 48705-48710 Sept. 17, 1993) and with International Recommendation OIML R 126 Edition 1998 (E). Calibration solutions are traceable to NIST standard material 1828.

In addition the operational procedures and functions of this Intoxilyzer have also been tested and found to be in compliance with the Oklahoma Board of Tests for Alcohol and Drug Influence - Rules and Regulations OAC Title 40: 25-1-2 and 30-1-3 and is hereby certified for use as an evidential breath testing device within this state.

CERTIFICATION TYPE	DATE
Manufacturer Calibration	12-05-13
Board of Tests Bench Check	12-13-13
Board of Tests Bench Check	10-07-15
Board of Tests Bench Check	06-12-17



Presentation of this form certifies that the Intoxilyzer, SN: 80-006119, functioned properly at the time of the breath test and hereby certifies the breath test result as valid.

Date:

Re: Request for Log of Test on Intoxilyzer 8000

Name:

Dear Sir:

The Log of Test information from a breath test conducted on a Intoxilyzer 8000 is printed in Section 2 of the Affidavit and is no longer retained as a separate form as with the Intoxilyzer 5000. The 8000 uses a nitrogen/ethanol dry gas for the external control test. The dry gas cylinder has the capacity to perform approximately 325 tests and does not require maintenance until its expiration date or can no longer provide the required pressure to perform a control test. The cylinder has an expiration date of two years from date filled at the manufacturer and is programmed in the 8000 upon installation, both dates are printed in Section 2 of the Affidavit. The Intoxilyzer is programmed to disable itself when the date of expiration is reached or the pressure monitor indicates 25 lbs psi.

Due to the reliability and stability of dry gas, maintenance on the 8000 is not required every 30 days or 25 tests, therefore maintenance is no longer performed by officers in the field. The Board of Tests for Alcohol and Drug Influence certifies each Intoxilyzer 8000 for field use by performing a bench check. The bench check is designed to check the operation and procedures programmed in the 8000 to ensure it is functioning properly for use in Oklahoma. A bench check is performed every time a new dry gas cylinder is installed or the instrument has returned from the manufacturer for repair. Once the instrument passes the bench check, all that is required to set it up for field use is to plug it in and turn it on and does not require a permit to set up. The Board of Tests has the ability to monitor each instrument once it is set up in the field.

Each Intoxilyzer has a database of valid permit numbers that it checks prior to allowing a breath test to be performed. The database is updated on a regular basis, if a permit number is invalid that officer will not be able to conduct a breath test on the 8000. If the 8000 prints an affidavit the breath test operator listed in Section 2 had a valid permit at the time of the test and the breath test is valid.

If you have any further questions please feel free to call me.

Sincerely,



J. Robert Blakeburn (OBA11059)
State Director of Tests for Alcohol
and Drug Influence

**STATE OF OKLAHOMA - DEPARTMENT OF PUBLIC SAFETY
CMI INTOXILYZER BENCH CHECK PROCEDURE REPORT**

BENCH CHECK DATE 06/12/2017	TIME PROCEDURE BEGAN 09:52	INSTRUMENT LOCATION 80-006119
INSTRUMENT TYPE Intoxilyzer	INSTRUMENT MODEL I-8000	INSTRUMENT SERIAL NUMBER 80-006119

TEST RESULTS

Diagnostics
PASS

Mouth Alcohol Test
PASS

RFI Test
PASS

Abort Test
PASS

Improper Sample Test
PASS

Interferent Detect Test
PASS

No Sample Given Test
PASS

Insufficient Sample Test
PASS

Air Blank Contamination Test
PASS

Range Exceeded Test
PASS


Duplicate Breath Sampling Results

Value	Test 1	Test 2	Control Test
.02	0.015	0.015	0.016
.04	0.034	0.033	0.034
.10	0.097	0.096	0.095
.20	0.196	0.192	0.193

PASS

CERTIFYING TECHNICIAN

I swear under penalty of perjury that in regards to the above listed instrument, I have complied with the Intoxilyzer Bench Checking Procedure approved by the State of Oklahoma.

SIGNATURE 	DATE 06/12/2017
NAME SMITH, JOSHUA	PERMIT # 0000016837



3600 N. Martin Luther King Bldg #9
 Oklahoma City, OK 73111
 405-425-2460

BENCH CHECK TECHNICIAN'S REPORT

Date: 06/12/17 Start Time: 1000 End Time: 1130

INSTRUMENT

Intoxilyzer Model 8000 Serial # 80-006119

GAS CANISTER

LOT # 09117080A2 # 816 EXP Date 6-5-2019


REFERENCE

Simulator Model			
Guth 2100	Guth 2100	Guth 2100	Guth 2100
Serial #			
DR3591	DR3753	DR3594	DR3754
Concentration			
.02	.04	.10	.20
LOT #			
16400	16320	15220	15240
Manufactured Date			
12-12-2016	10-21-2016	09-28-2015	11-11-2015
Expiration Date			
12-12-2018	10-21-2018	09-28-2017	11-11-2017
Solution Commission Date			
<u>5-30-17</u>	<u>5-30-17</u>	<u>6-8-17</u>	<u>6-8-17</u>

REASON FOR BENCH CHECK

- POST REPAIR CYLINDER REPLACEMENT
 TROUBLESHOOTING ROUTINE MAINTENANCE

COMMENTS


 Technician Signature / Permit # 16837
 Printed Name Josh Smith

Certificate of Analysis – Dry Gas

06/02/2017	.080
Date of Analysis	Labelled target value (g/210L)
016	.0793
Cylinder #	Average test result
 BOT Technician Name and Signature	



ILMO
specialty gases

BRT Form 5.0

ISO/IEC 17025:2005 Accredited Laboratory

In accordance with BOT Policy and Procedure Statement BRT-2.1.0, Maintenance, the above referenced dry gas cylinder is suitable for use as an external reference with the Intoxilyzer 8000.

ALPHA INSTRUMENT 2
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-004923
06/02/2017 13:24

Test	g/210L	Time
Air Blank	0.000	13:25
Cal Check	0.080	13:25
Air Blank	0.000	13:26
Cal Check	0.080	13:26
Air Blank	0.000	13:27
Cal Check	0.080	13:27
Air Blank	0.000	13:27
Cal Check	0.080	13:28
Air Blank	0.000	13:28
Cal Check	0.080	13:29
Air Blank	0.000	13:29
Cal Check	0.079	13:29
Air Blank	0.000	13:30
Cal Check	0.080	13:30
Air Blank	0.000	13:31
Cal Check	0.079	13:31
Air Blank	0.000	13:32
Cal Check	0.080	13:32
Air Blank	0.000	13:32
Cal Check	0.080	13:33
Air Blank	0.000	13:33
Cal Check Stats		
Average	0.0798	
Std Dev	0.0004	
Rel Std Dev(%)	0.5283	

Certificate of Analysis


Certificate ID: 10144 ✓
 Part #: BAC105L080T
 Cylinder Size: 105L
 Lot Number: 09117080A2 ✓
 Expiration: 6/5/2019 ✓

REC'D
06/01/17


0.080 BAC (For the calibration of instruments used to determine breath alcohol concentration)

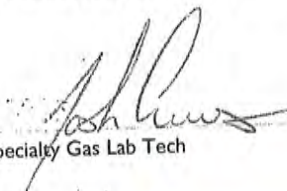
Contents: 105 Liters @ 1000 psig 70°F (21°C)

Component:	Concentration:	Accuracy:	Method:
Ethanol	208 ppm	+/- 0.002 or 2%	NDIR
Nitrogen	Balance	BAC whichever is greater	


Operator's Signature 16837

*NIST Standard Reference Material
Cylinder No. CC 274507 / Job No. 09160309
Certified 362.2 µmol/mol Ethanol in Nitrogen
for ILMO Products Co., Jacksonville, IL

Store in dry area, away from sources of heat, ignition and direct sunlight. Do not allow storage area to exceed 52 °C (125 °F).


Specialty Gas Lab Tech

04/25/17
Date



Distributed by:

CMI Inc.
316 East Ninth Street
Owensboro, KY 42303
Phone 866-835-0690
www.alcoholtest.com

OFFICER'S AFFIDAVIT AND NOTICE OF REVOCATION/DISQUALIFICATION

Section 1

ARREST DATE 06/12/2017		ARREST TIME 1000		COUNTY # 55	CITY # 55	CITATION # MOCK	
ARREST LOCATION MOCK				CITY MOCK		COUNTY OKLAHOMA	
SUBJECT NAME MOCK, SUBJECT TEST				DATE OF BIRTH 08/09/1961	SEX M	RACE W	HEIGHT 510
ADDRESS 3600 MARTIN LUTHER KING A				CITY OKLAHOMA CITY		STATE OK	ZIP CODE 73133
DRIVER LICENSE # T083454871		EXPIRATION DATE 06/30/2019	STATE OK	CLASS A	RESTRICTIONS NONE	ENDORSEMENTS NONE	
VEHICLE MAKE MOCK		MODEL MOCK	TAG # MOCK		STATE OK	CMV N	PLACARD N

On the above date, time, and location, the above named person was arrested, and I had reasonable grounds to believe the person was driving or in actual physical control of a motor vehicle upon a public road, highway, street, turnpike, public parking lot, or other public place, or upon a private road, street, alley or lane which provides access to one or more single or multi-family dwelling within this state while under the influence of alcohol and/or other intoxicating substances as prohibited by law.

(Describe driving behavior or circumstances):

MOCK SUBJECT TEST. ROUTINE MAINTENANCE, NEW GAS CANISTER INSTALL.

Describe person's condition (odor, actions, etc.):

NONE

THE PERSON WAS INFORMED OF THE IMPLIED CONSENT TEST REQUEST

Section 2

BREATH-ALCOHOL ANALYSIS RECORD, REPORT OF THE PERSON NAMED IN SECTION 1 ABOVE, and LOG OF TEST AND MAINTENANCE RECORD.

EQUIPMENT I-8000		SERIAL NUMBER 80-006119	SPECIALIST BOARD OF TESTS																																					
STANDARD Dry		TARGET VALUE 0.080	PERMIT NUMBER 16837																																					
MANUFACTURER ILMO		LOT # 09117080A2	DATE CYLINDER INSTALLED 06/12/2017	DATE CYLINDER EXPIRES 06/05/2019																																				
OBSERVATION START 1025	TEST DATE 06/12/2017	MOUTHPIECE G	<p style="text-align: center;">A 15-minute (or greater) deprivation/observation period was observed prior to and continuing to the time of test administration.</p> <p style="text-align: center;">I administered said Breath-Alcohol Test in compliance with the applicable operating procedure of the Board of Tests.</p> <p style="text-align: center;">_____ SIGNATURE OF OPERATOR</p>																																					
<table border="1" style="width: 100%; border-collapse: collapse;"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.00</td><td>11:05</td></tr> <tr><td>Subject Test</td><td>0.00</td><td>11:06</td></tr> <tr><td> Breath Volume</td><td>2.324 LITERS</td><td></td></tr> <tr><td>Air Blank</td><td>0.00</td><td>11:06</td></tr> <tr><td>Wait</td><td></td><td>11:08</td></tr> <tr><td>Air Blank</td><td>0.00</td><td>11:08</td></tr> <tr><td>Subject Test</td><td>0.00</td><td>11:09</td></tr> <tr><td> Breath Volume</td><td>1.753 LITERS</td><td></td></tr> <tr><td>Air Blank</td><td>0.00</td><td>11:09</td></tr> <tr><td>Cal Check</td><td>0.07</td><td>11:09</td></tr> <tr><td>Air Blank</td><td>0.00</td><td>11:10</td></tr> </tbody> </table>					Test	g/210L	Time	Air Blank	0.00	11:05	Subject Test	0.00	11:06	Breath Volume	2.324 LITERS		Air Blank	0.00	11:06	Wait		11:08	Air Blank	0.00	11:08	Subject Test	0.00	11:09	Breath Volume	1.753 LITERS		Air Blank	0.00	11:09	Cal Check	0.07	11:09	Air Blank	0.00	11:10
Test	g/210L	Time																																						
Air Blank	0.00	11:05																																						
Subject Test	0.00	11:06																																						
Breath Volume	2.324 LITERS																																							
Air Blank	0.00	11:06																																						
Wait		11:08																																						
Air Blank	0.00	11:08																																						
Subject Test	0.00	11:09																																						
Breath Volume	1.753 LITERS																																							
Air Blank	0.00	11:09																																						
Cal Check	0.07	11:09																																						
Air Blank	0.00	11:10																																						
RESULT TABLE			<div style="border: 1px solid black; padding: 5px; width: fit-content; margin: auto;"> <p>Test Result BrAC g/210L 0.00</p> </div>																																					
<table border="1" style="width: 100%; border-collapse: collapse;"> <tr> <td>NAME SMITH, JOSHUA</td> <td>PERMIT # 0000016837</td> </tr> <tr> <td>BADGE # 0000000435</td> <td>AGENCY BOARD OF TESTS</td> </tr> </table>					NAME SMITH, JOSHUA	PERMIT # 0000016837	BADGE # 0000000435	AGENCY BOARD OF TESTS																																
NAME SMITH, JOSHUA	PERMIT # 0000016837																																							
BADGE # 0000000435	AGENCY BOARD OF TESTS																																							

Section 3

EFFECTIVE THIRTY (30) DAYS FROM DATE SERVED. Notice of Revocation is given that your driver's license is revoked or denied for a period of one hundred eighty (180) days or more, if you refused or failed the state's test(s). In addition, your commercial driving privilege in this state will be disqualified if you refused or failed the state's test(s) (operating a CMV or CDL holder) OR a test result of .04 or more while operating a commercial motor vehicle. To appeal the driver's license revocation/disqualification, the Department of Public Safety must receive your written request for an Administrative Hearing **within fifteen (15) days from the date you received notice of revocation or disqualification.** Submit your **written** request in person at 3600 N. MLK Ave. or by mail to Legal Division, Department of Public Safety, P.O. Box 11415, OKC, OK, 73136

SERVING OF NOTICE: I personally hand-delivered a copy of the above containing the Notice of Revocation to the person named in Section 1 above.

Date served: **06/12/2017** Signature of officer: _____ Badge #: _____

In accordance with Title 12 O.S. Section 426, "I state under **penalty of perjury** under the laws of Oklahoma that the foregoing is true and correct." Section 4

Date: **06/12/2017** Place (location when signed): _____ Signature of arresting officer: _____

NAME SMITH, JOSHUA		AGENCY BOARD OF TESTS	
BADGE 0000000435	AGENCY ADDRESS INTER-AGENCY	ZIP CODE 00000	PHONE "4054252460"

Section 5

OTHER WITNESSES:

1	NAME	TITLE	ADDRESS	PHONE
2	NAME	TITLE	ADDRESS	PHONE

No Temporary License Issued:
Test Below Legal Limit

OFFICER'S AFFIDAVIT AND NOTICE OF REVOCATION/DISQUALIFICATION

Section 1

ARREST DATE 06/12/2017		ARREST TIME 1000		COUNTY # 55	CITY # 55	CITATION # MOCK	
ARREST LOCATION MOCK				CITY MOCK		COUNTY OKLAHOMA	
SUBJECT NAME MOCK, SUBJECT TEST				DATE OF BIRTH 08/09/1961	SEX M	RACE W	HEIGHT 510
ADDRESS 3600 MARTIN LUTHER KING A				CITY OKLAHOMA CITY		STATE OK	ZIP CODE 73133
DRIVER LICENSE # T083454871				EXPIRATION DATE 06/30/2019	STATE OK	CLASS A	RESTRICTIONS NONE
VEHICLE MAKE MOCK				MODEL MOCK	TAG # MOCK	STATE OK	CMV N
						PLACARD N	ENDORSEMENTS NONE

On the above date, time, and location, the above named person was arrested, and I had reasonable grounds to believe the person was driving or in actual physical control of a motor vehicle upon a public road, highway, street, turnpike, public parking lot, or other public place, or upon a private road, street, alley or lane which provides access to one or more single or multi-family dwelling within this state while under the influence of alcohol and/or other intoxicating substances as prohibited by law.

(Describe driving behavior or circumstances):

MOCK SUBJECT TEST. ROUTINE MAINTENANCE, NEW GAS CANISTER INSTALL.

Describe person's condition (odor, actions, etc.):

NONE

THE PERSON WAS INFORMED OF THE IMPLIED CONSENT TEST REQUEST

BREATH-ALCOHOL ANALYSIS RECORD, REPORT OF THE PERSON NAMED IN SECTION 1 ABOVE, and LOG OF TEST AND MAINTENANCE RECORD.

Section 2

EQUIPMENT I-8000		SERIAL NUMBER 80-006119	SPECIALIST BOARD OF TESTS																																					
STANDARD Dry		TARGET VALUE 0.080	PERMIT NUMBER 16837																																					
MANUFACTURER ILMO		LOT # 09117080A2	DATE CYLINDER INSTALLED 06/12/2017	DATE CYLINDER EXPIRES 06/05/2019																																				
OBSERVATION START 1025	TEST DATE 06/12/2017	MOUTHPIECE G		A 15-minute (or greater) deprivation/observation period was observed prior to and continuing to the time of test administration.																																				
<table border="1" style="width: 100%; border-collapse: collapse;"> <thead> <tr> <th>Test</th> <th>g/210L</th> <th>Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.00</td><td>11:05</td></tr> <tr><td>Subject Test</td><td>0.00</td><td>11:06</td></tr> <tr><td> Breath Volume</td><td>2.324 LITERS</td><td></td></tr> <tr><td>Air Blank</td><td>0.00</td><td>11:06</td></tr> <tr><td>Wait</td><td></td><td>11:08</td></tr> <tr><td>Air Blank</td><td>0.00</td><td>11:08</td></tr> <tr><td>Subject Test</td><td>0.00</td><td>11:09</td></tr> <tr><td> Breath Volume</td><td>1.753 LITERS</td><td></td></tr> <tr><td>Air Blank</td><td>0.00</td><td>11:09</td></tr> <tr><td>Cal Check</td><td>0.07</td><td>11:09</td></tr> <tr><td>Air Blank</td><td>0.00</td><td>11:10</td></tr> </tbody> </table>					Test	g/210L	Time	Air Blank	0.00	11:05	Subject Test	0.00	11:06	Breath Volume	2.324 LITERS		Air Blank	0.00	11:06	Wait		11:08	Air Blank	0.00	11:08	Subject Test	0.00	11:09	Breath Volume	1.753 LITERS		Air Blank	0.00	11:09	Cal Check	0.07	11:09	Air Blank	0.00	11:10
Test	g/210L	Time																																						
Air Blank	0.00	11:05																																						
Subject Test	0.00	11:06																																						
Breath Volume	2.324 LITERS																																							
Air Blank	0.00	11:06																																						
Wait		11:08																																						
Air Blank	0.00	11:08																																						
Subject Test	0.00	11:09																																						
Breath Volume	1.753 LITERS																																							
Air Blank	0.00	11:09																																						
Cal Check	0.07	11:09																																						
Air Blank	0.00	11:10																																						
SIGNATURE OF OPERATOR _____ NAME SMITH, JOSHUA BADGE # 0000000435 PERMIT # 0000016837 AGENCY BOARD OF TESTS																																								



Test Result BrAC g/210L <h2 style="margin: 0;">0.00</h2>
--

Section 3

EFFECTIVE THIRTY (30) DAYS FROM DATE SERVED. Notice of Revocation is given that your driver's license is revoked or denied for a period of one hundred eighty (180) days or more, if you refused or failed the state's test(s). In addition, your commercial driving privilege in this state will be disqualified if you refused or failed the state's test(s) (operating a CMV or CDL holder) OR a test result of .04 or more while operating a commercial motor vehicle. To appeal the driver's license revocation/disqualification, the Department of Public Safety must receive your written request for an Administrative Hearing **within fifteen (15) days from the date you received notice of revocation or disqualification.** Submit your **written** request in person at 3600 N. MLK Ave. or by mail to Legal Division, Department of Public Safety, P.O. Box 11415, OKC, OK, 73136

SERVING OF NOTICE: I personally hand-delivered a copy of the above containing the Notice of Revocation to the person named in Section 1 above.

Date served: **06/12/2017** Signature of officer: _____ Badge #: _____

In accordance with Title 12 O.S. Section 426, "I state under penalty of perjury under the laws of Oklahoma that the foregoing is true and correct."

Date: **06/12/2017** Place (location when signed): _____ Signature of arresting officer: _____

NAME SMITH, JOSHUA		AGENCY BOARD OF TESTS	
BADGE 0000000435	AGENCY ADDRESS INTER-AGENCY	ZIP CODE 00000	PHONE "4054252460"

Section 4

No Temporary License Issued:
Test Below Legal Limit

ARRESTEE'S COPY
BOT Form AFF02 07/2008

Certificate of Analysis - Wet Bath	
05-05-17	.020
Date of Analysis	Labelled target value (g/210L)
16400	.0195
Lot Number	Average test result (g/210L)
Jost Smirk <i>[Signature]</i>	
BOT Technician Name and Signature	



ORIES, INC.
 LG, PA 17111-4511 ☎ TELEPHONE: 717-564-5470

BRT Form 6.0

In accordance with BRT - 2.1.0, Maintenance, the above referenced simulator solution is suitable for use as an external reference in maintenance of the Intoxilyzer 8000.

CERTIFICATE OF ANALYSIS

INTOXILYZER 8000
 Instrument Initialization
 8133.12
 14:48 05/05/2017

Certified Alcohol Reference Solution for Simulator

*Rec'd
 4/01/2016
 #16837*

ALPHA INSTRUMENT 2
 Intoxilyzer - Alcohol Analyzer
 Model 8000
 05/05/2017

SN 80-004923
 14:56

Random Samples of Lot Number 16400 of Alcohol Reference Solution for Simulator were analyzed by gas chromatography on December 13, 2016, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain

Test	g/210L	Time
Air Blank	0.000	14:57
Cal Check	0.020	14:58
Air Blank	0.000	14:58
Cal Check	0.020	14:59
Air Blank	0.000	14:59
Cal Check	0.019	15:00
Air Blank	0.000	15:01
Cal Check	0.020	15:01
Air Blank	0.000	15:02
Cal Check	0.020	15:03
Air Blank	0.000	15:03
Cal Check	0.020	15:04
Air Blank	0.000	15:04
Cal Check	0.019	15:05
Air Blank	0.000	15:06
Cal Check	0.019	15:06
Air Blank	0.000	15:07
Cal Check	0.019	15:08
Air Blank	0.000	15:08
Cal Check	0.019	15:09
Air Blank	0.000	15:09
Cal Check Stats		
Average	0.0195	
Std Dev	0.0005	
Rel Std Dev(%)	2.7026	

0.0241% (w/vol) ethyl alcohol. The expiration date for this lot number is December 12, 2018 at 11:59 PM.


When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of 0.020 g/210L +/- 3%.

The alcohol and water used in this solution were free of test interfering substances.

Ted L. Pauley
 Ted L. Pauley, President
 GUTH LABORATORIES, INC.

#16837

*NIST Traceability:
 Testing was conducted using Cerilliant Reference Standard lot number FN11241503 whose values are traceable to NIST.
 All balances are calibrated annually by an outside agency using NIST traceable weights.
 Calibration verification is done prior to each use utilizing NIST traceable weights.*

Certificate of Analysis - Wet Bath	
Date of Analysis	3-8-2017
Labelled target value (g/210L)	.040
Lot Number	16320
Average test result (g/210L)	.0409
 BOT Technician Name and Signature	




GUTH LABORATORIES, INC.
 BURG, PA 17111-4511 • TELEPHONE: 717-564-5470

BRT Form 6.0

In accordance with BRT - 2.1.0, Maintenance, the above referenced simulator solution is suitable for use as an external reference in maintenance of the Intoxilyzer 8000.

CERTIFICATE OF ANALYSIS

Certified Alcohol Reference Solution for Simulator

REC'D
 03/01/17

 #16837

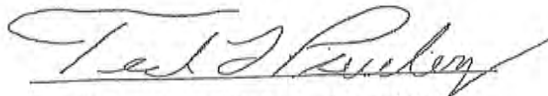
Random Samples of Lot Number 16320 of Alcohol Reference Solution for Simulator were analyzed by gas chromatography on October 25, 2016, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain 0.0495% (w/vol) ethyl alcohol. The expiration date for this lot number is **October 21, 2018** at 11:59 PM.

ALPHA INSTRUMENT 2
 Intoxilyzer - Alcohol Analyzer
 Model 8000
 03/08/2017
 SN 80-004923
 11:16

Test	g/210L	Time
Air Blank	0.000	11:16
Cal Check	0.041	11:17
Air Blank	0.000	11:17
Cal Check	0.041	11:18
Air Blank	0.000	11:19
Cal Check	0.040	11:19
Air Blank	0.000	11:20
Cal Check	0.041	11:21
Air Blank	0.000	11:21
Cal Check	0.041	11:22
Air Blank	0.000	11:22
Cal Check	0.041	11:23
Air Blank	0.000	11:24
Cal Check	0.041	11:24
Air Blank	0.000	11:25
Cal Check	0.040	11:26
Air Blank	0.000	11:26
Cal Check	0.042	11:27
Air Blank	0.000	11:27
Cal Check	0.041	11:28
Air Blank	0.000	11:29
Cal Check Stats		
Average	0.0409	
Std Dev	0.0005	
Rel Std Dev(%)	1.3878	

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of 0.040 g/210L +/- 3%.

The alcohol and water used in this solution were free of test interfering substances.


 Ted L. Pauley, President
 GUTH LABORATORIES, INC.


 #16837

NIST Traceability:
 Testing was conducted using Cerilliant Reference Standard lot number FNI218501 whose values are traceable to NIST.
 All balances are calibrated annually by an outside agency using NIST traceable weights.
 Calibration verification is done prior to each use utilizing NIST traceable weights.



GUTH LABORATORIES, INC.
 590 NORTH 6TH STREET • HARRISBURG, PA 17111-4511 • TELEPHONE: 717-564-5470

CERTIFICATE OF ANALYSIS

*REC'D
 11-4-2015
 [Signature]*

Certificate of Analysis - Wet Bath	
Date of Analysis	12-12-2016
Labelled target value (g/210L)	.100
Lot Number	15220
Average test result (g/210L)	0.0999
BOTT Technician Name and Signature JOSH SMITH <i>[Signature]</i> 16837	

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number 15220 of Alcohol Reference Solution for Simulator were analyzed by gas chromatography on **September 30, 2015**, using a Perkin Elmer Gas Chromatograph. Autosystem XL S/N: 610N9030209, and found to contain **0.1214%** (w/vol) ethyl alcohol. The expiration date for this lot number is **September 28, 2017** at 11:59 PM.

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of **0.100 g/210L +/- 3%**.

The alcohol and water used in this solution were free of test interfering substances.

[Signature]
 Ted L. Pauley, President
 GUTH LABORATORIES, INC.

Intoxilyzer - Alcohol Analyzer
 Model 8000 SN 80-002591
 12/12/2016 14:04

Test	g/210L	Time
Air Blank	0.000	14:04
Cal Check	0.100	14:05
Air Blank	0.000	14:05
Cal Check	0.100	14:06
Air Blank	0.000	14:07
Cal Check	0.100	14:07
Air Blank	0.000	14:08
Cal Check	0.100	14:09
Air Blank	0.000	14:09
Cal Check	0.099	14:10
Air Blank	0.000	14:10
Cal Check	0.100	14:11
Air Blank	0.000	14:12
Cal Check	0.100	14:12
Air Blank	0.000	14:13
Cal Check	0.100	14:14
Air Blank	0.000	14:14
Cal Check	0.100	14:15
Air Blank	0.000	14:16
Cal Check	0.100	14:16
Air Blank	0.000	14:17
Cal Check Stats		
Average	0.0999	
Std. Dev	0.0003	
Rel. Std. Dev (%)	0.3165	

NIST Traceability:
 Testing was conducted using Cerilliant Reference Standard lot number FN08051301 whose values are traceable to NIST.
 All balances are calibrated annually by an outside agency using NIST traceable weights.
 Calibration verification is done prior to each use utilizing NIST traceable weights.

212-16



BRT Form 5.0

In accordance with BRT - 2.1.0 Maintenance, the above referenced simulator solution is suitable for use as an external reference in maintenance of the Intoxilyzer 8000.



GUTH LABORATORIES, INC.

590 NORTH 67th STREET • HARRISBURG, PA 17111-4511 • TELEPHONE: 717-564-5470

REC'D
11-4-2015
J.P. # 16837

CERTIFICATE OF ANALYSIS

Certificate of Analysis - Wet Bath	
Date of Analysis 12-12-16	Labelled target value (g/210L) .200
Lot Number 15240	Average test result (g/210L) .2006
BOTT Technician Name and Signature Ted L. Pauley #16837	

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number 15240 of Alcohol Reference Solution for Simulator were analyzed by gas chromatography on November 12, 2015, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain 0.2418% (w/vol) ethyl alcohol. The expiration date for this lot number is November 11, 2017 at 11:59 PM.

ALPHA INSTRUMENT
Intoxilyzer - Alcohol Analyzer
Model 8000
12/12/2016

SN 80-002591
13:42

Test	g/210L	Time
Air Blank	0.000	13:42
Cal Check	0.202	13:43
Air Blank	0.000	13:43
Cal Check	0.201	13:44
Air Blank	0.000	13:45
Cal Check	0.201	13:45
Air Blank	0.000	13:46
Cal Check	0.201	13:47
Air Blank	0.000	13:47
Cal Check	0.200	13:48
Air Blank	0.000	13:49
Cal Check	0.200	13:49
Air Blank	0.000	13:50
Cal Check	0.201	13:51
Air Blank	0.000	13:51
Cal Check	0.200	13:52
Air Blank	0.000	13:52
Cal Check	0.199	13:53
Air Blank	0.000	13:54
Cal Check	0.201	13:54
Air Blank	0.000	13:55
Cal Check Stats		
Average	0.2006	
Std Dev	0.0008	
Rel Std Dev(%)	0.4203	

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of 0.200 g/210L +/- 3%.

The alcohol and water used in this solution were free of test interfering substances.

Ted L. Pauley
Ted L. Pauley, President
GUTH LABORATORIES, INC.

NIST Traceability:

Testing was conducted using Cerilliant Reference Standard lot number FN03211401 whose values are traceable to NIST.

All balances are calibrated annually by an outside agency using NIST traceable weights.

Calibration verification is done prior to each use utilizing NIST traceable weights.

12-12-16